



RNS® System Programming Manual

**For the RNS® Tablet Model 5000 and the
Patient Data Management System (PDMS) Model 4340**

R_XOnly

© 2021 NeuroPace, Inc.

DN 1018685 Rev 3
Rev Date: 2021-05

ABOUT THIS MANUAL

This manual includes instructions for use for:

- The RNS® Tablet model 5000 with software version 1.8
- The NeuroPace® Patient Data Management System model 4340

Note: The RNS® Tablet model 5000 is compatible with the RNS® Neurostimulator model RNS-300M and model RNS-320.

Note: The term “programmer” as used in this manual is a generic term that refers to either the RNS® Tablet or the NeuroPace® Programmer (laptop computer).

Note: Images in this manual are representative and may vary in detail from what a particular user experiences.

FCC INFORMATION

The following is communications regulation information on the neurostimulator models RNS-300M and RNS-320, and wand model W-02.

- Neurostimulator model RNS-300M FCC ID: WBWRF300
- Neurostimulator model RNS-320 FCC ID: WBWRF320
- Wand FCC ID: WBW5200 or WBW902

These components comply with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) These devices may not cause harmful interference, and (2) these devices must accept any interference received, including interference that may cause undesired operation.

IMPORTANT: Changes or modifications to these components not expressly approved by NeuroPace, Inc. could void the FCC Certification, and negate your authority to operate them.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

SYMBOLS

Explanation of symbols on product or package labeling



Caution



Do Not Resterilize



MR Conditional



MR Unsafe



Rx Only

Prescription Only



Non-Pyrogenic



Single Use



Sterilized Using Ethylene Oxide



Temperature limits during use

Use



Temperature limits during storage or transport

Storage




Manufacturer



Serial number



Type BF applied part

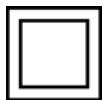


Class II electrical protection

5 V, 0.5 A

Direct current (DC) 5 V (volts), 0.5 A (amperes)

IP22

Ingress protection ratings: Level 2 for solid objects, which means testing confirms the device enclosure prevents ingress (entry) of items greater than 12.5 mm (~1/2 inch), such as fingers or similar objects; Level 2 for liquids, which means testing confirms vertically dripping water shall have no harmful effect when the enclosure is tilted at an angle up to 15° from its normal position.

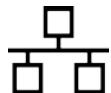

Ethernet Connection (Network Connection)

Proposition 65, a State of California voter initiative, requires the following notice:



WARNING: This product can expose you to chemicals including ethylene oxide, which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.

TABLE OF CONTENTS

ABOUT THIS MANUAL.....	II
FCC INFORMATION.....	III
SYMBOLS	IV
INTRODUCTION	4
DEVICE DESCRIPTION	4
RNS® TABLET	4
NEUROPACE® PATIENT DATA MANAGEMENT SYSTEM	4
PRODUCTS AND ACCESSORIES	4
CONTACTING NEUROPACE	5
TYPOGRAPHIC CONVENTIONS	5
WARNINGS AND CAUTIONS – PROGRAMMER	6
STERILIZATION AND STORAGE	8
Product Storage	8
Explant and Disposal.....	8
Return / Disposal of Programmer and Wand.....	9
Wand Maintenance, Cleaning and Sterilization.....	9
Electrostatic Discharge (ESD) / Static Electricity.....	9
RNS® TABLET INSTRUCTIONS.....	10
POWER ON THE TABLET.....	10
CONNECT AND LOGIN OVERVIEW	10
FUNCTIONALITY WITH AND WITHOUT A PDMS CONNECTION	11
CONNECT TO THE INTERNET FOR PDMS ACCESS.....	11
Mobile Broadband Works Automatically.....	11
Steps to Connect to a Wi-Fi Network	12
If You Lose Network Connectivity.....	13
LOG IN TO THE RNS® TABLET	14
NEUROPACE UPDATE CENTER	15
Update Types.....	15
Update Options	15
Update Process and Power Requirements.....	16
LOG IN TO THE PDMS	16
INTERROGATE THE RNS® NEUROSTIMULATOR	18
DATA RETENTION AND ACCESS ON THE NEUROSTIMULATOR, PROGRAMMER AND PDMS .	19
Programmed Settings.....	19
Storage of ECoGs and Neurostimulator Activity (Diagnostics)	19
HOME SCREEN.....	21
NEUROSTIM INFO SCREEN.....	23
MRI Mode for RNS Neurostimulator Model RNS-320	25

Electrode Impedance and Battery Voltage	26
ACTIVITY SCREEN.....	28
Activity Histogram Viewing Options.....	29
View Neurostim Settings	30
View Event List.....	33
Types of Reports.....	33
Advanced Views.....	34
ECoG LIBRARY SCREEN	34
LIVE ECoGs SCREEN.....	37
CREATE THE RECORDING MONTAGE	38
Identifying Ports, Leads and Electrodes on the Neurostimulator.....	38
Creating the Montage.....	39
Recommended Initial Montage Settings.....	41
Program the Neurostimulator with the New Settings.....	43
SET UP ECoG CAPTURE.....	44
Selecting ECoG Capture Settings	44
Recommended Initial ECoG Capture Settings	47
Program the Neurostimulator with the New Settings	49
RECOMMENDATION FOR REVIEW OF INITIAL ECoG CAPTURE SETTINGS.....	49
CONFIGURE PATTERN DETECTION: DEFINE AND MODIFY DETECTION SETTINGS	50
Phase 1: Select a Default Detection Set upon First Use	51
Phase 2: Define a Pattern Detection Set.....	54
Phase 3: Customize Pattern Detection.....	58
DETECTION TOOLS AND DETECTORS.....	62
Detection Tools and Detectors	62
Power Change Detector	63
Rhythmic Activity Detector.....	65
Spike Activity Detector.....	66
Bandpass Detection Tool	67
Area Detection Tool.....	68
Advanced Controls	69
SHARE AND/OR DOWNLOAD A DETECTION SET	70
Steps to Share a Detection Set	70
Steps to Download a Detection Set to the Tablet.....	72
OVERVIEW OF RESPONSIVE THERAPY	75
CONFIGURE RESPONSIVE THERAPY	76
Optional Advanced Settings	80
Test Stimulation Before Enabling	83
Enable Therapy: Review & Program	84
Recommended Initial Responsive Therapy Settings	85
Recommended Therapy Modifications After Observing Clinical Response.....	86
EXPORT AND IMPORT DEVICE SETTINGS FOR NEUROSTIMULATOR REPLACEMENT	87
Transfer Compatibility.....	87
Export Neurostimulator Settings.....	87
Import Neurostimulator Settings	88

PATIENT LIST ON TABLET AND THE PDMS ON A PERSONAL COMPUTER	90
PDMS ACCESS ON A PERSONAL COMPUTER.....	91
PDMS PROFILE AND SETTINGS	91
PATIENT & NEUROSTIM INFO	94
NEUROSTIM SETTINGS	94
ACTIVITY.....	96
ECoG LIBRARY	96
LIVE ECoGs.....	96
PATIENT FOLLOW-UP ACTIVITIES.....	97
TROUBLESHOOTING	98
Network Connectivity Problems.....	98
Abnormal Lead Impedance (greater than 3500 Ohms or less than 250 Ohms)	98
Insufficient Charge	101
Noise, Artifacts, Poor Signal Displayed, or No Signal Displayed in Live ECoG.....	101
Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet.....	102
Tablet Does Not Turn On	103
Tablet Freezes or Does Not Turn Off	103
Tablet Shuts Down	103
Impedance Measurement Was Rejected / Test Request Was Rejected	103
SPECIFICATIONS AND CHARACTERISTICS	104
RNS® TABLET	104
WAND (MODEL W-02)	104
TABLET WIRELESS	105
RNS® SYSTEM WIRELESS.....	105
ELECTROMAGNETIC EMISSIONS AND IMMUNITY.....	105
Hospital or Medical Environments	106
Home, Work or Public Environments.....	106
Guidance and Manufacturer's Declaration	107
Emissions and Immunity Information.....	108
GLOSSARY.....	112
INDEX.....	117

INTRODUCTION

DEVICE DESCRIPTION

- RNS® Tablet model 5000 with wand model W-02 and software version 1.8
- NeuroPace® Patient Data Management System model 4340

RNS® TABLET

The RNS® Tablet includes a tablet computer that runs proprietary NeuroPace software, and uses custom telemetry components to communicate with the RNS® Neurostimulator by holding the wand over it. The wand is applied to the patient (an applied part). Clinicians can use the tablet to program how the neurostimulator operates. Settings include, but are not limited to, the implanted system configuration, detection settings adapted to the patient's ECoG patterns, settings for ECoG record storage, and settings of the therapies. Clinicians can also use the tablet to review previously retrieved neurostimulator activity information, perform detection analysis, and communicate with the NeuroPace® Patient Data Management System (PDMS) via the Internet.

Note: The term “programmer” as used in this manual is a generic term that refers to either the RNS® Tablet or the NeuroPace® Programmer (laptop computer).

NEUROPACE® PATIENT DATA MANAGEMENT SYSTEM

The Patient Data Management System (PDMS) maintains patient and product data obtained from the programmer and remote monitor. Authorized users may access the PDMS via the Internet using a personal computer. An electronic signature in the form of a username and password are required for user authentication.

PRODUCTS AND ACCESSORIES

Note: Manuals and other literature are not listed below.

RNS® Tablet Kit Model 5001	
Item	Model Number
RNS® Tablet (with attached case)	5000
Stylus	
AC Adapter	
Carrying Bag	
Used With	
Wand	W-02

This manual contains information applicable to the RNS® Tablet, which is used to program how the RNS® Neurostimulator operates. It includes information about the secure Patient Data Management System (PDMS) database, which is accessible using the tablet or a web browser on a personal computer.

Refer to the RNS® System manual for additional device descriptions, device specifications, and indications for use, contraindications, warnings, cautions, and instructions for use.

All NeuroPace® manuals are available at www.NeuroPace.com or by contacting NeuroPace, Inc. (see **Contacting NeuroPace** on page 5).

Refer to the clinical summary booklet for information on the clinical study results of the RNS® System and adverse event data.

CONTACTING NEUROPACE

All questions or concerns regarding the NeuroPace® RNS® System should be forwarded to:

NeuroPace, Inc.
455 N. Bernardo Ave.
Mountain View, CA 94043

Customer Support: 866-726-3876 (Toll Free in the United States)

Website: www.NeuroPace.com

TYPOGRAPHIC CONVENTIONS

This manual uses the following typographic conventions.

WARNING: WARNING TITLE

Warnings alert you to serious adverse events and potential safety hazards and situations that may cause injury.

Caution: Caution Title

Cautions alert you to exercise special care.

Note: Notes provide additional information.

1. Numbered lists are used to identify a sequence of steps.

This format is used to identify figure titles and descriptions

BOLD SMALL CAPS indicate on-screen text, like the names of screens, buttons and fields.

Note: Software version numbers found on screenshots in this manual are provided for illustrative purposes only and may not be the same as the version number you see on screen.

WARNINGS AND CAUTIONS – PROGRAMMER

Note: See the RNS® System manual for warnings and cautions regarding other components of the RNS® System.

WARNING: PHYSICIAN AND CENTER ACCESS TO THE RNS® SYSTEM

The RNS® System should only be implanted by neurosurgeons with adequate experience in the implantation of subdural and stereotactic implantation of intraparenchymal electrodes and in the surgical treatment of intractable epilepsy. The RNS® System should only be used by neurologists or neurosurgeons with adequate experience in the management of intractable epilepsy and in the localization of epileptic foci, including the use of scalp and intracranial electrodes.

Neurologists and neurosurgeons using the RNS® System must have completed the NeuroPace® RNS® System training program. To qualify to manage patients with the RNS® System, physicians must demonstrate specific expertise related to epilepsy, video-EEG monitoring, interpretation of electrocorticograms (ECoGs), the pharmacology of antiepileptic medications and selection of patients for epilepsy surgery. Implantation of the RNS® System should be performed only by qualified neurosurgeons at centers capable of providing comprehensive epilepsy care, i.e. "Comprehensive Epilepsy Centers." These centers should have the expertise to provide diagnostic services that include video-EEG monitoring with scalp and intracranial electrodes and neuroimaging, and are experts in the treatment of epilepsy with antiepileptic medications, epilepsy surgery, and devices.

WARNING: MANAGEMENT OF PATIENTS WITH THE RNS® SYSTEM BY PHYSICIANS AT CENTERS THAT DO NOT PROVIDE THE SERVICES PROVIDED AT COMPREHENSIVE EPILEPSY CENTERS

In some instances, post-implant programming may be conducted by neurologists meeting the experience and certification requirements for neurologists at Comprehensive Epilepsy Centers, but who are not practicing in such centers. This situation might occur if the patient is not able to travel to a Comprehensive Epilepsy Center for regular follow-up (e.g. because of distance from the Center or limited access to transportation). These neurologists will be qualified by NeuroPace to provide post-implant programming. After NeuroPace® RNS® System training is complete, the qualified programming neurologist may receive external NeuroPace products (programmer, remote monitor).

WARNING: POTENTIAL SHOCK

Submerging any part of the programmer, or operating the programmer in or near a wet environment, may result in an electrical shock.

The programmer must be disconnected from the electrical outlet prior to cleaning to avoid the potential of electrical shock.

Electrical shock may occur if the programmer AC adapter and power cord are not properly connected to a grounded power source.

Do not attempt to modify or repair the wand or programmer. Contact NeuroPace Customer Support for assistance.

WARNING: RADIO FREQUENCY IDENTIFICATION (RFID) INTERFERENCE

Sources of RFID can result in signals that appear as ECoG activity to the neurostimulator. Signals that appear as ECoG activity could also result in delivering the programmed stimulation to the patient (per the device detection programming). The physician should be aware of possible sensing artifacts when assessing the ECoG recordings. Potential sources of RFID may occur in a health care environment, retail stores, public libraries, airports and business environments.

Refer to ***Electromagnetic Emissions and Immunity*** on page **105** for more information.

WARNING: SECURITY AND ELECTRONIC TRACKING SYSTEMS

Security screening devices (such as theft detectors, security tag deactivators, and airport security screening devices) can result in signals that appear as ECoG activity to the neurostimulator. Signals that appear as ECoG activity could also result in delivering the programmed stimulation to the patient (per the device detection programming). Such devices may be found at retail stores, public libraries and airports. The physician should be aware of possible sensing artifacts when assessing the ECoG recordings. Patients should be instructed to walk through the center of such security screening units without stopping, when possible, and exit the area of the screening device as soon as possible.

WARNING: NEUROPACE COMPONENTS

Use of accessories, transducers, and cables other than those provided by NeuroPace could result in increased electromagnetic emissions or decreased electromagnetic immunity of the RNS® System and result in improper operation. Do not plug the wand into equipment other than the programmer because it could damage the wand. Do not use a USB cable extension from the programmer to the wand.

WARNING: PORTABLE AND MOBILE RADIO FREQUENCY (RF) COMMUNICATIONS EQUIPMENT

Portable and mobile RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 12 inches (30 cm) to any part of the RNS® System, including cables. Otherwise, degradation of the performance of the RNS® System could result.

WARNING: NEUROPACE® EQUIPMENT PLACEMENT

Use of NeuroPace® equipment (for example, remote monitor or programmer) adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, the NeuroPace® equipment and other equipment should be observed to verify that they are operating normally.

Caution: Telemetry Artifact

Telemetry may produce an electrographic artifact. If responsive therapy is enabled with a sensitive detection set, detection of the electrographic artifact may occur resulting in therapy delivery. The physician should be aware of possible sensing artifacts when assessing the ECoG recordings.

Caution: Afterdischarge Activity

If evidence of afterdischarge activity resulting from stimulation is seen either on stored ECoGs or during test stimulation delivery, stimulation parameters should be adjusted to prevent such occurrence.

Caution: Programmer Failure

As with any electronic device, the programmer may be damaged or malfunction if the programmer AC adapter and power cord are not properly connected to a grounded power source.

Caution: Incompatibility of Programmer with Other Medical Devices

The effects of using the programmer to interrogate other electronic, programmable devices such as pacemakers, defibrillators, cochlear implants, and other neurostimulators or CPAP machines are unknown. It could result in reprogramming of the other device and therefore, the physicians familiar with each device should check the programmed parameters of each device before the patient is discharged and after each programming session of either device.

Caution: Electronic Interference

Communications between the programmer and the implanted neurostimulator may be interrupted by emissions from nearby electronic devices. Examples of sources of EMI are lithotripsy, computer monitors, cellular telephones, motorized wheel chairs, x-ray equipment and other monitoring equipment. Interruption of telemetry can result in incomplete communication. If EMI disrupts programming, move the programmer away from the likely source of EMI. Refer to **Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet** (page 102) for more information.

Caution: Placement of the Programmer Power Cord

Make sure nothing rests on the programmer power cord and that the cord is not located where it can be tripped over or stepped on.

Caution: Heating

The programmer's AC adapter may become hot during normal operation. Use care when handling during or immediately after operation.

STERILIZATION AND STORAGE

Product Storage

Components should be stored in a clean and secure area with a room temperature of approximately 14 to 28 degrees Celsius.

Explant and Disposal

Program all detection and therapy functions to DISABLED prior to explanting and shipping the RNS® Neurostimulator. Return the explanted neurostimulator and leads to NeuroPace. NeuroPace will provide shipping containers if requested.

DO NOT incinerate the neurostimulator; explosion can occur if the neurostimulator is exposed to incineration or cremation temperatures.

Return / Disposal of Programmer and Wand

The programmer and wand should be returned to NeuroPace if you are no longer using them. They contain electrical parts that need to be disposed of in accordance with local regulations.

Wand Maintenance, Cleaning and Sterilization

The wand is not user-serviceable and requires no maintenance. The wand can be cleaned by wiping with water. It can be placed in a sterile bag for use in the sterile field. DO NOT STERILIZE the wand.

Electrostatic Discharge (ESD) / Static Electricity

The ports on the tablet may be sensitive to electrostatic discharge / static electricity. Handle the tablet ports carefully. If exposed to electrostatic discharge, the tablet may experience telemetry artifacts or errors, or may freeze. In the event of a tablet freeze, refer to **Tablet Freezes** on page 103.

RNS® TABLET INSTRUCTIONS

This chapter primarily addresses use of the tablet during an interactive programming session, that is, when the clinician uses the wand for communication between the neurostimulator and tablet. For information on using the PDMS apart from an interactive programming session, either using the tablet or using a web browser on a personal computer, see **Patient List on Tablet and the PDMS on a Personal Computer** on page 90.

POWER ON THE TABLET

When interacting with a patient, the tablet should be used on battery power only. Make sure the tablet is sufficiently charged for use before powering on.

Press and hold the power button on the upper right side of the tablet. After 2-3 seconds, the tablet beeps to indicate it is powering on.

CONNECT AND LOGIN OVERVIEW

Access to all of the functionality of the tablet requires an Internet connection and two logins, in this order:

1. **Connect to the Internet** via Wi-Fi or mobile broadband on the tablet. Mobile broadband works automatically. To connect to Wi-Fi, see **Steps to Connect to a Wi-Fi Network** on page 12. Alternatively, the tablet computer can also support use of a docking station that has an Ethernet connector (wired connection). Contact NeuroPace customer service if you prefer this alternative.
2. **Login to the tablet.** NeuroPace provides the initial username and password for this login.
3. **Login to the Patient Data Management System (PDMS)**, which requires the Internet connection established in the first step. NeuroPace also provides the initial username and password for the PDMS. (For PDMS access using a personal computer, select the PDMS link at www.neuropace.com and log in. For further instructions, see **PDMS Access on a Personal Computer** on page 91.)

Note: Programming functions that are typically performed in the operating room do not require an Internet connection to the PDMS.

You can access all programming functions without connection to the PDMS except the following:

- View stored information such as ECoGs, neurostimulator activity and settings for any patient
 - Make a new detection set or adjust parameters of a detection set

For details, see **Functionality with and without a PDMS Connection** on page 11.

FUNCTIONALITY WITH AND WITHOUT A PDMS CONNECTION

Note: Programming functions that are typically performed in the operating room do not require an Internet connection to the PDMS.

The table below identifies the tablet functions that do not require an Internet connection to the PDMS and those that do.

Function	Internet Connection to the PDMS Required
Interrogate the neurostimulator	No
View neurostimulator information and settings downloaded during current session, including battery voltage, impedance, stimulation settings, recording montage	No
Create or change the recording montage	No
Set up or adjust ECoG capture settings	No
Select and program a default detection set or a custom detection set previously saved to the tablet from the PDMS	No
Select and program all stimulation settings	No
View stored information such as ECoGs, neurostimulator activity and settings for any patient	Yes
Make a new detection set or adjust parameters of a detection set	Yes

CONNECT TO THE INTERNET FOR PDMS ACCESS

When connected to the Internet and logged in to the Patient Data Management System (PDMS), the RNS® Tablet provides a look-through to the information stored on the PDMS, such as ECoGs, neurostimulator settings and activity. The PDMS is a patient database integrated with the RNS® Tablet; the PDMS is also available using a browser on any Internet-connected computer. The tablet ships with Wi-Fi enabled, so it can detect available Wi-Fi networks out of the box.

Note: If you attempt to connect to a Wi-Fi network that requires you to go to an Internet website to log in, as in some hotels and coffee shops, you will not be able to log in because the tablet does not have an Internet browser program. As a medical device, its Internet connectivity is restricted to the PDMS only.

The tablet ships with Wi-Fi and mobile broadband enabled.

Mobile Broadband Works Automatically

Mobile broadband is pre-configured by NeuroPace to be on and working, so no further effort is necessary and the tablet will connect to mobile broadband automatically if a signal is available.

Steps to Connect to a Wi-Fi Network

Follow these steps to connect via Wi-Fi to the Internet on the RNS® Tablet.

1. The tablet login screen appears first upon startup. Select a user. The screen prompts you for a username and password and displays the keyboard, but ignore these for now.
2. Touch to select the small symbol that indicates network status. (If necessary, touch outside the login fields to dismiss the keyboard and see the whole screen.)


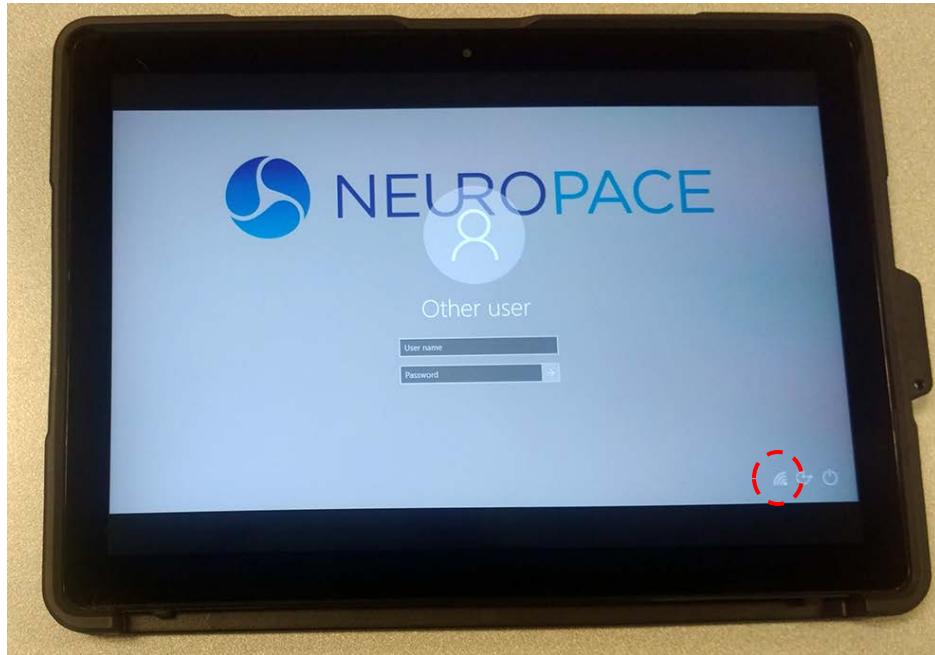


Figure 1 Touching the network icon at lower left



When a wireless network is available, it looks like the signal strength symbol with an asterisk above it. A different symbol is present if there is no wireless network available. Select it in any case to open the **NETWORKS** panel on the right side of the screen.

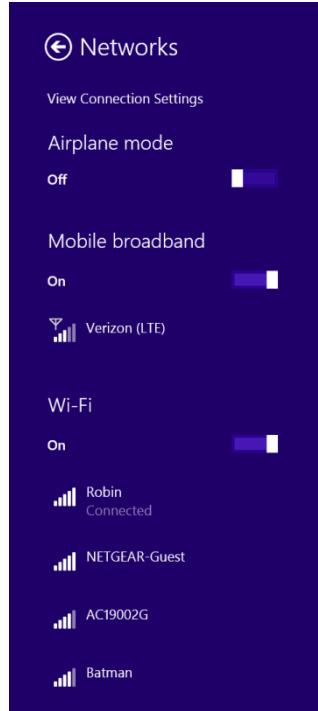

3. From this point, Wi-Fi access is the same as on a smartphone or laptop computer:
 - a. Select the desired network and then select **CONNECT**.
 - b. Enter the network security key and select **NEXT**. (The security key is your Wi-Fi password.) The tablet connects if the security key is correct. Internet connection is complete.

If You Lose Network Connectivity

If you lose network connectivity, you will not lose information you already programmed into the neurostimulator or saved to the PDMS, but you will lose information not already programmed into the neurostimulator or saved to the PDMS.

If you are already logged in to the tablet and have lost network connectivity, first log off to return to the login screen: return to the Home screen and touch the circle-X (exit) button at upper right, and then select **EXIT** on the **PROGRAMMER EXIT** dialog. Then follow the **Steps to Connect to a Wi-Fi Network** on page 12. If you still cannot connect, contact NeuroPace for further assistance (see **Contacting NeuroPace** on page 5).

LOG IN TO THE RNS® TABLET

This is the first of two logins. The second one gives you access to the PDMS.

1. The tablet login screen appears first upon startup. Select a user. The screen prompts you for a username and password and displays the keyboard.



2. Enter your username and password (initially provided by NeuroPace) and select **Enter**. Login completes if the password is correct. The tablet automatically checks for software updates.
 - a. If updates are available, the **NEUROPACE UPDATE CENTER** opens. Follow the instructions below to use the update center.
 - b. If there are no updates, the **NEUROPACE LOG IN** screen appears, which lets you **LOG IN TO THE PDMS** (see page 17).

NEUROPACE UPDATE CENTER

After logging into the tablet, the **NEUROPACE UPDATE CENTER** opens if any updates are available for installation. As described below, the Update Center identifies the update type(s) and allows you to select a convenient time to perform the updates.

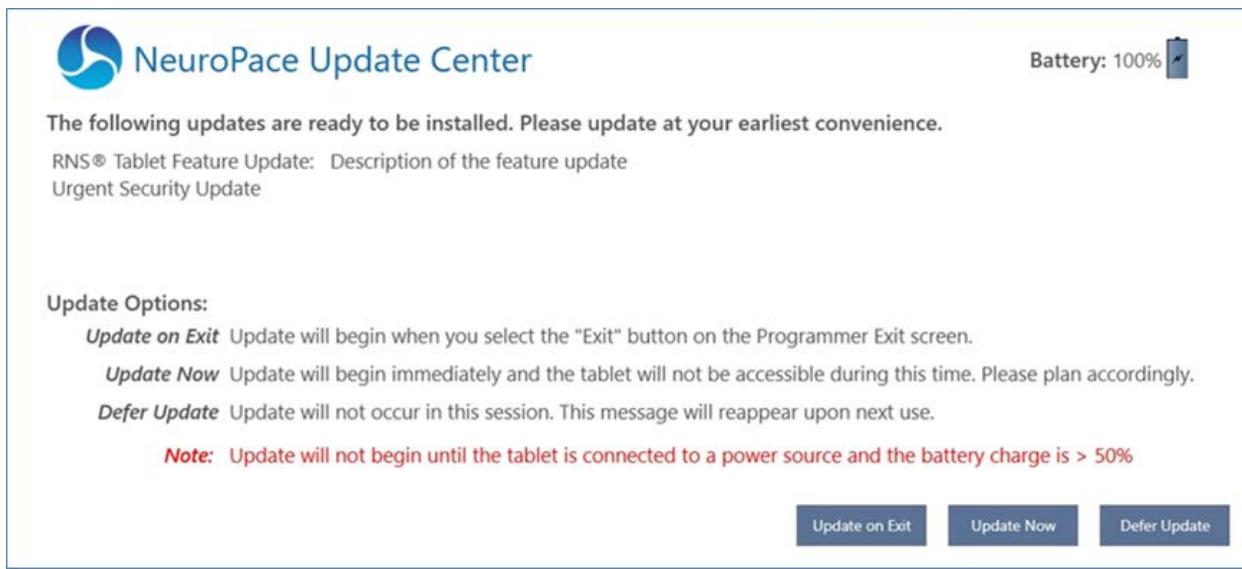


Figure 2 NeuroPace Update Center

Update Types

1. **RNS Tablet Feature Update:** These updates contain tablet interface changes or new feature updates.
2. **Urgent Security Update:** These updates are required to ensure cybersecurity for the tablet and should be done as soon as possible. They do not include any tablet interface changes or new features.

Note: Some updates require the tablet to restart. Under these circumstances, upon completion of the update, the system will restart automatically and will require you to log in again.

Update Options

It is important to install updates, however the time it takes for an update to complete can vary from less than a minute to up to an hour. Accordingly, be sure to perform an update when you have enough time. Once an update starts it cannot be cancelled; interrupting an update by powering off the tablet can render the tablet unusable.

1. **UPDATE ON EXIT:** This option postpones the update until you select **Exit** on the **Programmer Exit** dialog. To maintain tablet availability within the next hour (i.e. for a patient clinic visit), we recommend selecting **UPDATE ON EXIT**.
2. **UPDATE NOW:** This option starts the update immediately and the tablet will not be available during this time. Do not use this option if the tablet is needed within the next hour.

3. **DEFER UPDATE:** This option postpones the update until the next use of the tablet, at which time the **NEUROPACE UPDATE CENTER** will appear again.

Update Process and Power Requirements

When an update begins, the screen name shows **NEUROPACE UPDATE CENTER: INSTALLING UPDATE(S)** and reports update progress.

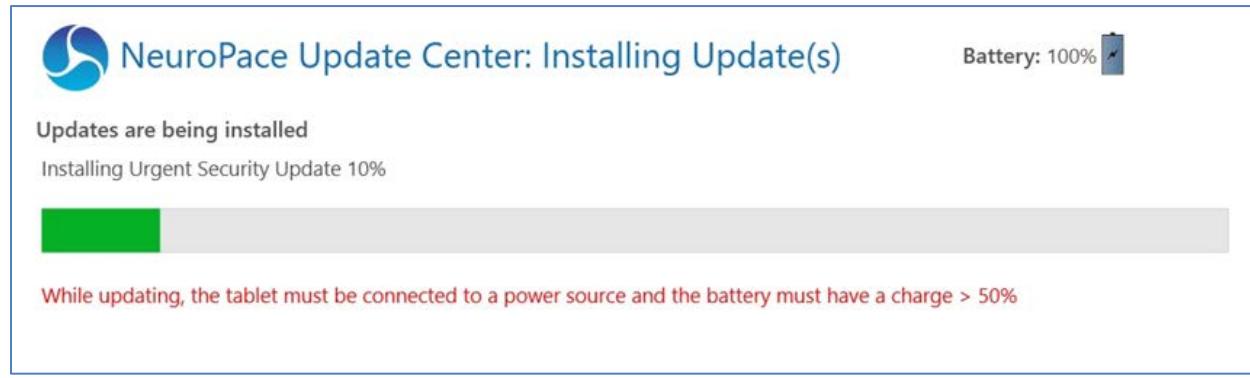


Figure 3 Installing Updates Progress

- **Power requirements:** **Updates cannot proceed until tablet is connected to a power source and the battery charge is greater than 50%.** These power requirements ensure that the tablet remains powered on throughout the update, so it is not interrupted. It is recommended to keep the tablet plugged in until the update completes.

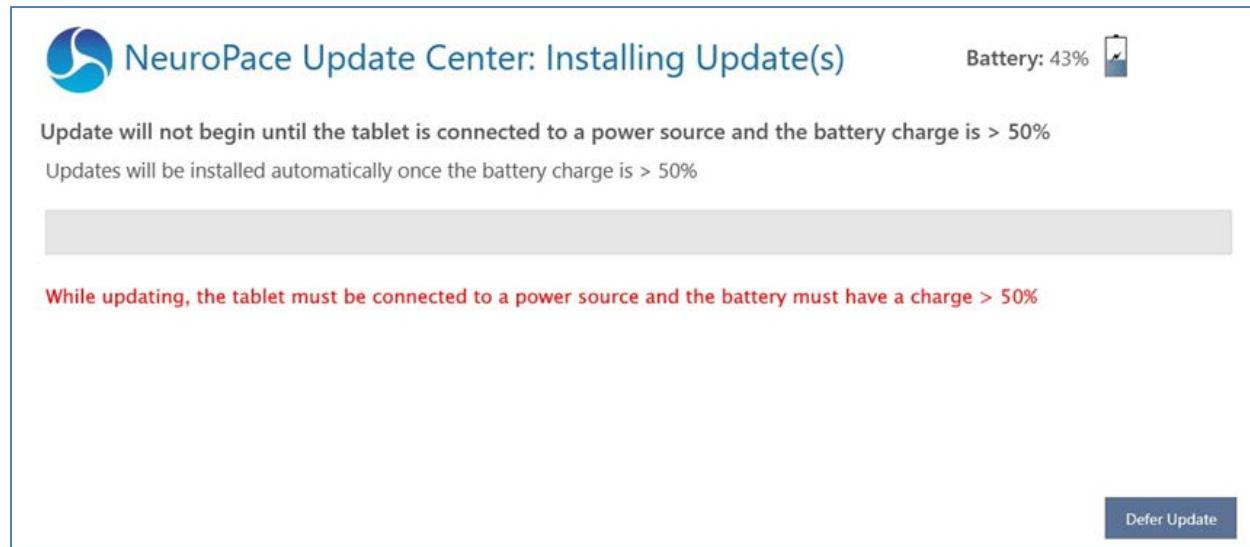


Figure 4 Power requirements must be met before updates will start

If the power requirements are not met, instructions to address the unmet power requirements will be displayed and the progress status bar will not commence (see **Figure 4**). Follow the instructions on

screen as applicable to either plug in the tablet or wait for the battery to reach the required level. When the power requirements are met, the update will begin automatically.

If the power requirements have not been met and the update has not begun, you can select the option to **DEFER UPDATE** until the next time you log in to the tablet, at which time the **NEUROPACE UPDATE CENTER** will appear again.

LOG IN TO THE PDMS

The **NEUROPACE Log In** screen appears after tablet login, and whenever you attempt to access functions that require PDMS access, like the **PATIENT LIST** and **ECoG LIBRARY**.

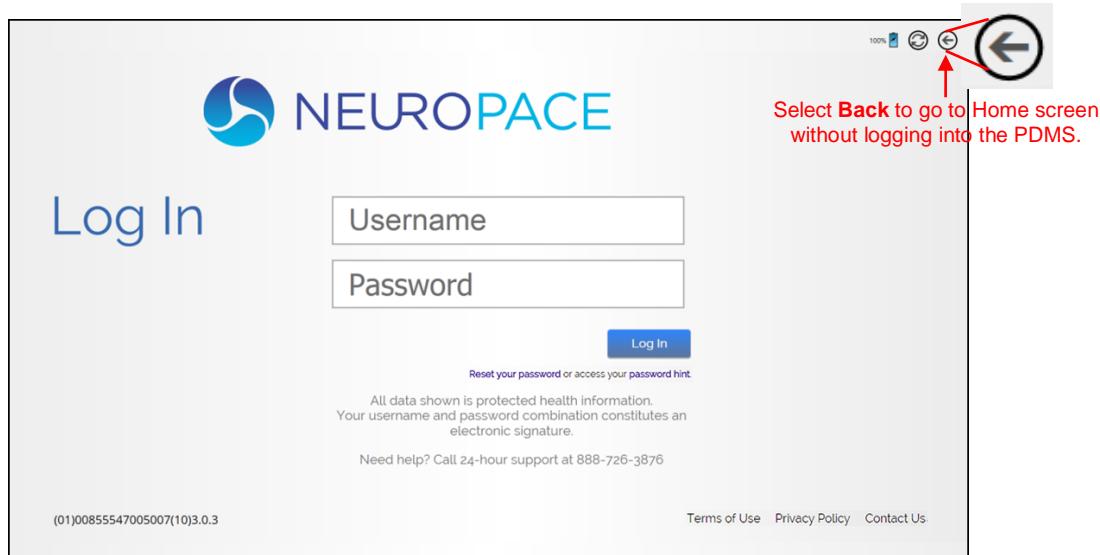


Figure 5 NeuroPace Log In screen—for PDMS login

To skip PDMS login, select the **BACK** button in the upper right corner. The Start screen opens.

1. Select the fields to enter the username and password and then select **Log IN**. Login completes and the Start screen opens (**Figure 3**).



Figure 6 Start screen

INTERROGATE THE RNS® NEUROSTIMULATOR

To interrogate the RNS® Neurostimulator means to identify it, perform routine tests and measurements, and retrieve currently programmed settings and stored data from it. Follow these steps to interrogate:

1. Connect the wand to the USB port on the tablet.

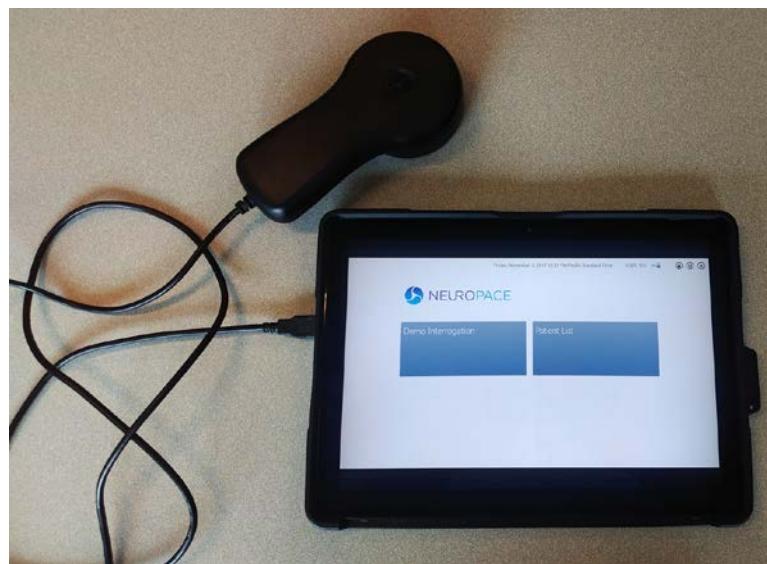


Figure 7 Tablet with connected wand

2. Place the wand within approximately 1 inch of the neurostimulator, concave side of wand facing the neurostimulator. Then select **INTERROGATE** from the Start screen (left tile—see **Figure 3**). You can also select the green **INTERROGATE** tile at upper right in the **HOME** screen (see **Figure 6**.) The wand must be held in place over the neurostimulator for interrogation to succeed.

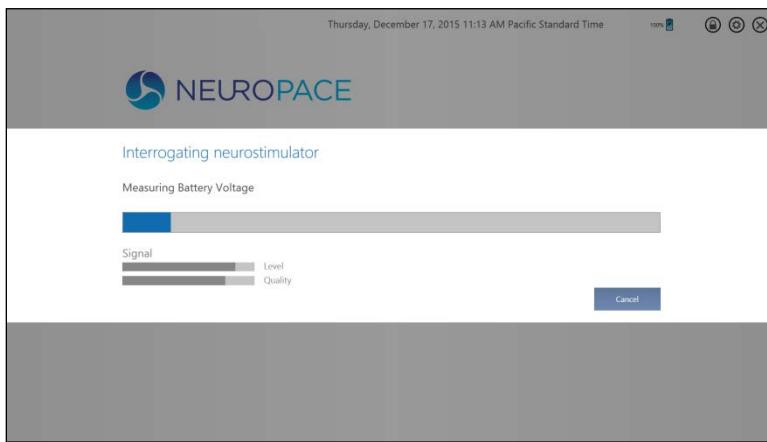


Figure 8 Interrogating neurostimulator dialog—interrogation in progress

The **INTERROGATING NEUROSTIMULATOR** dialog shows progress and the **HOME** screen appears when complete. It takes from 10 seconds up to 5 minutes, depending on the amount of data to be retrieved. Interrogation does not significantly drain the neurostimulator battery.

The dialog has bars that show the **LEVEL** and **QUALITY** of the **SIGNAL**. The **LEVEL** measures the strength of the communication between the neurostimulator and the tablet. The **QUALITY** measures the absence of electronic interference in the environment. Poor signal level or quality can result in **Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet**—see page **102** for troubleshooting information.

DATA RETENTION AND ACCESS ON THE NEUROSTIMULATOR, PROGRAMMER AND PDMS

Programmed Settings

The neurostimulator retains its programmed settings until it is programmed again. These include the recording montage, ECoG capture triggers and reservations, and the detection and stimulation settings.

*Note: ECoG capture triggers are inactive during programming. **Programming clears activity information (diagnostics) and ECoGs stored in the neurostimulator.** If there is any new information, be sure to interrogate the neurostimulator prior to programming.*

Storage of ECoGs and Neurostimulator Activity (Diagnostics)

Between interrogations by the patient or a clinician, the neurostimulator stores ECoGs and neurostimulator activity events up to its storage limit. (Neurostimulator activity events are also known as diagnostics, which include, for example, detection and stimulation events and magnet swipes.) Once the memory for ECoGs or diagnostics is full, the neurostimulator will store new data by

overwriting the oldest ECoG or diagnostics. Reserved ECoGs are overwritten only by newer ECoGs of the same ECoG capture trigger type.

- The limit to the number of ECoGs stored on the neurostimulator depends on the capture duration you select in the ECoG capture settings (see **Set up ECoG Capture** on page 44 for details).
- You can make no adjustments that affect the number of neurostimulator activity events that can be stored. However, the number of events that can be stored is large compared to the number of ECoGs, because the data size of each event is relatively small.

At each interrogation with the remote monitor or programmer, the neurostimulator transfers all stored ECoGs and diagnostics to the programmer or remote monitor, at which point the successfully transferred data is cleared from the neurostimulator.

- If the programmer or remote monitor is not connected to the PDMS (via the Internet), the transferred data remains on the programmer or remote monitor hard drive. You cannot view ECoGs or diagnostics while stored on the programmer or remote monitor; it must be transferred to the PDMS before you can view it.
- When the programmer is connected to the PDMS (via the Internet), the transferred data is automatically uploaded to the PDMS and cleared from the programmer memory. When the patient synchronizes using the remote monitor, data is transferred to the PDMS and cleared from the remote monitor. Data is retained on the PDMS.

Note: If automatic upload to the PDMS does not complete successfully, data is not lost. Data is cleared from the programmer or remote monitor only after a successful upload.

HOME SCREEN

The Tablet Home screen is the launch point for all programming activities.


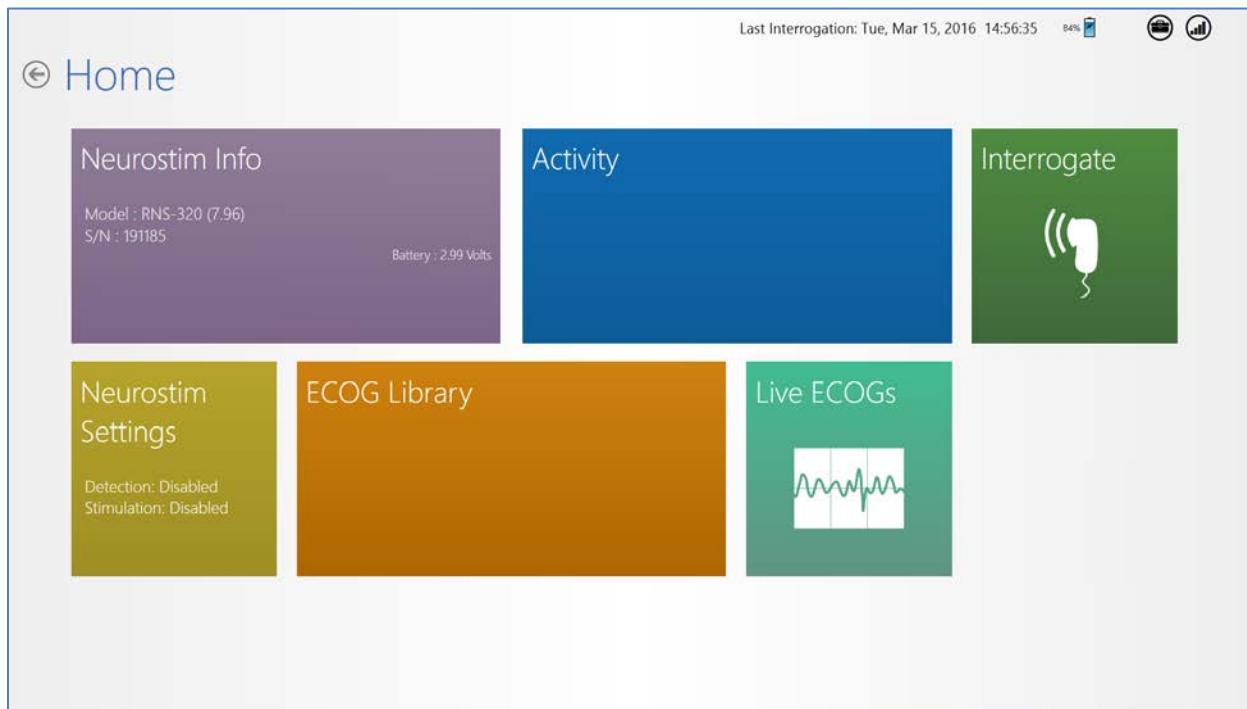


Figure 9 Tablet Home screen

Note: When logged in to the PDMS, changes you make with the tablet are immediately saved to the PDMS. If you lose connectivity, you will not lose information you already programmed into the neurostimulator or saved to the PDMS, but you will lose information not already programmed into the neurostimulator or saved to the PDMS.

The elements in the far upper right of the screen are present on all tablet screens.

- **LAST INTERROGATION** date and time (24-hour time)
- **LOGGED IN** username
- Tablet battery status (%)

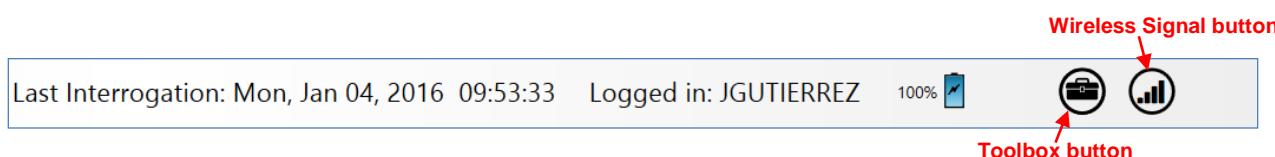

- **WIRELESS SIGNAL** button to access the **PDMS COMMUNICATION LOG**

PDMS Communication Log					
Time	Task	Status	Message	% Complete	Error Code
1/8/2016 12:35:05 PM	Background	DBMGR_SYNC	Synchronization completed successfully.	100	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Transferring Data...	66	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Transferring Data...	33	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Transferring Data.	0	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Transferring Data	0	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Transferring Data	0	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Connecting to PDMS.	100	
1/8/2016 12:35:05 PM	Background	DBMGR_UPLOADING	Connecting to PDMS.	0	

Sync Cancel

- **TOOLBOX** button to access **UTILITIES** dialog to **IMPORT DEVICE** and **EXPORT DEVICE SETTINGS**. These features are designed to support neurostimulator replacement. For instructions, see **Export and Import Device Settings for Neurostimulator Replacement** on page 87.



Figure 10 Utilities dialog

NEUROSTIM INFO SCREEN

Select **NEUROSTIM INFO** from the **HOME** screen to open the **NEUROSTIM INFO** screen.


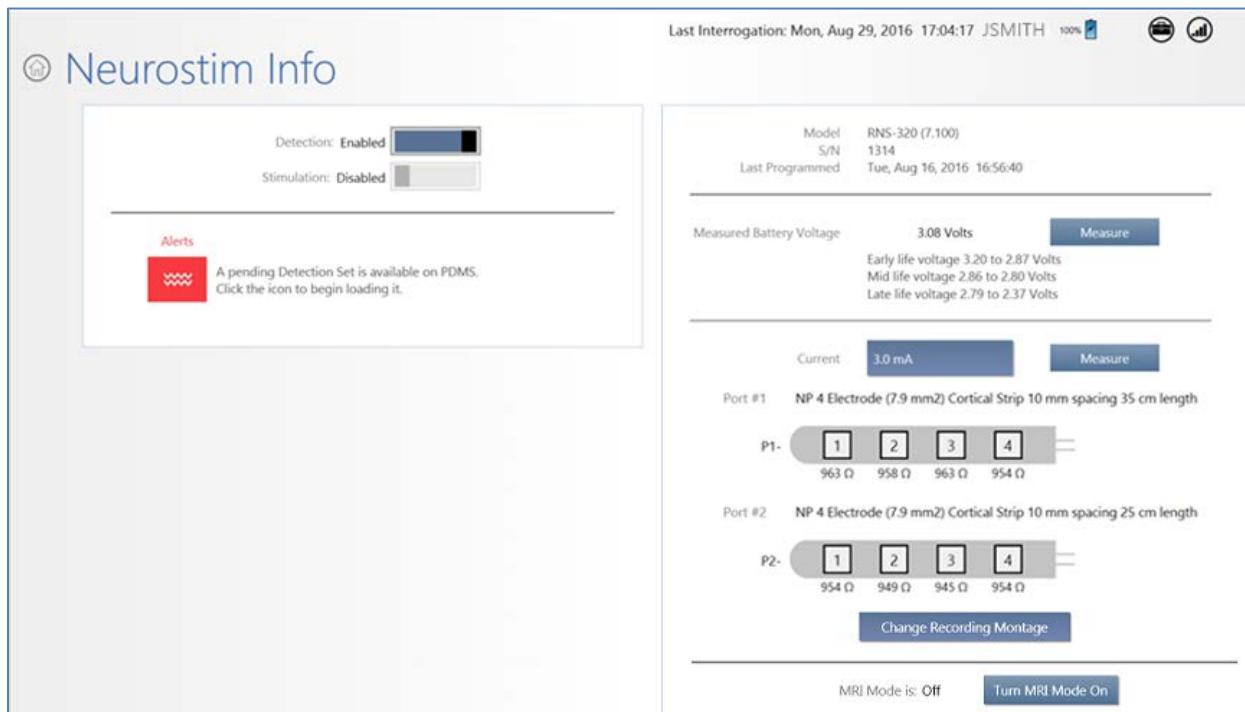



Figure 11 Neurostim Info screen

The **NEUROSTIM INFO** screen provides information about the neurostimulator, including:


- At upper left, whether **DETECTION** and **STIMULATION** are **ENABLED** or **DISABLED** and Alerts, if any

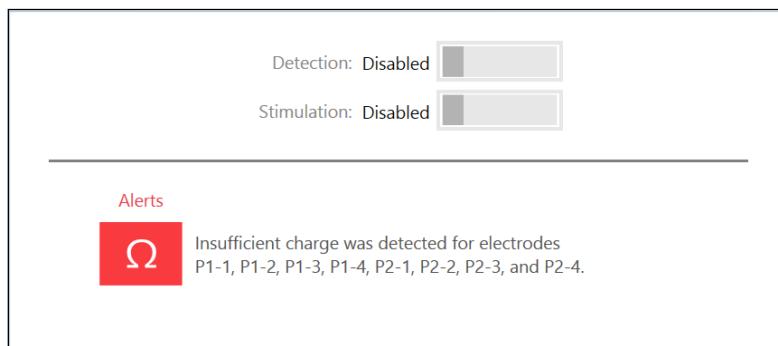


- At upper right, **MODEL**, serial number (**S/N**) and **LAST PROGRAMMED** date and time (24-hour time)

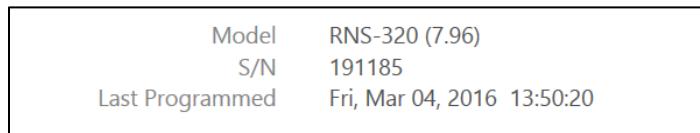

- At middle right, **MEASURED BATTERY VOLTAGE** with **MEASURE** button. The early, mid and late life voltage ranges are displayed to give context to the measured voltage. See **Significance of Battery Voltage** on page 27.

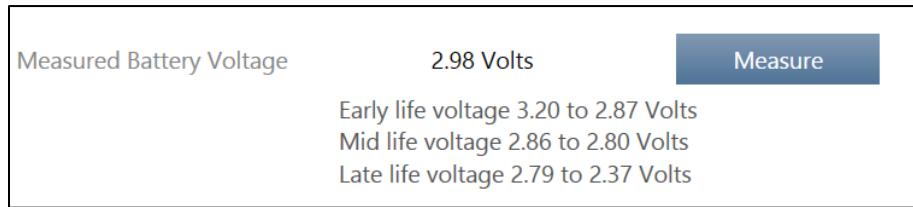

- At lower right, measured impedance for electrodes on the **PORT #1** and **PORT #2** leads

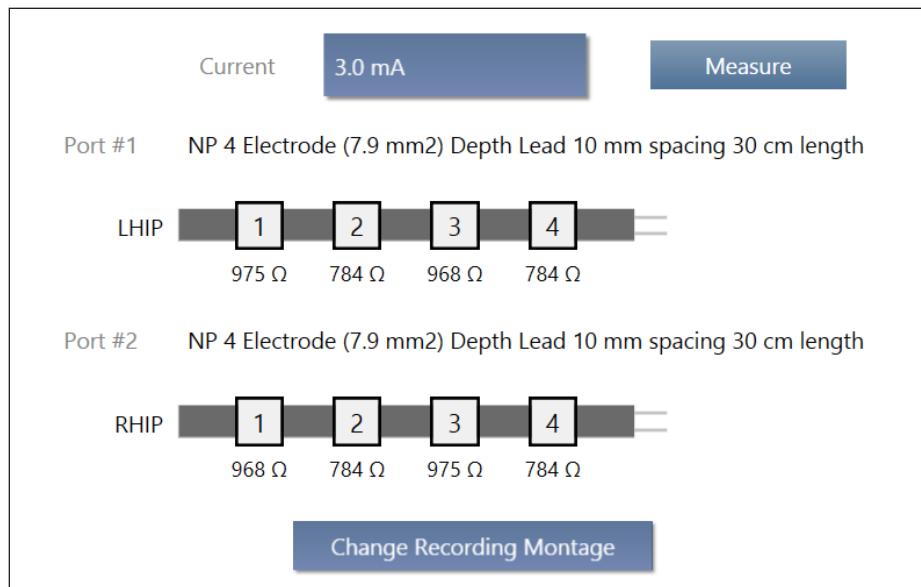

- Buttons to **MEASURE** impedance manually using the amount of **CURRENT** selected. (It is measured routinely during interrogation.)
- A button to **CHANGE RECORDING MONTAGE**, as explained under **Creating the Montage** on page 39

MRI Mode for RNS Neurostimulator Model RNS-320

Note: MRI Mode software buttons and features are available only when you interrogate an RNS® Neurostimulator model RNS-320, which is MR Conditional. Refer to the *MRI Guidelines for the RNS® System* to determine patient eligibility for MRI and for the specific conditions required to safely perform an MRI scan on patients implanted with an RNS® Neurostimulator model RNS-320. Scanning under different conditions may result in device damage or malfunction and serious patient risks including permanent brain damage which may cause severe injury, coma, or death. All NeuroPace® manuals are available at www.NeuroPace.com or by contacting NeuroPace, Inc.



RNS® Neurostimulator model RNS-300M of the RNS® System is MR Unsafe. Having an MRI scan with a model RNS-300M neurostimulator implanted may result in serious injury or possible death.


RNS® System External Components: All external components and accessories such as the Magnet, RNS® Tablet, NeuroPace® Programmer, NeuroPace® Remote Monitor, and Wand are MR Unsafe and can pose a projectile hazard in the MR environment, and therefore, must be kept out of the MRI scanner room.

Instructions to Turn MRI Mode ON


Using an RNS® Tablet, **MRI Mode** must be turned **On** prior to a patient implanted with an RNS® Neurostimulator model RNS-320 receiving an MR scan. This step must be completed outside the MRI scanner room, because the RNS® Tablet and Wand are MR Unsafe. Log in to the RNS® Tablet and interrogate the RNS® Neurostimulator. Select **Neurostim Info** from the **Home** screen to open the **Neurostim Info** screen. At bottom right, the **Neurostim Info** screen reports **MRI Mode is: Off** and provides the corresponding button to **Turn MRI Mode On**.



To turn MRI Mode on, place the wand within approximately 1 inch of the neurostimulator, concave side of wand facing the neurostimulator, and select **Turn MRI Mode On**. The tablet immediately sends the programming signal to turn MRI Mode on. A telemetry dialog indicates progress. Upon completion, the **Neurostim Info** screen displays an alert at upper left indicating that MRI Mode is on.

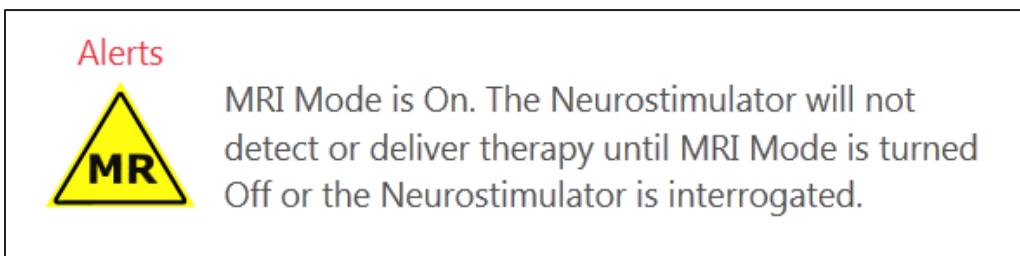

At bottom right, the screen reports **MRI Mode is: On**, and the button toggles to **Turn MRI Mode Off**.

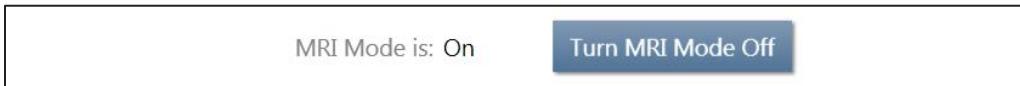

Note: It is strongly recommended that the neurostimulator remain in MRI Mode only for as long as necessary. While in MRI Mode, the neurostimulator is not detecting or delivering therapy to the patient. Additionally, the neurostimulator uses more battery power in MRI Mode than in normal operating mode. Note that the neurostimulator can be in MRI Mode for up to approximately two days per year without affecting battery longevity. Battery longevity estimates are specified in the RNS® System Physician Manual.

Note: Device status, including MRI Mode status, will propagate to the Patient Data Management System (PDMS) whenever the RNS® Tablet is connected to the internet, and therefore become available to all authorized users.

Instructions to Turn MRI Mode OFF

Using an RNS® Tablet, **MRI MODE** must be turned **OFF** after a patient implanted with an RNS® Neurostimulator model RNS-320 receives an MR scan. This step must be completed outside the MRI scanner room, because the RNS® Tablet and Wand are MR Unsafe.

Note: Please verify that the MRI scan was performed or cancelled prior to turning off MRI Mode.

Method 1: If the **NEUROSTIM INFO** screen remains open on the RNS® Tablet that was used to turn MRI Mode on, the button at bottom right of the **NEUROSTIM INFO** screen can be used to turn MRI Mode off. Select **TURN MRI MODE OFF** (while holding the concave side of the wand within approximately 1 inch of the neurostimulator) to turn MRI Mode off. This scenario may occur, for example, if a clinician with an RNS® Tablet is present immediately before the MR scan to turn MRI Mode on, and afterward to turn MRI Mode off.

Method 2: Interrogation of the RNS® Neurostimulator with any RNS® Tablet turns off MRI Mode. During interrogation of a neurostimulator with MRI Mode on, a notification on the RNS® Tablet says that continuing the interrogation will turn MRI Mode off. Select **CONTINUE** to proceed with interrogation and turn MRI Mode off.

Electrode Impedance and Battery Voltage

Interrogation includes measuring electrode impedance and battery voltage, and these measurements are reported on the **NEUROSTIM INFO** screen. The **NEUROSTIM INFO** tile and screen show an alert if the impedance measurement detects insufficient charge for one or more electrodes.

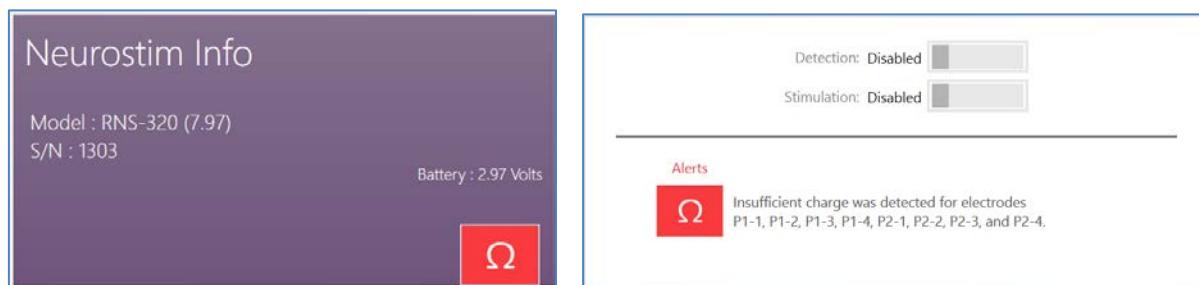


Figure 12 Insufficient charge alert shown on the **NEUROSTIM INFO** tile and screen

To measure these items again, place the wand within approximately 1 inch of the neurostimulator and select the applicable **MEASURE** button on that screen (select **NEUROSTIM INFO** from the **HOME** screen to

get there). For impedance measurements, use the **CURRENT** field to adjust the current delivered during measurement.

Note: The system cannot measure impedance during a detected episode or during therapy delivery. To measure impedance, either wait for the episode or therapy to complete, or disable detection and therapy.

Significance of Impedance

Impedance measurements are useful for verifying the integrity of a particular electrode. Each impedance measurement delivers a single pulse of current between an electrode and the neurostimulator canister (can) to measure the voltage that results as the current flows. Impedance measurements are reported in Ohms (voltage divided by current).

Impedance measurements greater than 3500 Ohms or less than 250 Ohms are considered abnormal. If impedance measurements are abnormally high, there may be insufficient charge to deliver therapy. When this is the case, alerts may appear on the Neurostim Info screen as shown below in **Figure 10**.

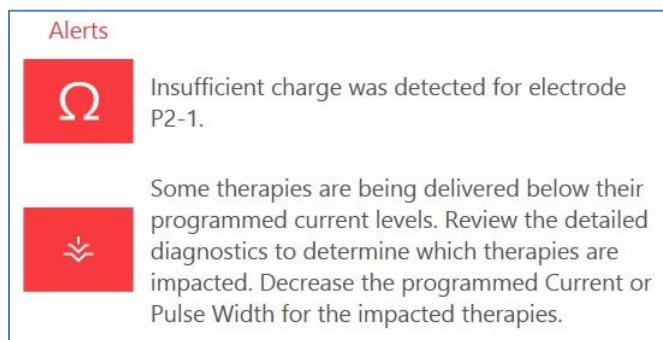


Figure 13 Alerts when there is insufficient charge to deliver therapies

For troubleshooting information, refer to **Abnormal Lead Impedance (greater than 3500 Ohms or less than 250 Ohms)** on page **98**; and **Insufficient Charge** on page **101**.

Significance of Battery Voltage

Normal RNS[®] Neurostimulator battery voltage measurements are between 2.38 and 3.20 Volts. When the battery voltage measurement drops to 2.37 Volts or less, neurostimulator replacement is recommended, and the tablet notifies you on screen by showing the elective replacement indicator (ERI).

ACTIVITY SCREEN

Select **ACTIVITY** from the **HOME** screen to open the **ACTIVITY** screen.

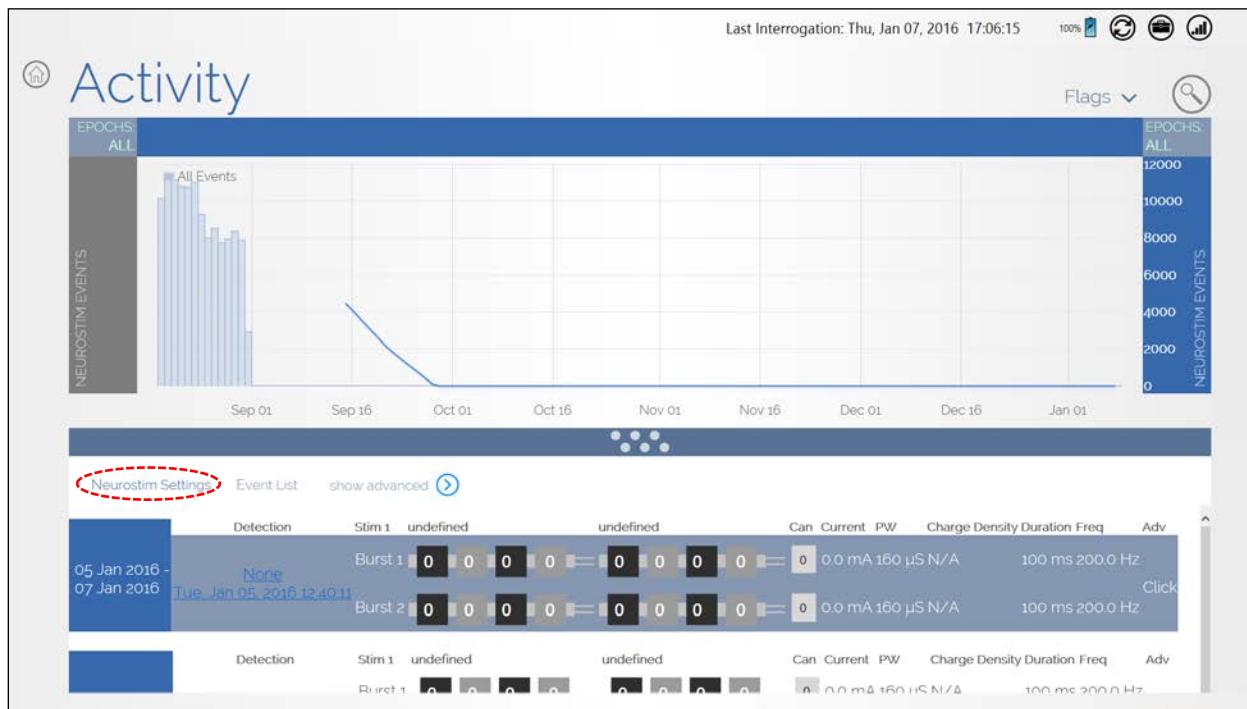


Figure 14 Activity screen


Note: You must be logged in to the PDMS to access the **ACTIVITY** screen. If you are not, the **NEUROPACE LOG IN** screen opens when you select **ACTIVITY**, prompting you to log in. (See **Log In to the PDMS** on page 16.) If you are not connected to the Internet, see **Connect to the Internet for PDMS Access** on page 11.

The upper part of the **ACTIVITY** screen is a histogram that shows the number of detection and stimulation events over time. Swipe left or right over the dates to view earlier or later activity. The lower part of the screen shows (by default) a list of the **NEUROSTIM SETTINGS**—detection and stimulation settings—used over time, one row for each period of time between programming changes, which is called a programming epoch. More view options are described below.


Activity Histogram Viewing Options

The interface provides many options for viewing activity.

- Select **EPOCHS** at upper left or right to see events from **ALL EPOCHS**, **DETECTION ONLY** or **STIMULATION ONLY**.
- Select **FLAGS** at upper right to see events flagged with **MAGNET SWIPE**, **OFFICE INTERROGATION** or **IMPLANTS**.

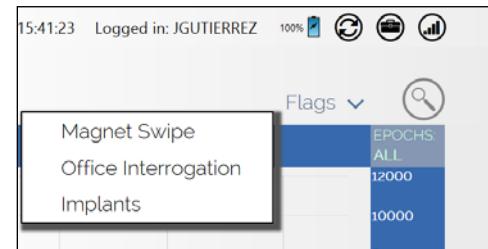

- Select **NEUROSTIM EVENTS** on the left or right sides of the histogram to select items for layering on the graph; you can select multiple items concurrently, up to all items in the upper part of the list. A color legend appears on the left to help you recognize what is on the histogram.

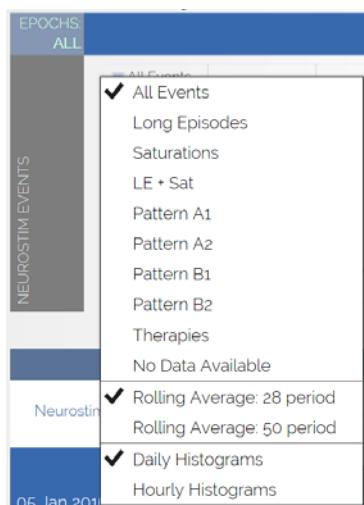

Select event types to layer on graph

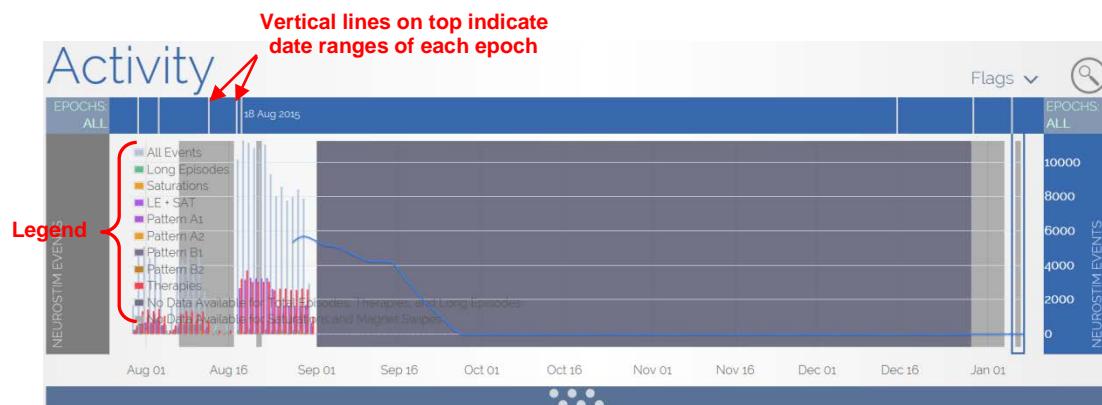


Figure 15 Activity histogram with all events showing and supporting legend on left

View Neurostim Settings

The lower part of the **ACTIVITY** screen details the history of detection and stimulation settings used—collectively called **NEUROSTIM SETTINGS**—between each iteration of programming (each programming epoch), which are in rows from most recent on top to oldest on the bottom. Drag up the gray separator bar to view more of these settings, which will cover the upper histogram portion. Swipe up on the setting rows to view more rows, if they are present.

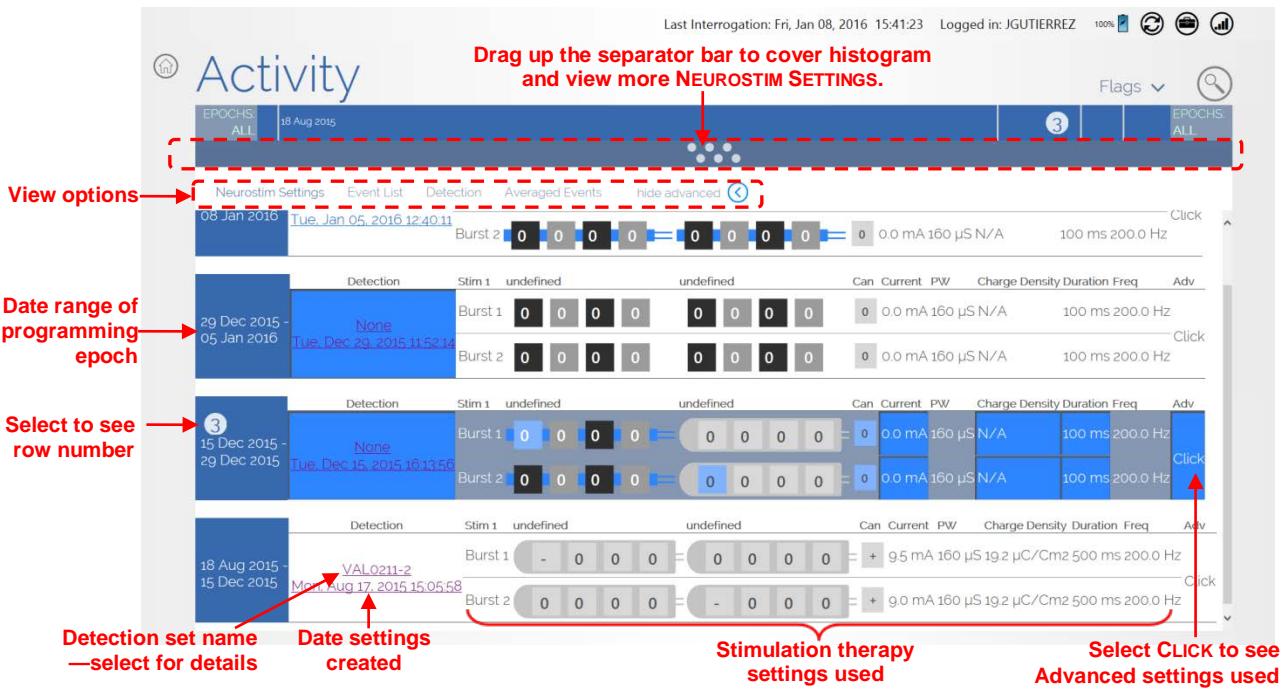


Figure 16 Activity screen—view detection and stimulation settings used between programming sessions

At the left end of each row is a blue tile showing the date range when those **NEUROSTIM SETTINGS** were used. Vertical lines above the histogram correspond to these dates (see **Figure 12** above). When you select a blue tile, its row number appears on the tile (see **Figure 13** above, row numbers increase as you go down).

Also the corresponding time period is bracketed in the histogram graph, and the row number and start date appear over the bracketed time period, if there is enough room on screen (see **FIGURE 14**). This helps to relate the programming changes to the device activity.

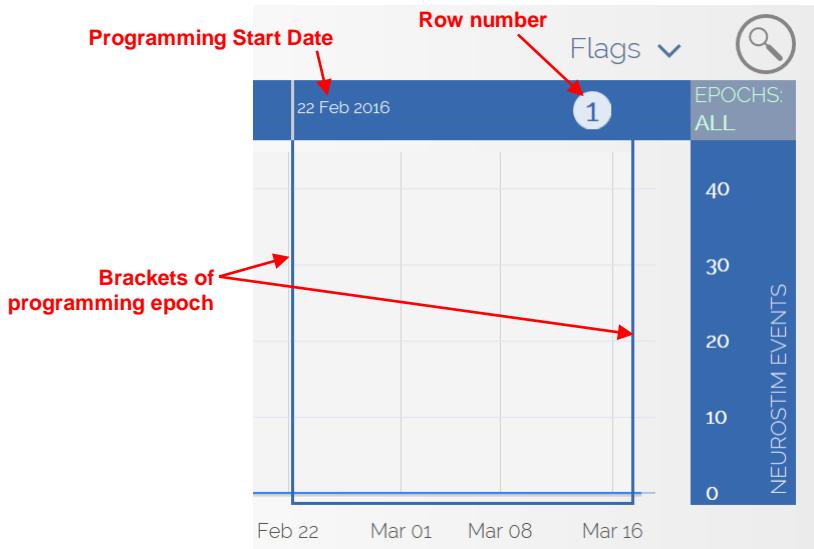


Figure 17 Activity histogram excerpt showing programming epoch 1 bracketed

Next to the blue tiles, the second column in **Figure 13** shows the name of the **DETECTION** set. Select the name to view its details in the screen that opens (**Figure 15**). Touch (X) at upper right to close.

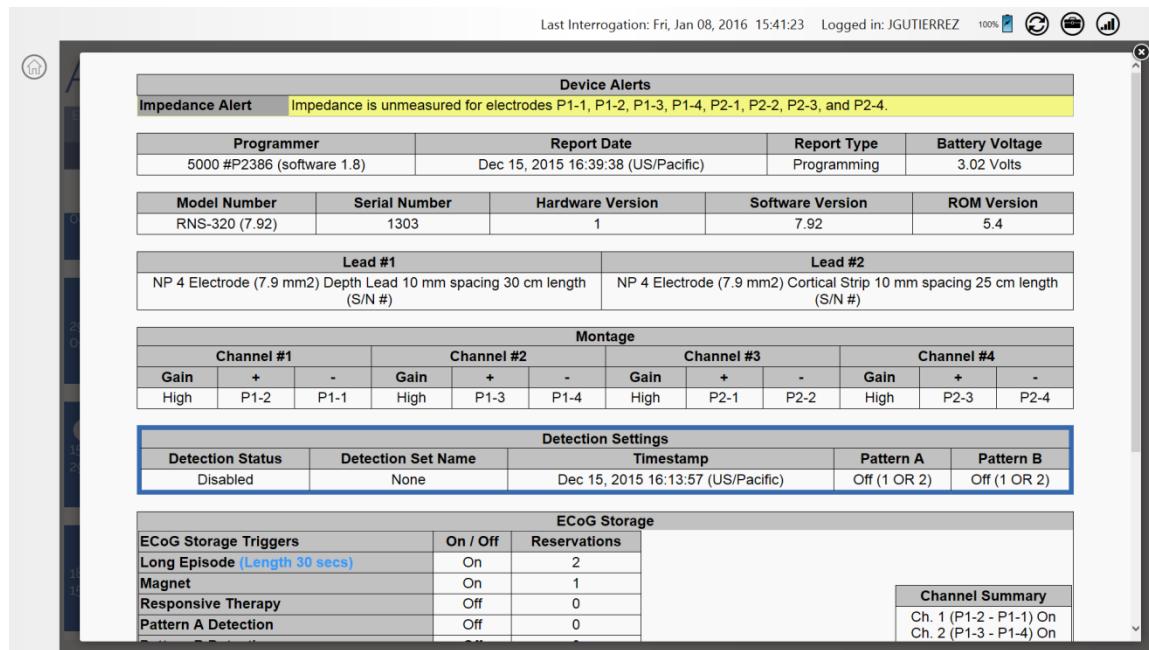


Figure 18 Detection set details

At the right end of each row in **Figure 13** is the ADV (Advanced) column. Touch CLICK below ADV to view any advanced settings used in the screen that opens (**Figure 16**). Touch (X) at upper right to close.

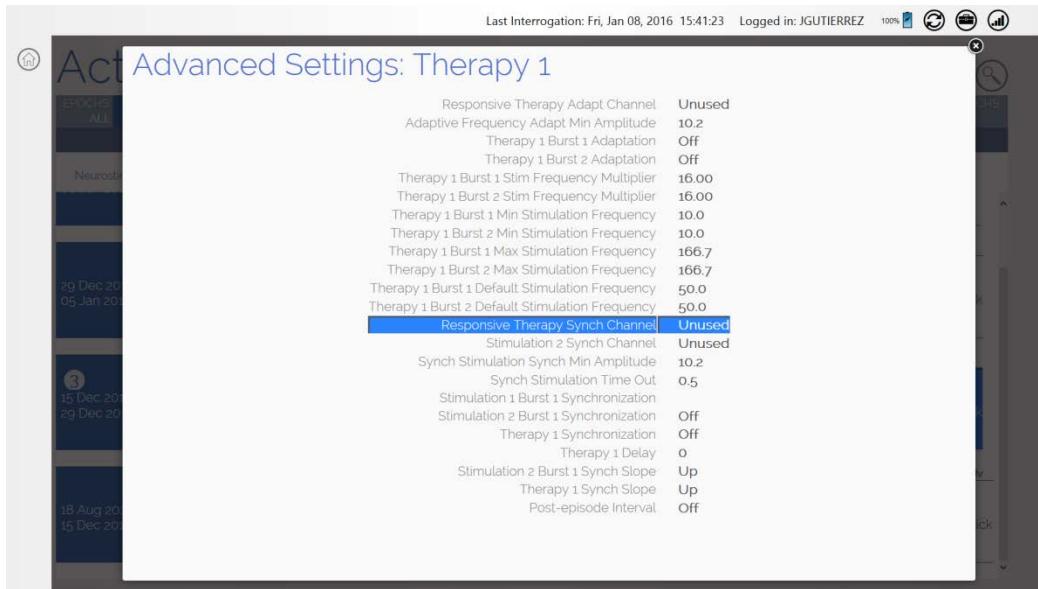


Figure 19 Advanced stimulation therapy settings used

View Event List

Select **EVENT LIST** next to **NEUROSTIM SETTINGS** above the left end of the row (**Figure 17**) to see a complete list of events for all rows. From there, select links under **REPORT TYPES** for detailed reports on each event. Swipe up or down in the blue column on the left to scroll between programming epochs. Swipe up or down within the event list itself to scroll the list of events.

Date	Time	Duration	PGM/DTR	Report Type
12/12/2015	17:46:00		5000	Battery Measurement
12/12/2015	17:45:55		5000	Battery Measurement
12/12/2015	17:45:44		5000	Electrode Impedance
12/12/2015	17:44:49		5000	Initial Interrogation
12/12/2015	17:44:49		5000	Electrode Impedance
12/12/2015	17:42:01		5000	Electrode Impedance
12/12/2015	17:41:11		5000	Initial Interrogation
12/12/2015	17:41:11		5000	Electrode Impedance
12/05/2015	17:21:24		5000	Electrode Impedance
12/05/2015	17:20:36		5000	Electrode Impedance
12/05/2015	17:19:28		5000	Electrode Impedance
12/05/2015	17:18:51		5000	Electrode Impedance
12/05/2015	17:18:35		5000	Programming
03/01/2015	00:00:00			Neurostimulator Implant IRNS-300ML - Serial Number: 1404

Figure 20 Event List selected—lists all events in a Neurostim Settings period, with links to reports

Types of Reports

There are six types of patient reports:

- **Interrogation** (Initial or Subsequent within the tablet session) ▶ Interrogation time/date
 - Programmed Neurostimulator settings
 - Event counters (# of detections, stimulations, etc.)
 - Event date/time and details
- **Programming**
 - Programming date/time
 - Programmed Neurostimulator settings
- **Electrode Impedance**
 - Impedance measurement date/time
 - Impedance measurement of each electrode (in ohms)
- **Battery Measurement**
 - Battery measurement date/time
 - Battery measurement value (in volts)

- **Therapy Testing** (shows therapy testing values)
- **ECoG**
 - ECoG trigger time/date
 - ECoG record
 - ECoG trigger reason
 - ECoG markers (detections, stimulations, etc.)

Advanced Views

DETECTION, **AVERAGED EVENTS** and **PROGRAMMING EPOCHS** are the advanced views. **DETECTION** shows a graph of detection occurrence percentage for each programmed detection pattern since the last change in detection settings. **AVERAGED EVENTS** shows a week by hour histogram of the normalized average number of detected episodes per hour for the last 12 weeks. **PROGRAMMING EPOCHS** shows the history of recent programming epochs, summarizes the settings and activity of each epoch, and highlights the changes in settings and activity from one epoch to the next.

ECoG LIBRARY SCREEN

It is recommended to review ECoGs stored in the library to assess how early and how often detections occur and the electrographic response to stimulation therapies delivered. Select **ECoG LIBRARY** from the **HOME** screen to open the **ECoG LIBRARY** screen, which displays thumbnails of all ECoGs stored for the current patient.

Note: You must be logged in to the PDMS to access the ECoG library. If you are not, the NEUROPACE LOG IN screen opens when you select ECoG LIBRARY, prompting you to log in. (See Log In to the PDMS on page 16.) If you are not connected to the Internet, see Connect to the Internet for PDMS Access on page 11.

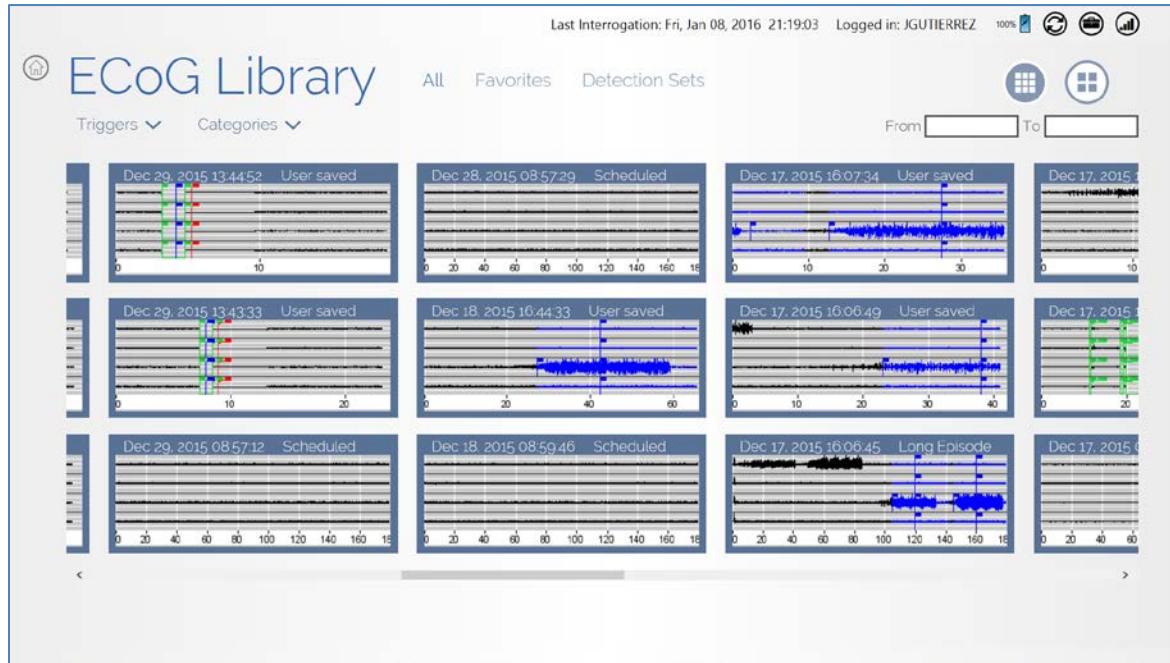


Figure 21 ECoG Library screen

Each ECoG is identified by date and time captured and by capture trigger. Newest ECoGs are on the left, older ones toward the right. Swipe left or right to view all ECoGs. You can filter which ECoGs are displayed using **TRIGGERS** and **CATEGORIES** at upper left (**Figure 18**). Use the **FROM** and **To** fields at upper right to view ECoGs by date range.

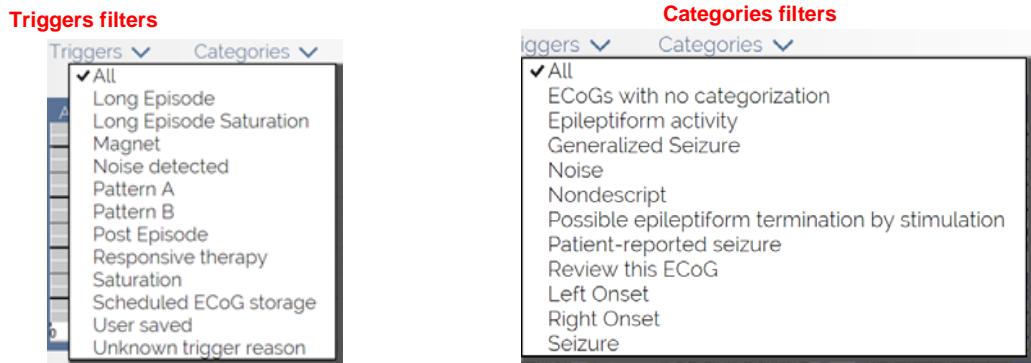


Figure 22 Filter stored ECoGs by triggers and categories

At top center (see **Figure 18**), you can select **ALL**, **FAVORITES** or **DETECTION SETS** to display ECoGs by those classifications.

Select a thumbnail to open it in a larger size, as shown below. ECoGs highlight detection and stimulation events. Detected activity appears in blue. Undetected activity appears in black. Colored flags mark specific points of detection or stimulation. You can **ZOOM** in on the ECoG and increase the **GAIN**. Use the **ASSIGN CATEGORY** field to assign the ECoG to a category. Select the **STAR** button at upper right to specify this ECoG as a favorite, for fast access later using the **FAVORITES** filter on the **ECoG LIBRARY** screen.

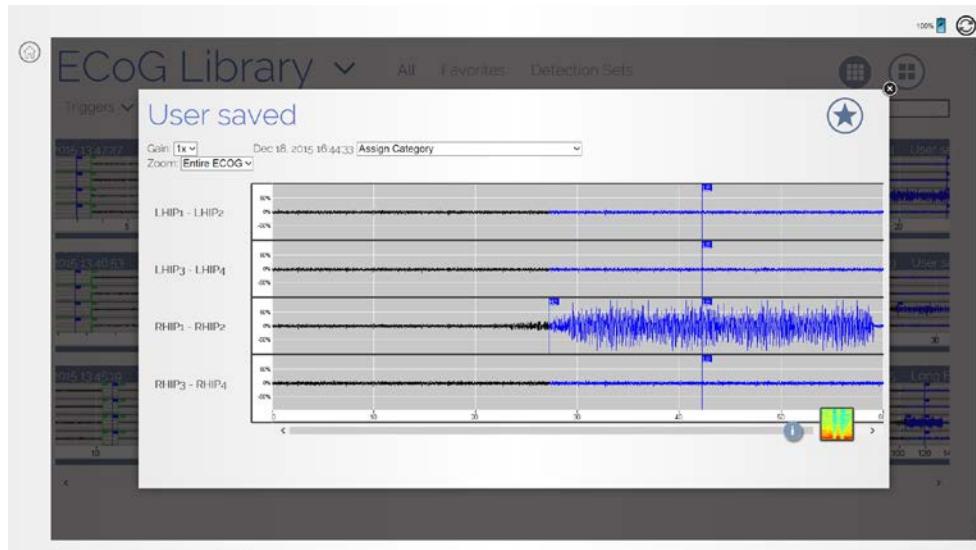


Figure 23 Viewing a User saved ECoG

Captured ECoGs reside in the RNS® Neurostimulator memory until they are transferred to the programmer at the next interrogation with the wand, and are then cleared from memory. If neurostimulator memory capacity is reached, captured ECoGs are overwritten by subsequent ECoGs. The number of ECoGs that can be stored depends on the duration of each ECoG, which you can specify in the **CAPTURE WINDOW** field in the **DURATION** area on the far upper right side of the **ECoG CAPTURE** screen. See **Selecting ECoG Capture Settings** on page 44 for details.

LIVE ECoGs SCREEN

Note: You can access the **LIVE ECoGs** screen only after interrogation of the neurostimulator, and the wand must be held within approximately 1 inch of the neurostimulator to acquire live ECoGs.

After interrogation, keep the wand in place and select **LIVE ECoGs** from the **HOME** screen to open the **LIVE ECoGs** screen. An ECoG is live when you open it. Buttons provide various options for viewing the live ECoG. Select **STOP/STORE** to save the ECoG in the patient's ECoG library. It remains on screen after you stop and store it until you start another live ECoG. Captured live ECoGs are stored in the **ECoG LIBRARY** and classified as *User saved*.



Figure 24 Live ECoGs screen

Live ECoG capture can be useful for immediate or later analysis of patient ECoGs, for adjusting amplifier gain settings, and for simulated detection as part of **Configure Pattern Detection: Define and Modify Detection Settings** (see page 50).

*Important Note: The recording montage must be programmed before live ECoGs can be viewed. See **Create the Recording Montage** on page 38.*

Note: The maximum length for a live ECoG is four minutes. If the live ECoG is allowed to scroll for more than four minutes, only the last four minutes are stored.

Note: Live ECoGs are stored on the tablet, not in the neurostimulator memory, and are uploaded to the PDMS if and when the tablet is logged into the PDMS.

Caution: **Telemetry Artifact**

Telemetry may produce an electrographic artifact. If responsive therapy is enabled with a sensitive detection set, detection of the electrographic artifact may occur, resulting in therapy delivery. The physician should be aware of possible sensing artifacts when assessing the ECoG recordings.

CREATE THE RECORDING MONTAGE

The RNS® Neurostimulator can sense brain activity (ECOGs) from up to four amplifier channels. The programmed configuration of leads and electrodes you assign to each channel constitute the recording montage. You create the montage using the tablet. **There can be only one current montage and corresponding detection settings. When the montage is altered, the detection settings created with the previous montage are discarded and the detection settings must be reprogrammed.**

Identifying Ports, Leads and Electrodes on the Neurostimulator

The figure below identifies neurostimulator **PORT #1** and **PORT #2** with their leads extending out of the ports.

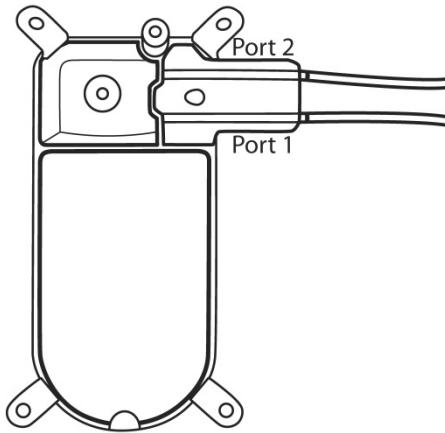


Figure 25 RNS® Neurostimulator port #1 and #2 locations. Port #1 is closer to the center of the neurostimulator and port #2 is closer to the edge.

For both depth leads and cortical strip leads, the electrodes are numbered 1 through 4 starting from the distal end of the lead.

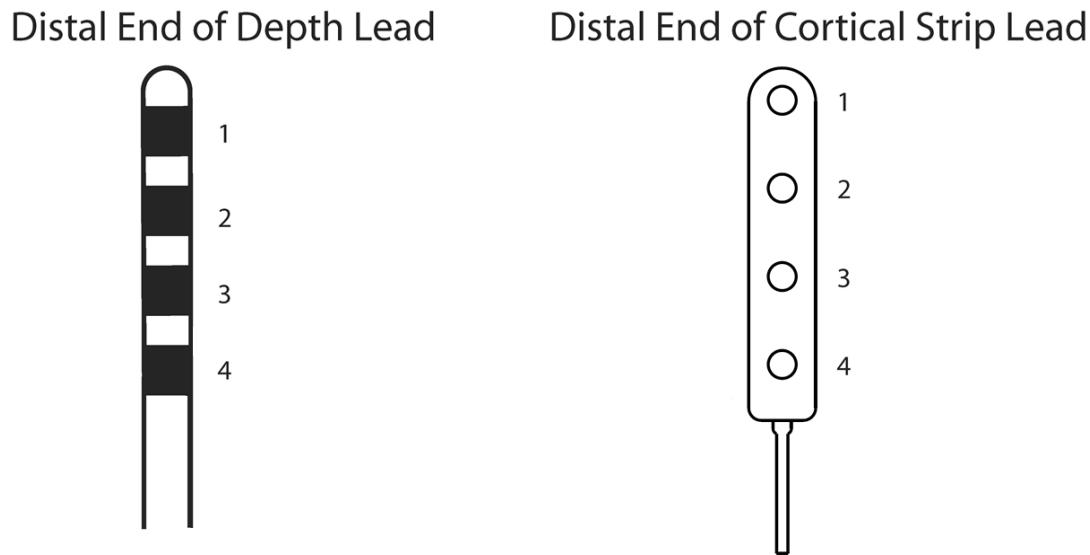



Figure 26 For both lead types, electrodes are numbered 1 through 4. Electrode 1 is most distal.

Creating the Montage



1. From the **HOME** screen, select **NEUROSTIM INFO**. The **NEUROSTIM INFO** screen opens.

Note: For all settings, changed settings are indicated by blue buttons. Gray buttons indicate settings have not changed.

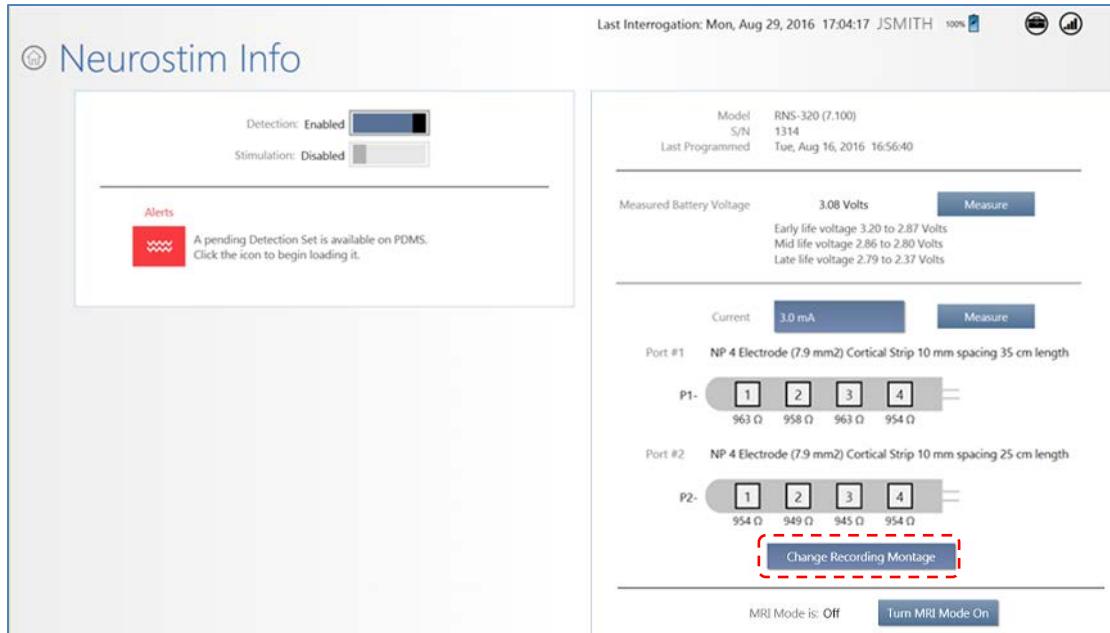

- On the NEUROSTIM INFO screen, select CHANGE RECORDING MONTAGE at lower right. The RECORDING MONTAGE screen opens.

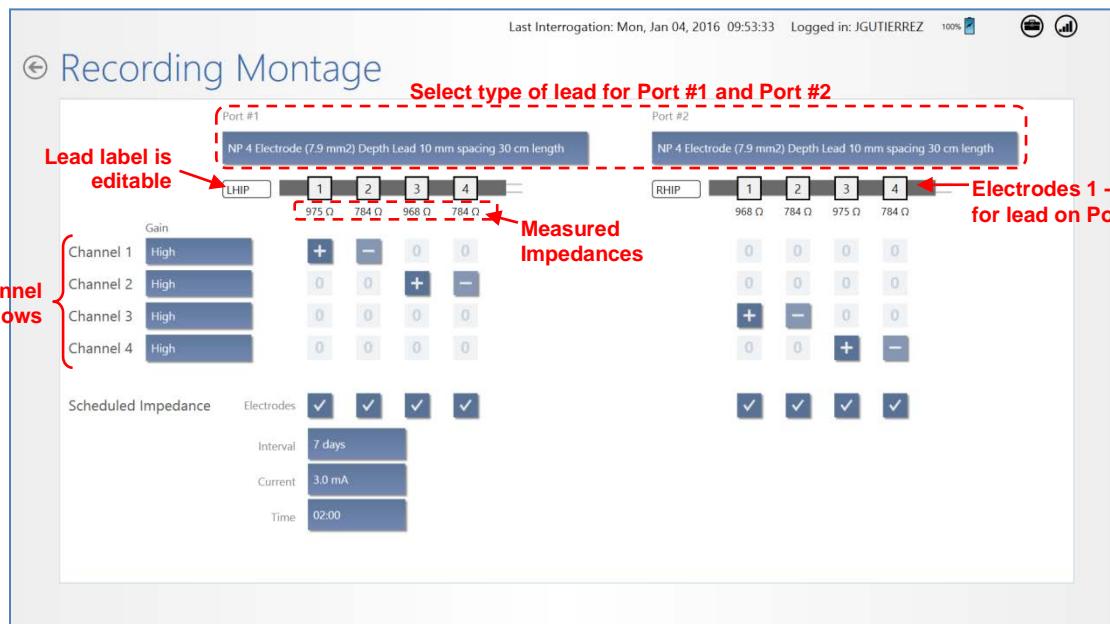


Figure 27 Recording Montage screen

- At Neurostimulator Implant and When Changing Leads – Select Lead Type: Upon first use of leads, you must select the lead type for PORT #1 and/or PORT #2. To do this, select the field that says NONE, just below the PORT #1 or PORT #2 names at the top, and select the type of lead installed in each port. When lead types are already selected for the ports, proceed to the next step.

4. **Assign Channels:** Use the channel rows and electrode columns to assign two electrodes to each **CHANNEL** (1 through 4). To do this, select the small, square field where the channel and electrode intersect and then select a positive [+] non-inverting input, a negative [-] inverting input, or zero [0]. Detection occurs across the selected electrodes for each channel.

Each electrode can be assigned to one or two channels. The tablet requires you to use a positive [+] input for electrodes assigned to two channels.

Note: The tablet does not permit use of a negative [-] input for electrodes assigned to two channels because it may result in ECoG artifacts.

5. **Set Gain:** Next to the channel name on the left, select the **GAIN** field for each channel to select **HIGH**, **MEDIUM-HIGH**, **MEDIUM-LOW**, or **LOW**. The default is **HIGH** and it is usually not changed. After programming, you can assess ECoG signal amplification by viewing the live ECoGs.
6. **Scheduled Impedance Settings:** These settings are at lower left. By default, the tablet schedules an impedance test on all electrodes every 7 days with 3.0 mA current at 2:00 AM. Select the fields next to **INTERVAL**, **CURRENT** and **TIME** to choose different settings. Normally there is no need to use other than the default settings.

Recommended Initial Montage Settings

- **Channel 1** = Lead 1: Electrode 1 [+], Electrode 2 [-]
- **Channel 2** = Lead 1: Electrode 3 [+], Electrode 4 [-]
- **Channel 3** = Lead 2: Electrode 1 [+], Electrode 2 [-]
- **Channel 4** = Lead 2: Electrode 3 [+], Electrode 4 [-]
- **GAIN** set to **HIGH** for all channels (default)
- **SCHEDULED IMPEDANCE** settings = Interval: 7 days, Current: 3.0 mA, Time: 2:00 AM (defaults)

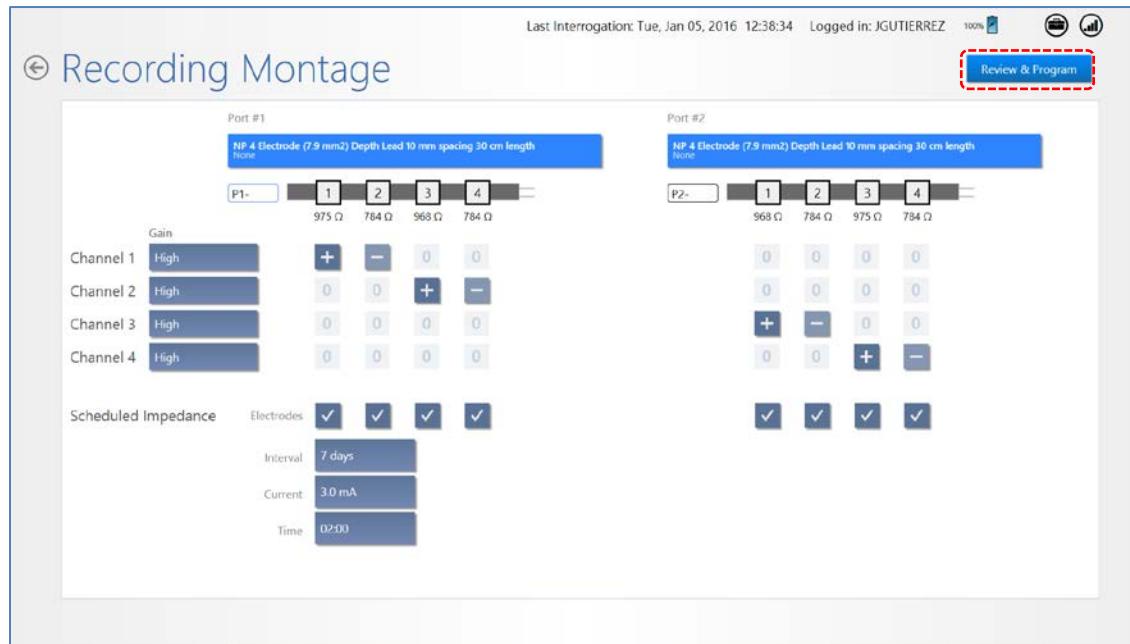



Figure 28 Recommended initial montage settings

Program the Neurostimulator with the New Settings

The new settings are not sent to the neurostimulator until you complete this step.

*Note: ECoG capture triggers are inactive during programming. **Programming clears activity information and ECoGs stored in the neurostimulator.** If there is any new information, be sure to interrogate the neurostimulator prior to programming.*

7. Select **REVIEW & PROGRAM** at upper right (**Figure 25**).

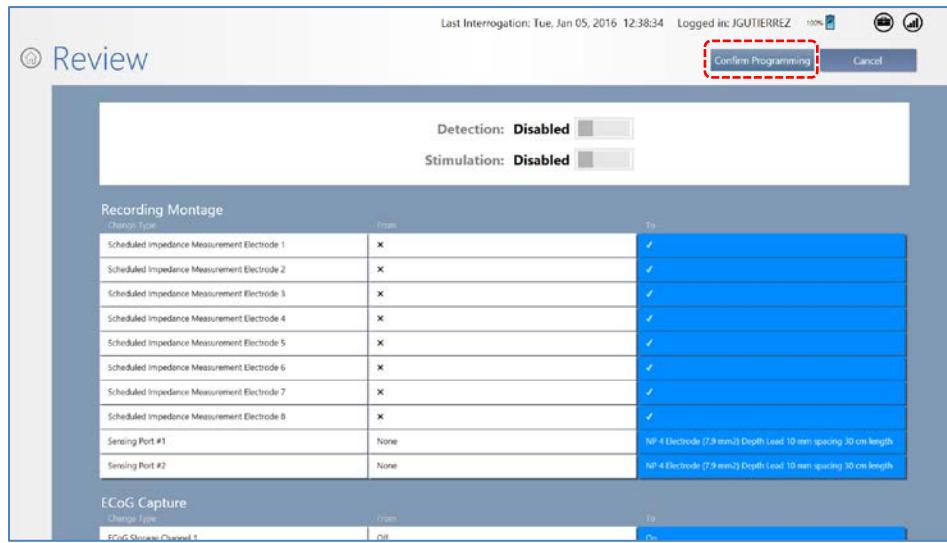

8. On the **REVIEW** screen, if you are satisfied with the new or changed settings, place the wand within approximately 1 inch of the neurostimulator and select **CONFIRM PROGRAMMING**. If not satisfied, select **CANCEL** to return to the previous screen.

The **PROGRAMMING NEUROSTIMULATOR** dialog shows programming progress until complete.


SET UP ECoG CAPTURE

Electrocorticograms (ECoGs) can be stored by the RNS® Neurostimulator when an ECoG capture trigger occurs. The maximum number of ECoGs that can be stored is affected by the number of channel(s) selected ON for storage, and the pre-trigger and post-trigger durations. The oldest stored ECoGs are continually replaced by the newly stored ECoGs, except as noted below.

The number of ECoGs stored for a particular ECoG capture trigger can be reserved. Reserved ECoGs are only overwritten by newer ECoGs of the same ECoG capture trigger type.

ECoG capture may also be scheduled at four different times on a 24-hour clock.

Selecting ECoG Capture Settings

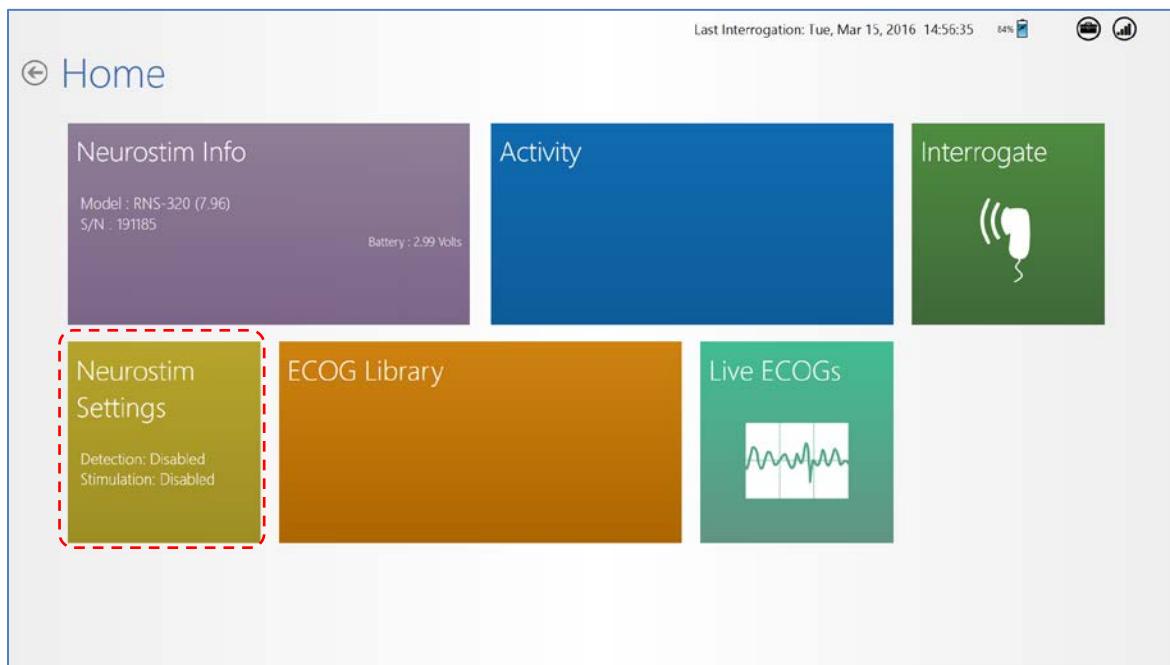

1. Select **NEUROSTIM SETTINGS** from the **HOME** screen.

2. On the NEUROSTIM SETTINGS screen, select ECOG CAPTURE at upper right. The ECOG CAPTURE screen opens.

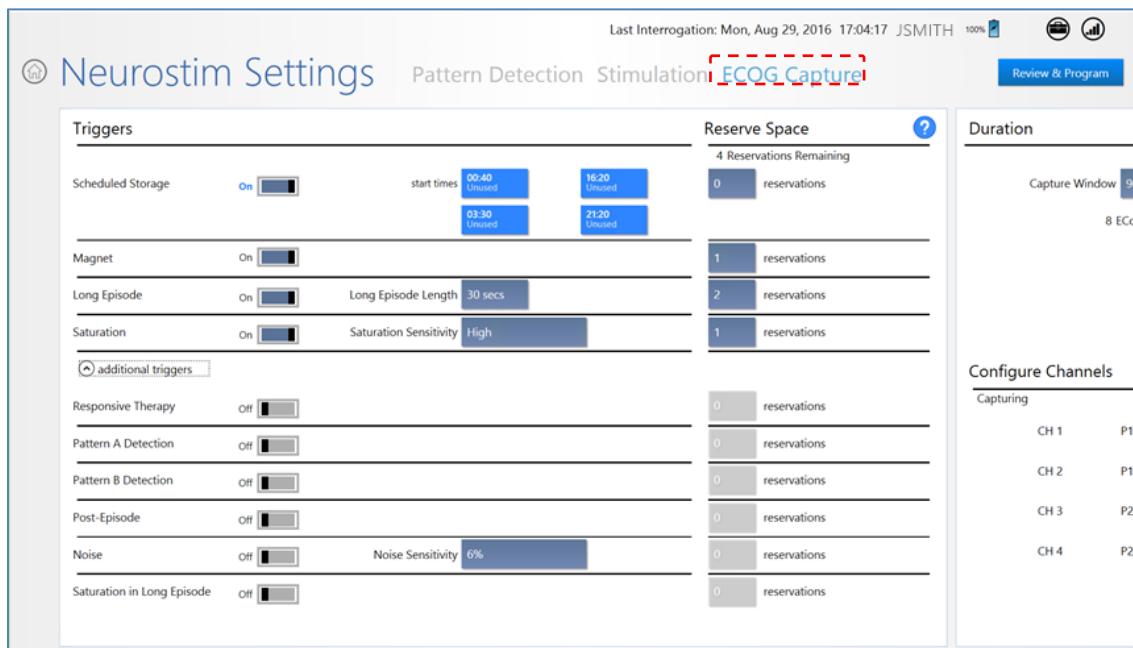


Figure 29 ECoG Capture screen, left side. Swipe left to see right side.

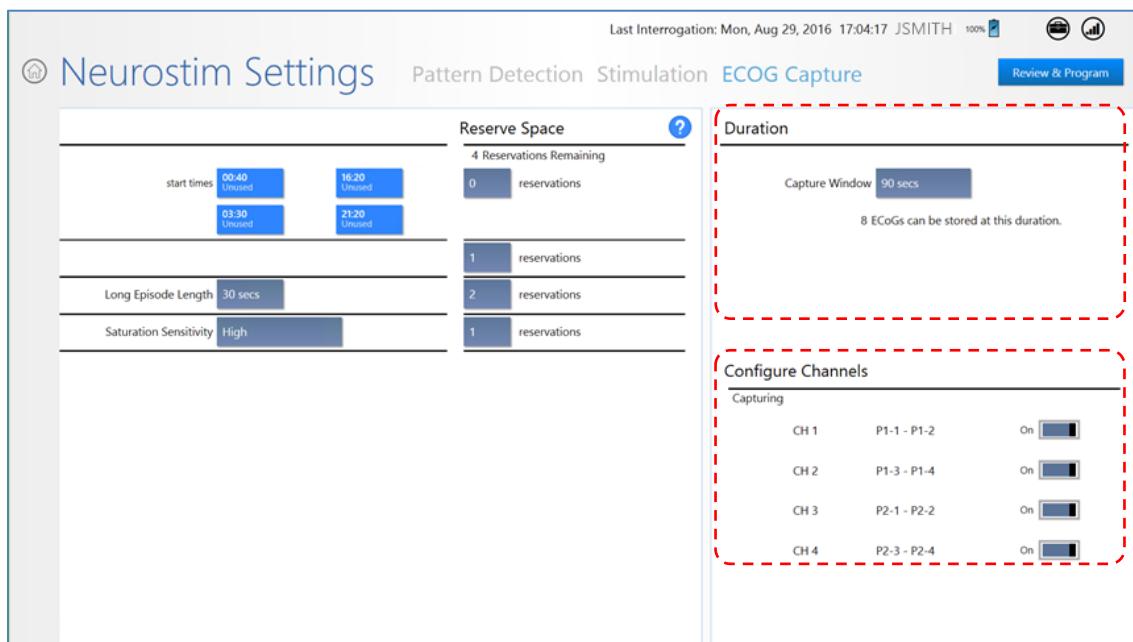


Figure 30 Right side of ECoG Capture screen

3. **Turn On Desired Channels:** In the **CONFIGURE CHANNELS** area at lower right, select the **ON/OFF** switch next to the desired **CHANNEL(S)** (**CH 1** through **CH 4**) to turn each channel on or off as desired. Only the channels selected **ON** will be recorded and stored.
4. **Set Capture Duration:** In the **DURATION** area at upper right (*Figure 27* above), select the **CAPTURE WINDOW** field and select one of the capture duration options. After your selection, the **DURATION** area shows the total duration (in seconds) and the number of ECoGs that can be stored at this duration.
5. **Set Capture Triggers:** In the **TRIGGERS** area at upper left (*Figure 26* above and *Figure 28* below), select the **ON/OFF** switch to turn **ON** the triggers that you wish to initiate ECoG capture.



Figure 31 Selecting ECoG Capture Triggers

There are three recommended ECoG triggers to be programmed at implant. The first is a magnet swipe by the patient, which could indicate that a clinical seizure has occurred. The others are saturation (high amplitude ECoG) and long episode (sustained changes in the ECoG), both of which could indicate an electrographic seizure.

Commonly Used Triggers: The more commonly used triggers are listed first.

- **SCHEDULED STORAGE:** Initiate ECoG capture at up to four specific start times, which you can select using the fields to the right. Options are available in ten minute intervals using 24 hour time.
- **MAGNET:** Initiates ECoG capture when the patient swipes the magnet over the neurostimulator. If the patient swipes the magnet when they experience clinical symptoms, selecting this trigger may be helpful in evaluating these events.
- **LONG EPISODE:** Initiates ECoG capture when a detected episode exceeds the duration you select in the **LONG EPISODE LENGTH** field to the right. Duration options range from **1** to **120** seconds.

- **SATURATION:** Initiates ECoG capture when the ECoG amplitude exceeds the threshold associated with the value you select in the **SATURATION SENSITIVITY** field to the right. The sensitivity options are **High**, **Medium**, **Low** and **Off**. Selecting this trigger may be helpful in troubleshooting the RNS® System.
- Additional Triggers:** Select the down arrow next to **ADDITIONAL TRIGGERS** at lower left to access the triggers below.
- **RESPONSIVE THERAPY:** Initiates ECoG capture when electrical stimulation therapy is delivered in response to a detected event. Selecting this trigger may be helpful in determining if therapy was delivered appropriately as well as evaluating the result of therapy delivered.
 - **PATTERN A DETECTION** and **PATTERN B DETECTION:** These are separate triggers; each initiates ECoG capture when a detected episode begins with the patterns designated as Pattern A and Pattern B on the **PATTERN DETECTION** screen. Selecting a pattern trigger may be helpful in evaluating the pattern detection settings.
 - **POST-EPIISODE:** Initiates ECoG capture at the end of a detected episode. Selecting this trigger may be helpful in evaluating activity after an episode.
 - **NOISE:** Initiates ECoG capture when the ECoG signal exceeds 60 Hz. Selecting this trigger may be helpful in troubleshooting the RNS® System.
 - **SATURATION IN LONG EPISODE:** Initiates ECoG capture when the ECoG amplitude exceeds the threshold selected in the **SATURATION SENSITIVITY** field during an episode that exceeds the duration selected in the **LONG EPISODE LENGTH** field.
6. **Set Reservations for Each Trigger:** In the **RESERVE SPACE** column to the right of each trigger (see *Figure 28* above), select the **reservations** field and select the number of reservations to make for each trigger. A reservation means that up to the specified number of ECoGs will be stored for that trigger, and can only be overwritten by newer ECoGs of the same trigger type. The total number of reservations available depends on the duration selected in the **DURATION** area at upper right.

Recommended Initial ECoG Capture Settings

- **Configure Channels for Capturing:** All four channels **ON**
- **Duration: 90 SECS**
- **Triggers On for ECoG Capture:**
 - **LONG EPISODE:** Reservations = **2**, Long Episode Length = **30 SECS**
 - **MAGNET:** Reservations = **1**
 - **SATURATION:** Saturation Sensitivity: **Low**, Reservations = **1**
 - All other triggers **OFF**

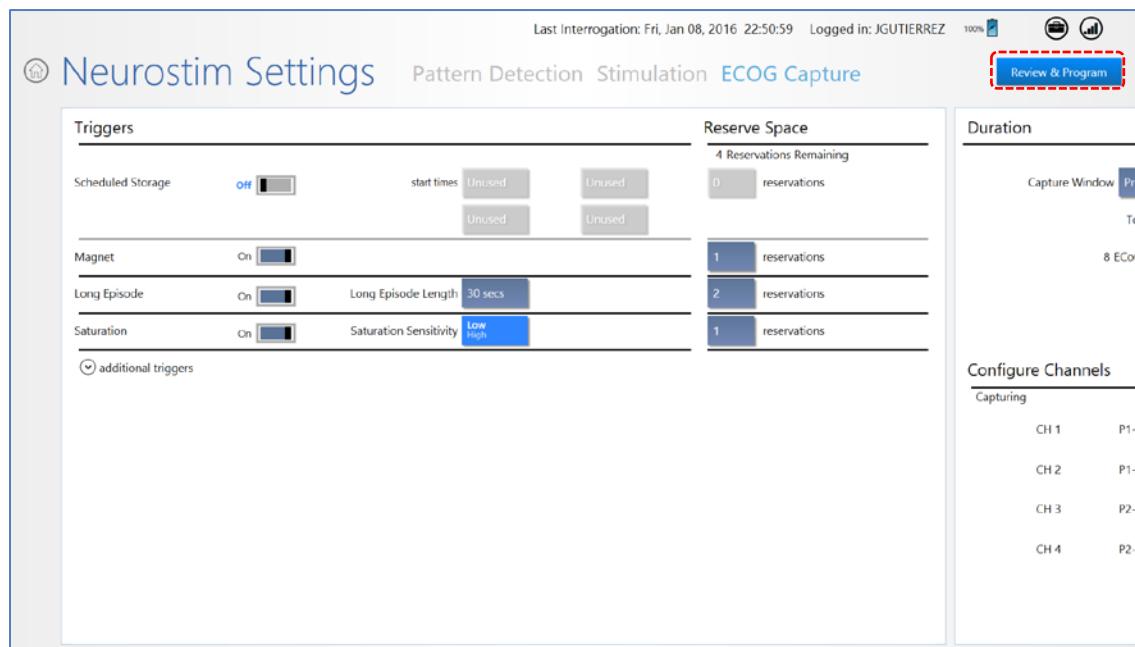



Figure 32 ECoG Capture showing recommended initial settings

Program the Neurostimulator with the New Settings

The new settings are not sent to the neurostimulator until you complete this step.

*Note: ECoG capture triggers are inactive during programming. **Programming clears activity information and ECoGs stored in the neurostimulator.** If there is any new information, be sure to interrogate the neurostimulator prior to programming.*

7. Select **REVIEW & PROGRAM** at upper right. The **REVIEW** screen opens. It lists all changes you are about to program.

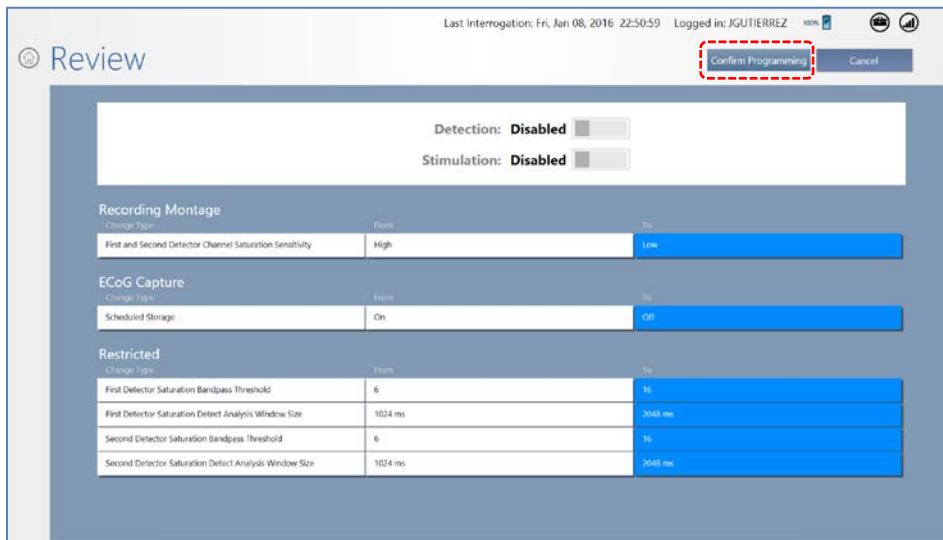

8. On the **REVIEW** screen, if you are satisfied with the new or changed settings, place the wand within approximately 1 inch of the neurostimulator and select **CONFIRM PROGRAMMING**. If not satisfied, select **CANCEL** to return to the previous screen.

The **PROGRAMMING NEUROSTIMULATOR** dialog shows programming progress until complete.

RECOMMENDATION FOR REVIEW OF INITIAL ECoG CAPTURE SETTINGS

ECoGs stored by the RNS® Neurostimulator during the first several days after implantation of the neurostimulator and leads may not be representative of the patient's typical baseline activity. ECoG capture settings should be reviewed to ensure that the desired activity is being stored. If it appears that ECoGs containing the desired activity are being overwritten by later events, instruct the patient to use the remote monitor and wand to interrogate the neurostimulator immediately after each clinical seizure, or to interrogate several times per day. Several interrogations per day during the first several days after implantation will not have a significant impact on battery life.

- **Recommendation:** Review ECoGs on the PDMS to determine whether detections and storage are optimal. If not, consider adjustments to the following triggers recommended to be **ON** initially.
 - **LONG EPISODE LENGTH:** Adjust the duration that triggers capture of a long episode.
 - **MAGNET:** Consider turning **OFF** the **MAGNET** trigger if the patient is not able to apply the magnet before or immediately after a clinical seizure, or if ECoGs captured by this trigger are not useful. Magnet swipes will still be counted and reported in the diagnostics.

- **DURATION:** Select a briefer **CAPTURE WINDOW** value to yield more stored ECoGs, or a lengthier value to yield fewer stored ECoGs.

If you make changes, be sure to follow the steps above to review and program the new settings into the neurostimulator.

CONFIGURE PATTERN DETECTION: DEFINE AND MODIFY DETECTION SETTINGS

Pattern detection settings are used by the RNS® Neurostimulator to recognize epileptiform activity as defined by the physician. To configure pattern detection means to select specific ECoG patterns to be detected, and then to program the system to detect these patterns. Configuring pattern detection is a prerequisite for electrical stimulation therapy, but it does not require you to enable therapy.

Note: A pattern defines the type of ECoG activity to be detected.

Note: When defining or modifying a detection set, the programmer and the PDMS are able to automatically simulate detection on all ECoGs stored for a patient. Simulation is displayed for groups of six ECoGs at a time, with the user able to navigate through all ECoGs for that patient.

Note: Although detection sets can be created and modified using the PDMS either in a web browser or on the programmer, only the programmer can be used to program the new settings into the neurostimulator. When you select a new detection set that was created on the PDMS, it is not necessary to re-run simulation on the programmer.

Configuring pattern detection has three phases:

- **Phase 1: Select a Default Detection Set upon First Use** on page 51, you select one of the three default detection sets to gather an initial collection of ECoGs. This is usually configured and programmed in the operating room immediately after implantation.
- **Phase 2: Define a Pattern Detection Set** on page 54, after initial ECoGs have been stored, you review the stored ECoGs and identify the specific ECoG patterns to be detected. Using the SimpleStart Feature, the system provides a starting point for this process. It is a tool to assist you as you exercise clinical judgment to define detection settings. It suggests a starting point only. It does not preempt nor replace your clinical judgment. This is usually programmed at the first visit post-implantation.
- **Phase 3: Customize Pattern Detection** on page 58, you repeat the process of reviewing stored ECoGs and adjusting the parameters of each pattern to optimize detection of ECoG activity for this patient. This is programmed at any visit after phase 2.

Note: For all settings, changed settings are indicated by blue buttons. Gray buttons indicate settings have not changed.

To put these three phases in context:

Immediately after implantation of the neurostimulator, you select a default detection set in **Phase 1**, in order to collect ECoGs with activity of interest specific to the patient, and ECoGs with baseline activity.

After collecting sufficient ECoGs, you identify the specific electrocorticographic patterns to be detected and then, in **Phase 2**, the system provides a starting point for establishing a patient-specific detection set. You may then modify the parameters of the suggested starting detection set to make detection more or less sensitive, or to detect earlier or later.

In future visits, after collecting additional ECoGs, you then, in **Phase 3**, iterate to improve detection (and, once responsive therapy has been enabled, to improve clinical response) by making small, incremental changes to detection parameters with the sliders or **MORE CONTROLS**.

Note: Typically, physicians repeat earlier phases (Phase 1 or Phase 2) only if they observe a new type of ECoG activity of interest that requires an additional detector; or if a good clinical response has not been obtained and they wish to start over.


Phase 1: Select a Default Detection Set upon First Use

Upon first use of the tablet for each patient, the patient has no saved ECoGs from which to select a pattern for detection, and so you must select one of three default detection sets. This is usually configured in the operating room immediately after implantation. The default sets are used to gather an initial set of ECoGs for a patient, which you can then use in a subsequent visit to define a pattern detection set specific to the patient. If you have already performed these steps once, proceed to the next section, **Phase 2: Define a Pattern Detection Set**.




1. Select **NEUROSTIM SETTINGS** from the **HOME** screen.

2. On the **NEUROSTIM SETTINGS** screen, select **PATTERN DETECTION** at upper center. A list of three **COMPATIBLE DETECTION SETS** appears.

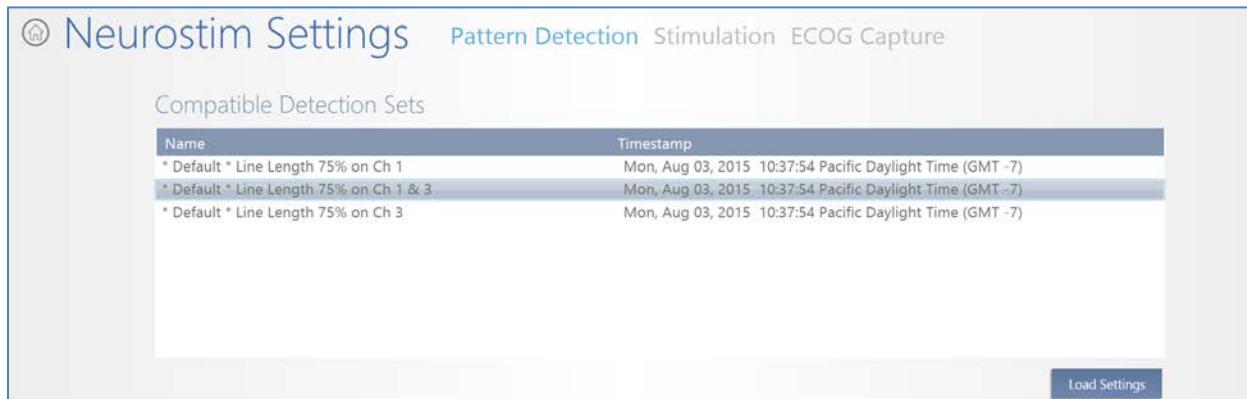

3. Select one of the three default detection settings. In the usual case, a patient has two active leads to be configured, and therefore the recommended default pattern is **POWER CHANGE 75% ON CH 1 & 3 (DEFAULT)**, the second item in the list. (If only one lead is implanted, this two-channel option will not appear in the list.) The other default patterns detect on only one lead.
 - Power change 75% is a detection setting that captures any large changes in the signal that clinicians can program and optimize as detection patterns.
4. Select **LOAD SETTINGS** at lower right, then **REVIEW & PROGRAM** at upper right on the next screen. The **REVIEW** screen opens.

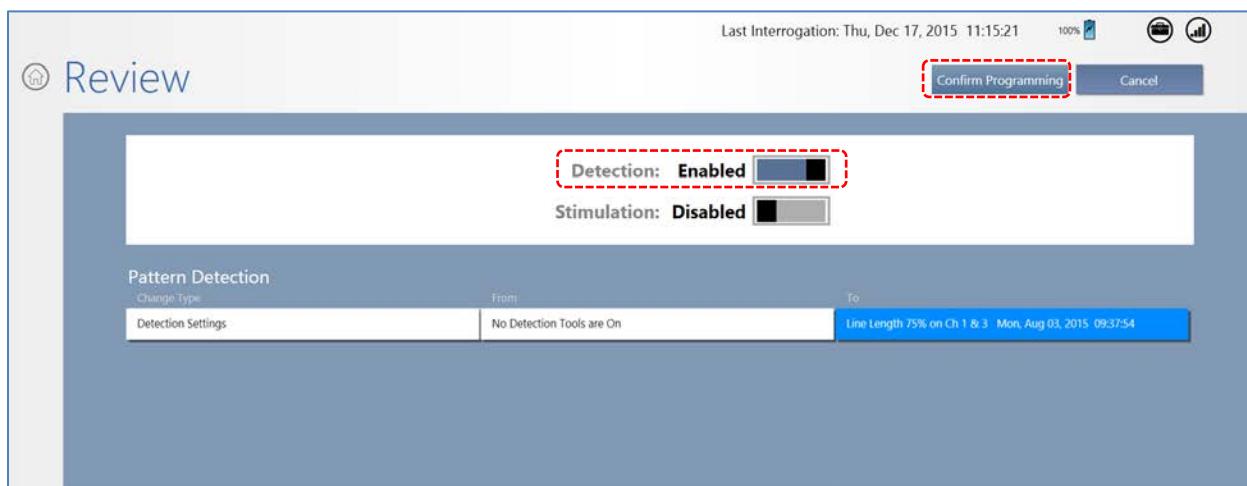

5. On the **REVIEW** screen, make sure **DETECTION** is **ENABLED**. Then, if you are satisfied with the new settings, place the wand within approximately 1 inch of the neurostimulator and select **CONFIRM PROGRAMMING**. If not satisfied, select **CANCEL** to return to the previous screen.

When you confirm programming, a **TELEMETRY INFORMATION** dialog reports that “Detection is Enabled and Responsive Therapy is disabled,” and prompts you to select **CONTINUE** or **CANCEL**.

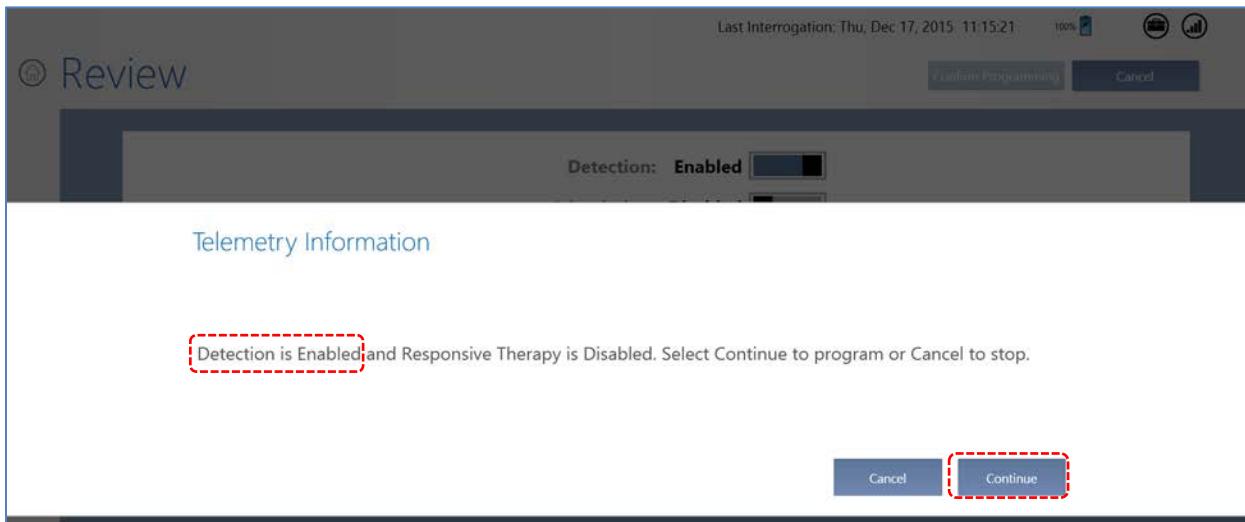


Figure 33 Telemetry Information dialog to confirm programming again

6. Make sure the message reports that **DETECTION** is **ENABLED**, keep the wand in place and select **CONTINUE** to proceed with programming, or **CANCEL** to cancel programming and return to the previous screen.

When you select **CONTINUE**, the **PROGRAMMING NEUROSTIMULATOR** dialog shows programming progress.

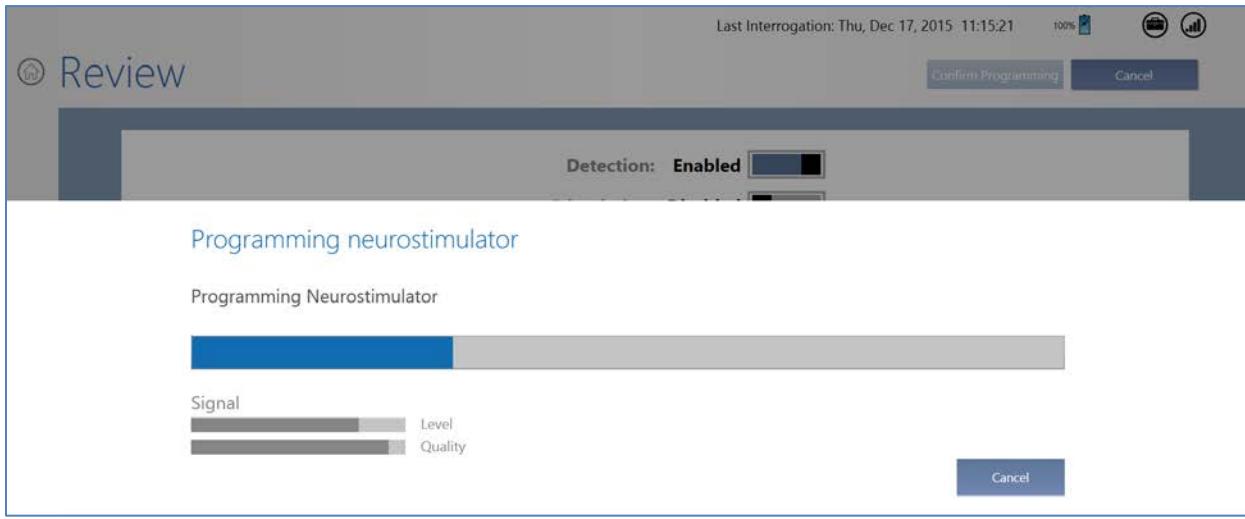


Figure 34 Programming neurostimulator dialog shows progress of programming

When complete, you return to the list of **COMPATIBLE DETECTION SETS**. Initial configuration is complete. (You can select a different default setting and redo initial configuration if desired.)

Phase 2: Define a Pattern Detection Set

Note: You must be connected to the PDMS (via the Internet) to define a pattern detection set. If not, you can program only a default detection set or a detection set previously defined and downloaded to the tablet. If you are not connected to the Internet, see **Connect to the Internet for PDMS Access** on page 11.

At this point, the patient should have ECoGs stored in the PDMS from which to select specific patterns for detection. A set of patterns you define and name is called a **pattern detection set** (or simply a detection set). A pattern defines the type of ECoG activity to be detected. You can define up to four patterns in a set, two patterns on each of two channels. The patterns defined on a specific channel can operate only in that selected ECoG channel. Follow these steps to define a pattern detection set.

1. Select **NEUROSTIM SETTINGS** from the **HOME** screen,
2. On the **NEUROSTIM SETTINGS** screen, select **PATTERN DETECTION** at upper center. The **PATTERN DETECTION** screen opens.

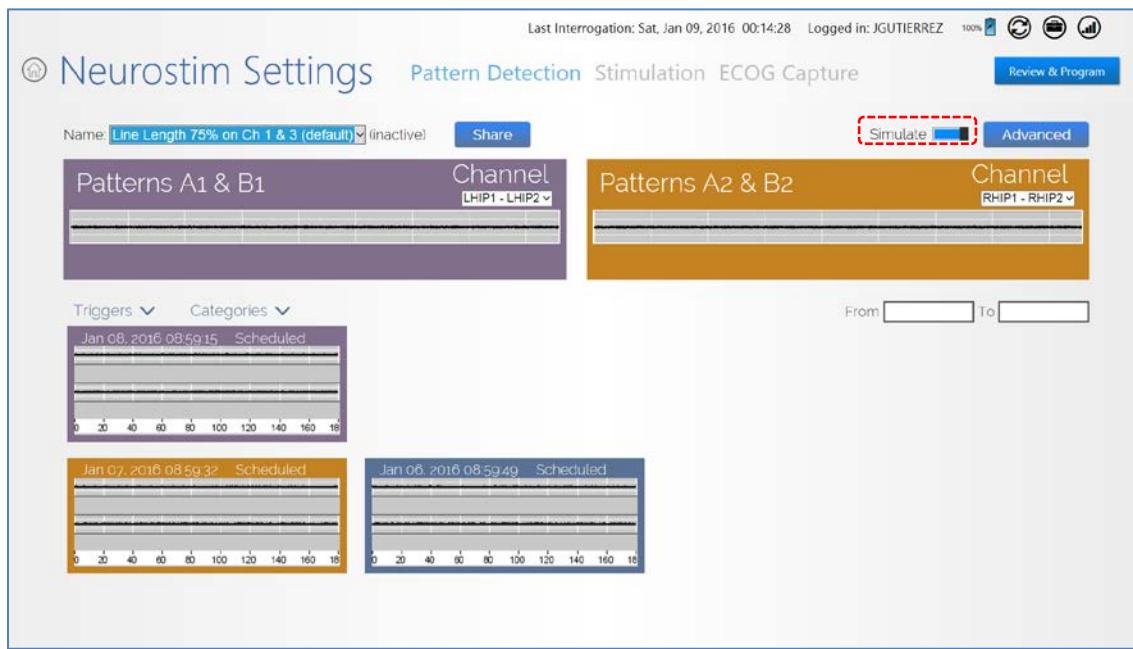


Figure 35 Pattern Detection screen

- The **NAME** field at upper left identifies the currently active detection set, as indicated by **(ACTIVE)** next to the field. You can edit this set or select an inactive set to program by selecting the **NAME** field.
- The upper half of the screen provides tiles to select the current detection patterns and channels; the tiles are named **PATTERNS A1 & B1** (upper left) and **PATTERNS A2 & B2** (upper right).
- The lower half of the screen shows thumbnails of all ECoGs (if no filters are on) stored in the PDMS that were captured using the current recording montage. (ECoGs stored in the PDMS were uploaded after interrogations from the remote monitor or programmer.) You can filter

- the ECoGs shown using the **TRIGGERS** and **CATEGORIES** menus at left center, and/or by date range using the **FROM** and **TO** fields at middle right.
- The **SIMULATE** button at upper right (next to **ADVANCED**) is switched **ON** (blue) by default, which enables simulations of detections to be shown on all ECoGs shown on this screen.
 - 3. To define or change a detection set, select the desired detection set in the **NAME** field at upper left.
 - The detection set may be one of the three default detection sets to be used as a starting point for configuration; or it may be a detection set created and saved previously.
 - If you don't want to modify the detection set you chose, you can program the detection set right away: Select **REVIEW & PROGRAM** at upper right and, while holding the wand within approximately 1 inch of the neurostimulator, follow the prompted steps to **CONFIRM PROGRAMMING** and **CONTINUE** sending the new settings to the neurostimulator.
 - 4. Use the **DEFINE PATTERNS** dialog to define or modify a detection set. There are two routes to the **DEFINE PATTERNS** dialog.
 - Select one of the two **PATTERN** tiles at the top to go there directly, or
 - Select one of the ECoG thumbnails at the bottom to open its ECoG dialog, which is named by its date and trigger type, and then select one of the two **MODIFY** buttons (**A1 & A2**, **B1 & B2**) at upper right. (You can also use the ECoG dialog to assign the ECoG to a category if desired.)

If You Selected a New ECoG First

5. If you selected a new ECoG first and selected a **MODIFY** button, a **SELECT CHANNEL** dialog opens. Touch within the channel that has a pattern to detect. Then the **DEFINE PATTERNS** dialog opens.

Using the Define Patterns Dialog

Note: A pattern defines the type of ECoG activity to be detected.

6. In the **DEFINE PATTERNS** dialog, the upper portion shows the ECoG in the selected channel. Touch on the point of the ECoG where you would like detection to occur. Starting from the point where you touch, a 3-second window of time (or region of interest, ROI) is bracketed on the ECoG. This touch feature to bracket a 3-second region of interest is called the SimpleStart Feature.

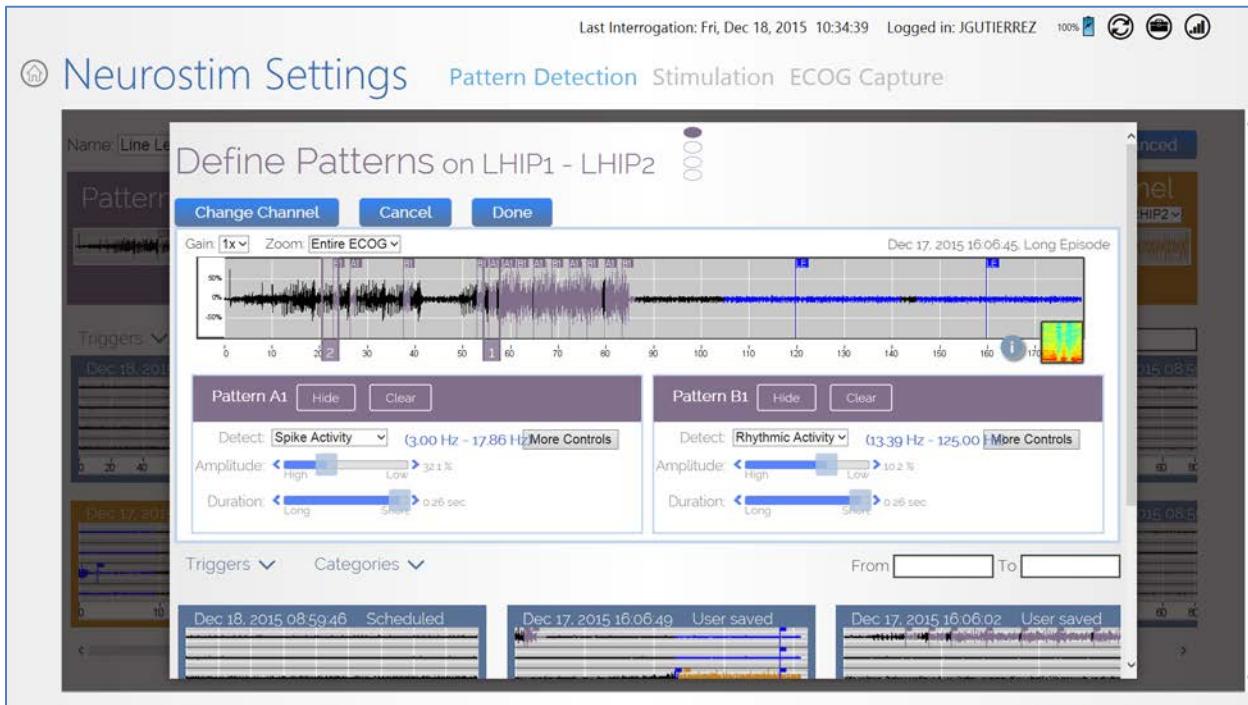


Figure 36 Define Patterns dialog

- The 3-second window bracketed on the ECoG is the pattern you are identifying for detection; it should contain the type of pattern you would like to detect.
- Each 3-second window is identified as one of the four possible patterns you can define, **A1, A2, B1 or B2**. (The 1 and 2 indicate the first and second ROI you select. Typically the A and B patterns are defined on different leads.)
- You can touch on the window identifier (**A1, A2, B1 or B2**) and slide the window left or right to place it where it gives the desired simulated detection results. Note that the software may apply a different suggested detection tool when you move the window.

Note: You cannot select a 3-second window within the first 5 seconds of an ECOG, because the software looks for differences in the ECOG data before the 3-second window and needs sufficient preceding data to determine this difference. If you touch before 5 seconds in the ECoG, the software places the window starting at the 5-second mark.

- Normally you want to touch near the beginning of the pattern in order to define a pattern that will be detected promptly after it begins. Detection settings should promptly detect ECoG activity and not detect on baseline activity.

The software automatically applies a suggested detector (shown in the **DETECT** field) and settings that fit the waveform pattern in that 3-second window (or region of interest, ROI), and simulates detection, showing the detection flags the chosen pattern would apply to this ECoG—and to all the patient's other stored ECoGs shown in smaller size below. It is recommended to review simulated detection activity on an ensemble of stored ECoGs before finalizing the detection settings.

- To hide the simulation on the ECoG, select **HIDE**. To remove the selection on the ECoG, select **CLEAR**.
- At any point, you can select a different ECoG by selecting an ECoG thumbnail along the bottom; swipe left or right to view all available thumbnails.
- By default, if you define both patterns, detections will occur for both patterns on this channel. The channel is indicated by the filled oval at top center of the dialog, from 1 to 4 top to bottom. Use the **CHANGE CHANNEL** button to change the channel on which to detect the pattern.

Note: If simulating performance of a 2-channel detection set, be sure to select ECoGs that contain ECoG activity patterns you wish to detect in both channels.

- Make adjustments to the parameters applicable to the detection tool using the sliders available and observe the effects on the simulation. For **SPIKE ACTIVITY** and **RHYTHMIC ACTIVITY**, the sliders adjust **AMPLITUDE (%)** and **DURATION** (seconds). For **POWER CHANGE**, the sliders adjust **SENSITIVITY (%)** and **DETECT** (seconds).
 - All these detection parameters affect sensitivity. Adjustments that increase sensitivity tend to increase how early and how often the pattern is detected; adjustments that decrease sensitivity tend to decrease how early and how often the pattern is detected.
 - Detection settings should be programmed for early detection of activity. Detection settings that result in early activity detection usually also result in numerous detections; this is expected and acceptable behavior.
 - The parameter values are shown to the right of each slider. For all four parameters, lower values result in greater sensitivity.
 - The **MORE CONTROLS** button gives access to additional adjustment options and **Advanced Controls** (see page 69).
7. Select **DONE** at upper right to close the **DEFINE PATTERNS** dialog. The detection set is not saved until you save it in the next step.

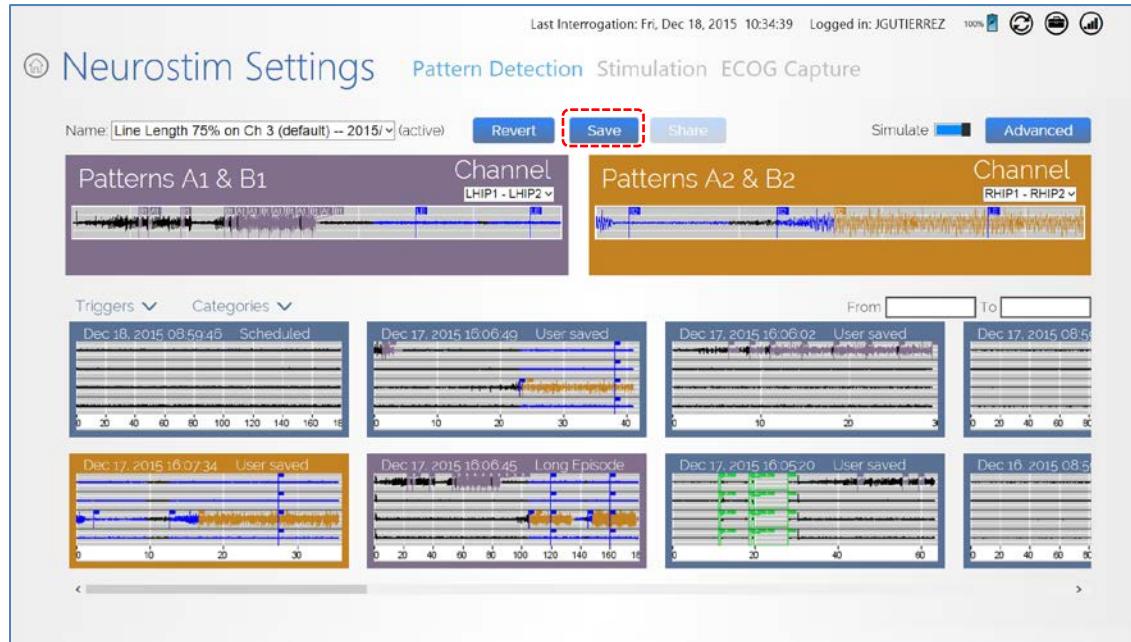


Figure 37 Detection settings ready to be saved

*Note: ECoG capture triggers are inactive during programming. **Programming clears activity information and ECoGs stored in the neurostimulator.** If there is any new information, be sure to interrogate the neurostimulator prior to programming.*

8. **Save and Program:** Back in the **NEUROSTIM SETTINGS** screen, select **SAVE** at upper center; name the detection set in the **SAVE DETECTION SET** dialog and select **SAVE**; then select **REVIEW & PROGRAM** at upper right and, while holding the wand within approximately 1 inch of the neurostimulator, follow the prompted steps to **CONFIRM PROGRAMMING** and **CONTINUE** sending the new settings to the neurostimulator.

Phase 3: Customize Pattern Detection

To customize pattern detection is to refine detection settings to promptly detect the specific patterns for this patient, and not detect on baseline activity, since detection serves as the basis for neurostimulator treatment. On visits subsequent to **Phase 2: Define a Pattern Detection Set**, the usual, systematic approach is to make small, incremental adjustments to individual parameters of a detection pattern. The adjustment you make may be intended to increase or decrease how early and how often the pattern is detected, based on your assessment of current detection activity as viewed on stored ECoGs and on the **Activity Screen** (see page 28).

Note: Typically, physicians repeat earlier phases (Phase 1 or Phase 2) only if they observe a new type of ECoG activity of interest that requires an additional detector; or if a good clinical response has not been obtained and they wish to start over.

1. Select **NEUROSTIM SETTINGS** from the **HOME** screen,
2. On the **NEUROSTIM SETTINGS** screen, select **PATTERN DETECTION** at upper center. The **PATTERN DETECTION** screen opens.

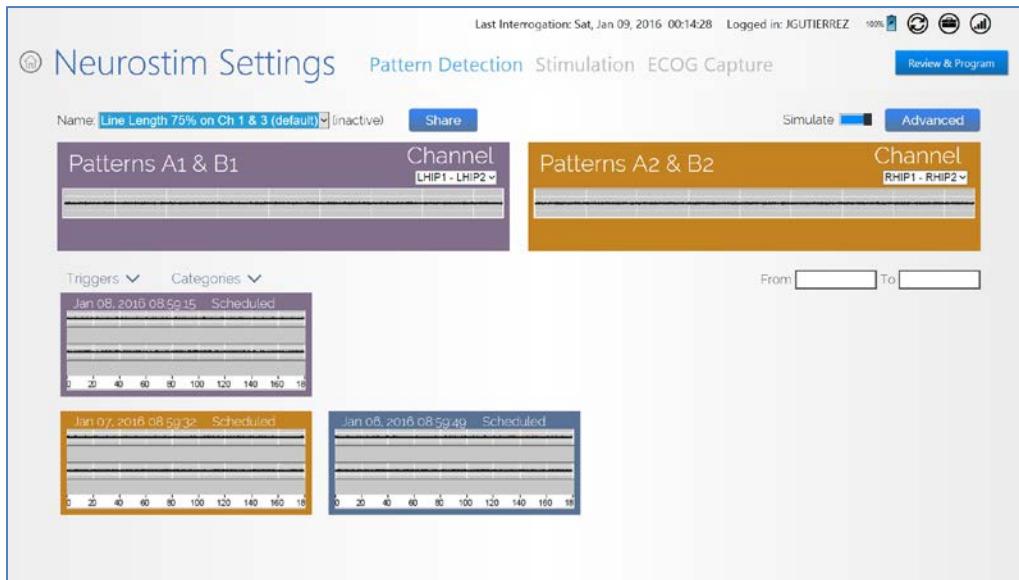


Figure 38 Pattern Detection screen

- The **NAME** field at upper left identifies the currently active detection set, as indicated by **(ACTIVE)** next to the field. You can edit this set or select an inactive set to program by selecting the **NAME** field.

*Note: To customize the current pattern detection set, do not change the active detection set shown by default in the **NAME** field at upper left.*

- The upper half of the screen provides tiles to select the current detection patterns and channels; the tiles are named **PATTERNS A1 & B1** (upper left) and **PATTERNS A2 & B2** (upper right).
 - The lower half of the screen shows thumbnails of all ECoGs (if no filters are on) stored in the PDMS that were captured using the current recording montage. You can filter the ECoGs shown using the **TRIGGERS** and **CATEGORIES** menus at left center.
 - The **SIMULATE** button at upper right (next to **ADVANCED**) is switched **ON** (blue) by default, which enables simulations of detections to be shown on all ECoGs shown on this screen.
3. Select one of the two **PATTERN** tiles at the top to access the **DEFINE PATTERNS** dialog.

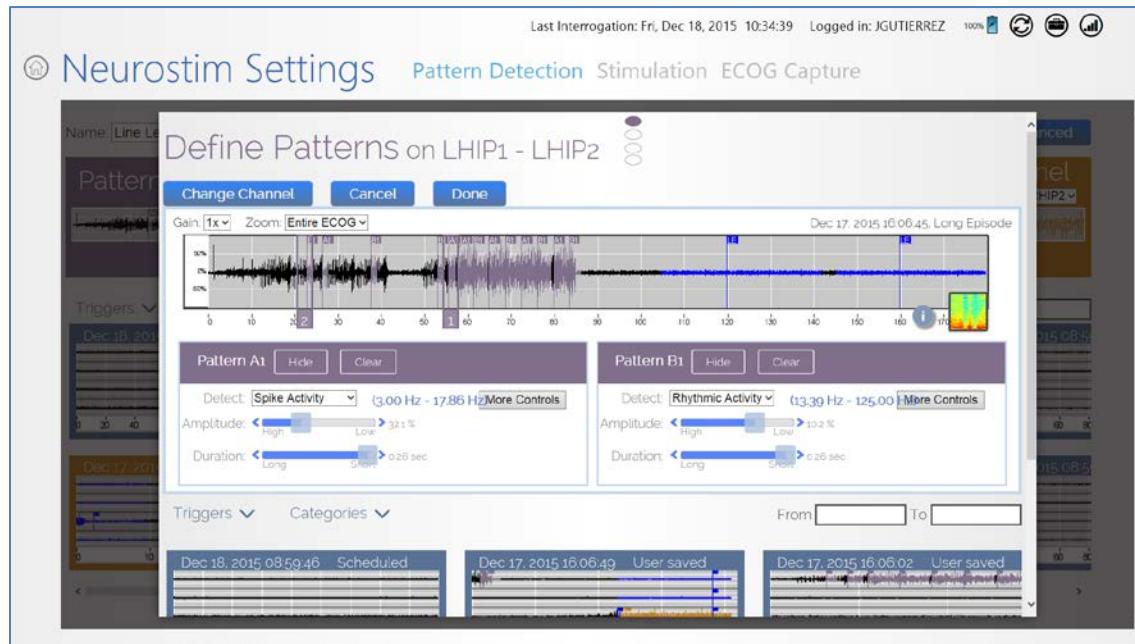


Figure 39 Define Patterns dialog

4. In the **DEFINE PATTERNS** dialog, select and adjust the desired parameter using the sliders so as to increase or decrease how early and how often the pattern is detected.

The software automatically simulates detection with the changed parameter as applied to this ECoG and to the patient's other stored ECoGs shown in smaller size below. It is recommended to review simulated detection activity on an ensemble of stored ECoGs before finalizing your adjustment.

- Swipe up to view the bottom of the dialog where the ECoGs appear with simulated detection; swipe left or right to see more ECoGs if available. You can filter the ECoGs using the **TRIGGERS** and **CATEGORIES** menus, and/or the **FROM** and **To** date range fields.

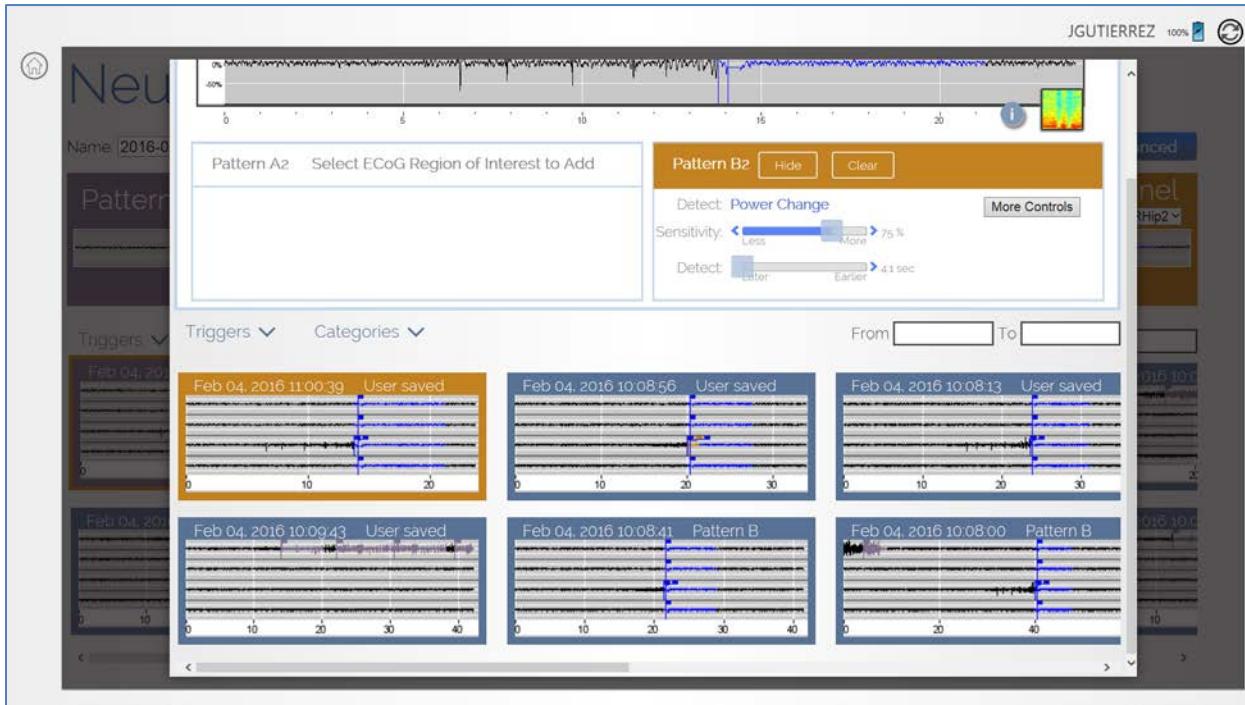


Figure 40 Detection simulation in all ECoGs

- To hide the simulation on the ECoG, select **HIDE**. To remove the pattern on the ECoG, select **CLEAR**.
- Detection parameters that affect sensitivity are key: adjustments that increase sensitivity tend to increase how early and how often the pattern is detected; adjustments that decrease sensitivity tend to decrease how early and how often the pattern is detected.
- To adjust sensitivity, the two most common parameters to adjust are the ones with sliders directly on the **DEFINE PATTERNS** dialog. For **SPIKE ACTIVITY** and **RHYTHMIC ACTIVITY**, the sliders adjust **AMPLITUDE (%)** and **DURATION** (seconds). For **POWER CHANGE**, the sliders adjust **SENSITIVITY (%)** and **DETECT** (seconds).

Note: The parameter values are shown to the right of each slider. For all four parameters, lower values result in greater sensitivity. Make adjustments and observe the effects on the simulation.

- The **MORE CONTROLS** button gives access to additional adjustment options and advanced settings.
- 5. Select **DONE** at upper right to close the **DEFINE PATTERNS** dialog. The detection set is not saved until you save it in the next step.
- 6. **Save and Program:** Back in the **NEUROSTIM SETTINGS** screen, select **SAVE** at upper center; rename the updated detection set if desired in the **SAVE DETECTION SET** dialog, and select **SAVE**. Then select **REVIEW & PROGRAM** at upper right and, while holding the wand within approximately 1 inch of the neurostimulator, follow the prompted steps to **CONFIRM PROGRAMMING** and **CONTINUE** sending the new settings to the neurostimulator.

- After changing detection settings, the **NEUROSTIM SETTINGS** screen has a **REVERT** button next to the **SAVE** button. Select **REVERT** if you want to remove all changes to the detection settings and revert to the current detection set.

DETECTION TOOLS AND DETECTORS

Note: Detection tools and detectors support you as you exercise clinical judgment to define detection settings. Detectors and settings initially suggested by the software (as described in Phase 2 on page 54) offer a starting point only. The software does not preempt nor replace your clinical judgment.

In previous versions of programmer software, the physician selected a set of ECoGs with specific patterns to be detected, chose a detection tool and settings, observed the simulated detections that would occur on that set of ECoGs, and then repeated the process iteratively to define a detection set.

In the current version, when you use the SimpleStart Feature to select a waveform pattern to detect, the software suggests a detector (one of **POWER CHANGE**, **RHYTHMIC ACTIVITY** or **SPIKE ACTIVITY**) and associated settings as a proposed fit for that pattern. You use this procedure initially (in **Phase 2: Define a Pattern Detection Set**, see page 54) to create a new detection set when you have first gathered ECoGs containing the type of activity you wish to detect. This detection set is a suggested starting point to assist you as you exercise clinical judgment to define detection settings. It suggests a starting point only. It does not preempt nor replace your clinical judgment. The software also simulates the suggested detector and settings on all of the patient's stored ECoGs. By reviewing the simulated detection on an ensemble of stored ECoGs, you can confirm that detection occurs when desired, or change the selected pattern, detectors and/or settings to make detection occur when desired. Having defined up to four detection patterns and associated settings, you then save them as a detection set. (Up to four detection patterns on up to two channels comprise a detection set.)

Detection Tools and Detectors

The RNS® Tablet uses three kinds of detection tools: **POWER CHANGE** (page 63), **BANDPASS** (page 67), and **AREA** (page 68).

- The **POWER CHANGE** detection tool (formerly line length) is used in the **POWER CHANGE** detector.
- The **BANDPASS** detection tool is used in two detectors, **RHYTHMIC ACTIVITY** and **SPIKE ACTIVITY**.
- The **AREA** detection tool is infrequently used and can be selected only manually from the **DEFINE PATTERNS** dialog (select **MORE CONTROLS > ADV CONTROLS**).

*Note: **RHYTHMIC ACTIVITY** and **SPIKE ACTIVITY** both are types of ECoG activity that can be detected by the **BANDPASS** tool. The **AREA** tool is not used to generate a suggested detector from an ROI, but can be manually programmed by selecting **MORE CONTROLS > ADV CONTROLS**.*

The table below summarizes this information.

Detection Tool	Detectors Used With
POWER CHANGE (FORMERLY LINE LENGTH)	Power Change
BANDPASS	Rhythmic Activity and Spike Activity
AREA	Area (configured manually only, not used as a suggested detector)

Power Change Detector

To understand how the power change detector works refer to **Figure 39**. Imagine the ECoG as a string that can be stretched so the length can be measured. ECoGs with higher amplitude or frequency oscillations will have longer length.

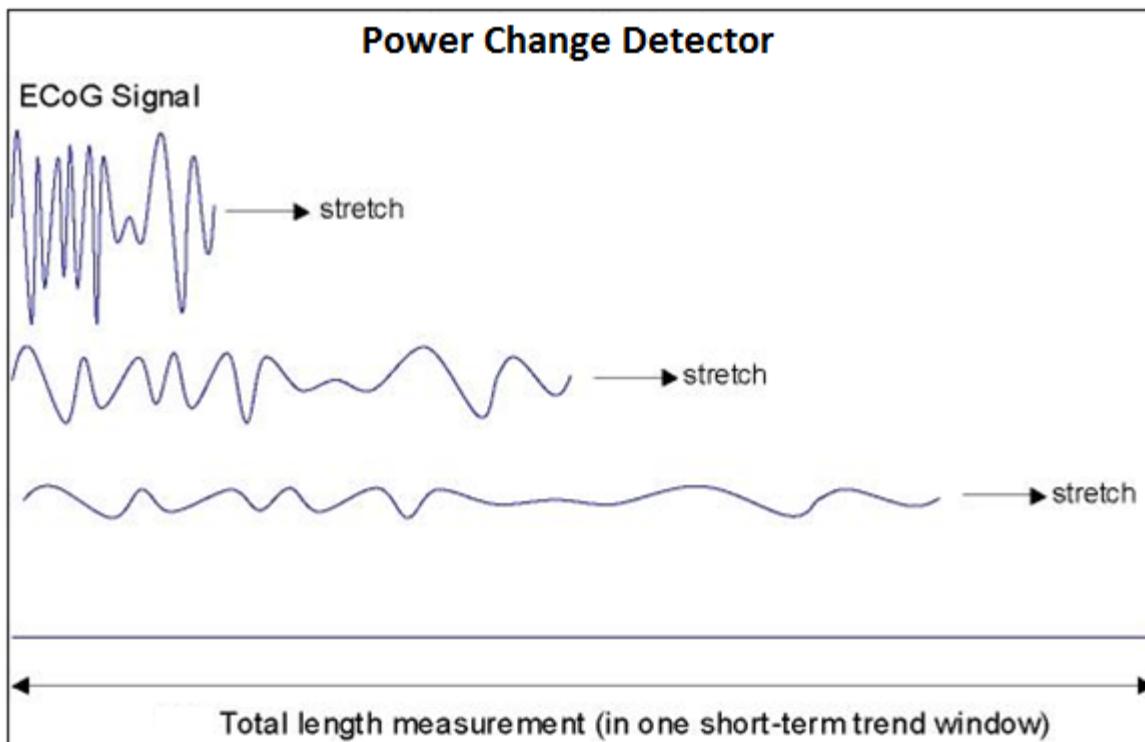


Figure 41 An illustration of the Power Change detector

The ECoG in the bottom of **Figure 39** below shows a developing seizure that corresponds in time with the short-term power change trend line (red) and the long-term power change trend line (blue). The power change detector works by comparing the normalized ECoG length (length per second) in a short-term trend window (typically 1 to 4 seconds, red line in graph) to the normalized ECoG length in a long-term trend window (typically 1 to 2 minutes, blue line in graph).



Figure 42 Power Change detector example

Detection occurs when the difference between the short-term and long-term trend windows exceeds a programmed **SENSITIVITY** percentage (typically 37.5% to 100%). Lower **SENSITIVITY** settings result in earlier and more frequent detections. The **DTECT** setting specifies the length of the short-term trend (in seconds). Shorter **DTECT** settings result in earlier and more frequent detections.

Rhythmic Activity Detector

The Rhythmic Activity detector provides highly specific and early detection of rhythmic activity within a defined bandwidth. It acts like a frequency filter and can therefore be used to detect activity within specific frequency bands (theta, alpha, beta, and gamma). To understand how the Rhythmic Activity detector works refer to **Figure 40**.

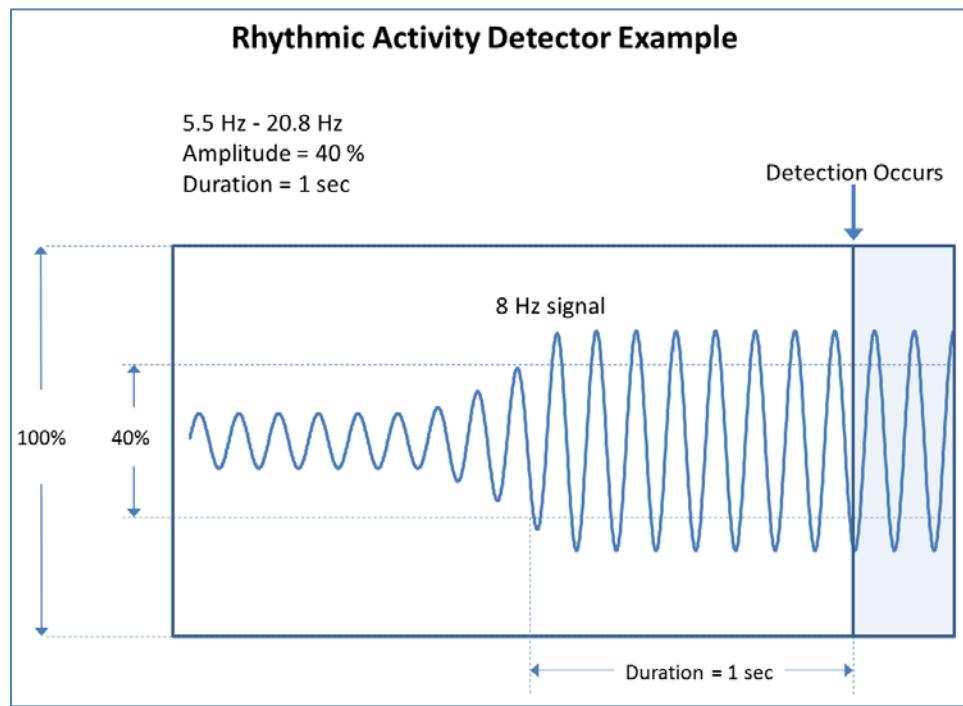


Figure 43 Rhythmic Activity detector example

The Rhythmic Activity detector analyzes the ECoG to determine when a signal is present within a specified bandwidth that meets specified amplitude and duration requirements. Four configurable parameters define the Rhythmic Activity detector operation: the **MINIMUM FREQUENCY** and **MAXIMUM FREQUENCY** that define the bandpass, the **AMPLITUDE**, and the **DURATION**. Detection occurs when a signal within the bandwidth exceeds the **AMPLITUDE** for the minimum **DURATION**. As **DURATION** is increased, the detector becomes less sensitive because it will take longer for a signal at a given frequency to be detected. As **DURATION** is decreased, the detector becomes more sensitive because it will take less time for a signal at a given frequency to be detected. As **AMPLITUDE** is increased, the detector becomes less sensitive because lower amplitude signals will be rejected. As **AMPLITUDE** is decreased, the detector becomes more sensitive because lower amplitude signals will be accepted.

Spike Activity Detector

The Spike Activity detector functions in a manner very similar to the Rhythmic Activity detector, but it is optimized for the faster transients that comprise spikes. To understand how the Spike Activity detector works refer to **Figure 41**.

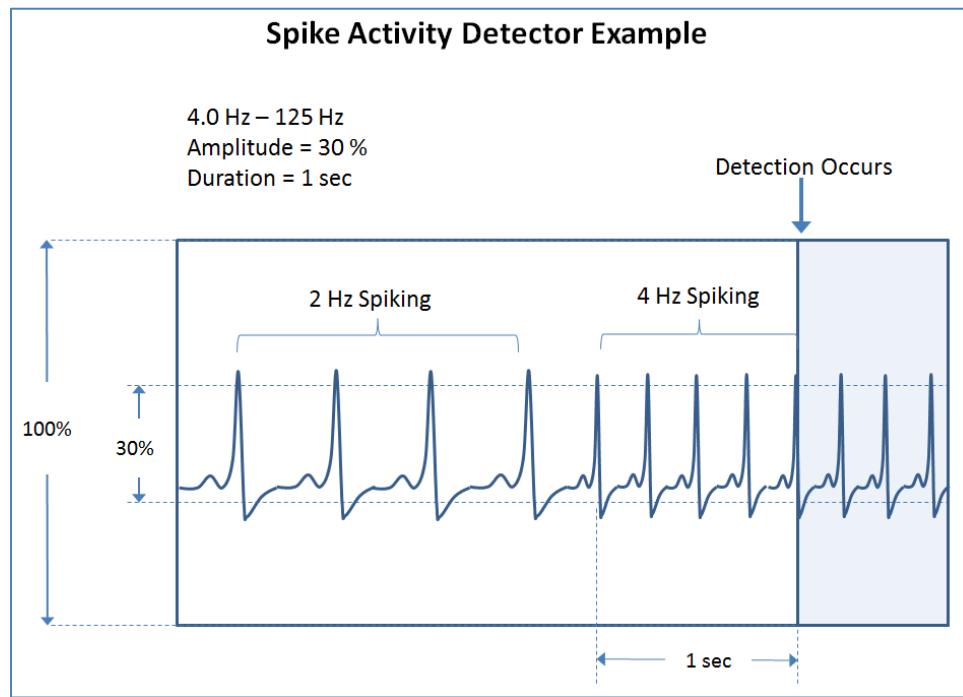


Figure 44 Spike Activity detector example

The Spike Activity detector analyzes the ECoG to determine when a signal is present within a specified bandwidth that meets specified amplitude and duration requirements. Four configurable parameters define the Spike Activity detector operation: the **MINIMUM FREQUENCY** and **MAXIMUM FREQUENCY**, the **AMPLITUDE**, and the **DURATION**. The **MINIMUM FREQUENCY** defines the slowest spiking rate that can be detected. The **MAXIMUM FREQUENCY** defines the upper passband for the detector. Because the signal frequency content of the fast transient portions of spikes are generally much higher than the spike rate, the **MAXIMUM FREQUENCY** usually needs to be set much higher than the spike rate. Detection will occur when spikes meeting the frequency requirements also exceed the **AMPLITUDE** for the **DURATION**. As **DURATION** is increased, the detector becomes less sensitive because it will take longer for detection to occur. As **DURATION** is decreased, the detector becomes more sensitive because it will take less time for a detection to occur.

Bandpass Detection Tool

The Bandpass detection tool provides highly specific and early detection of rhythmic and spiking activity. It acts like a frequency filter and can therefore be used to detect activity within specific frequency bands (for example, theta, alpha, beta, and gamma). The **RHYTHMIC ACTIVITY** and **SPIKE ACTIVITY** detectors use the Bandpass tool.

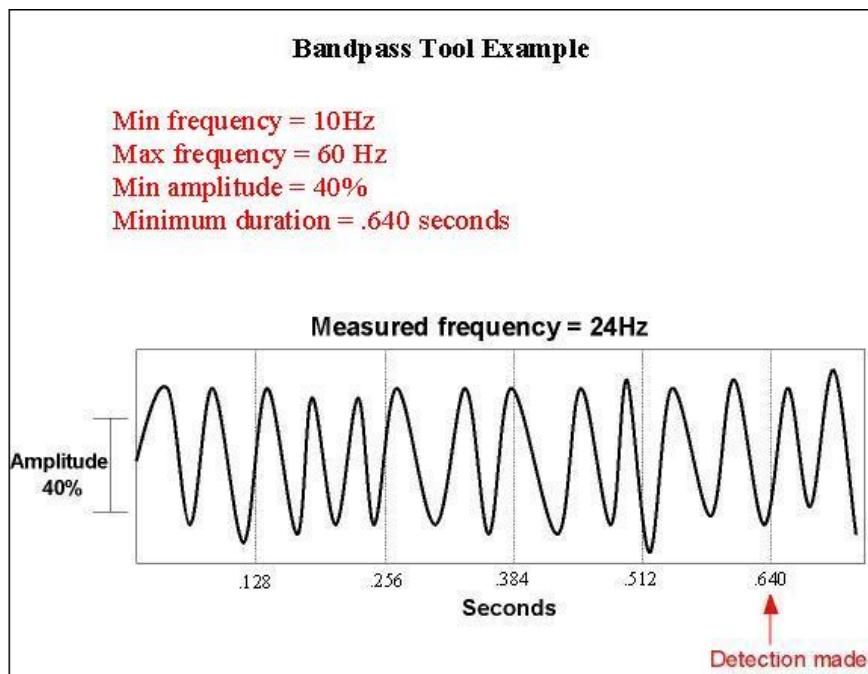


Figure 45 Bandpass detection tool example.

Area Detection Tool

The **AREA** detection tool can be selected only manually from the **DEFINE PATTERNS** dialog (select **MORE CONTROLS > ADV CONTROLS**). It is not used as a suggested detector by the software. The **AREA** detection tool measures the area between the ECoG and the time axis (see **Figure 43**).

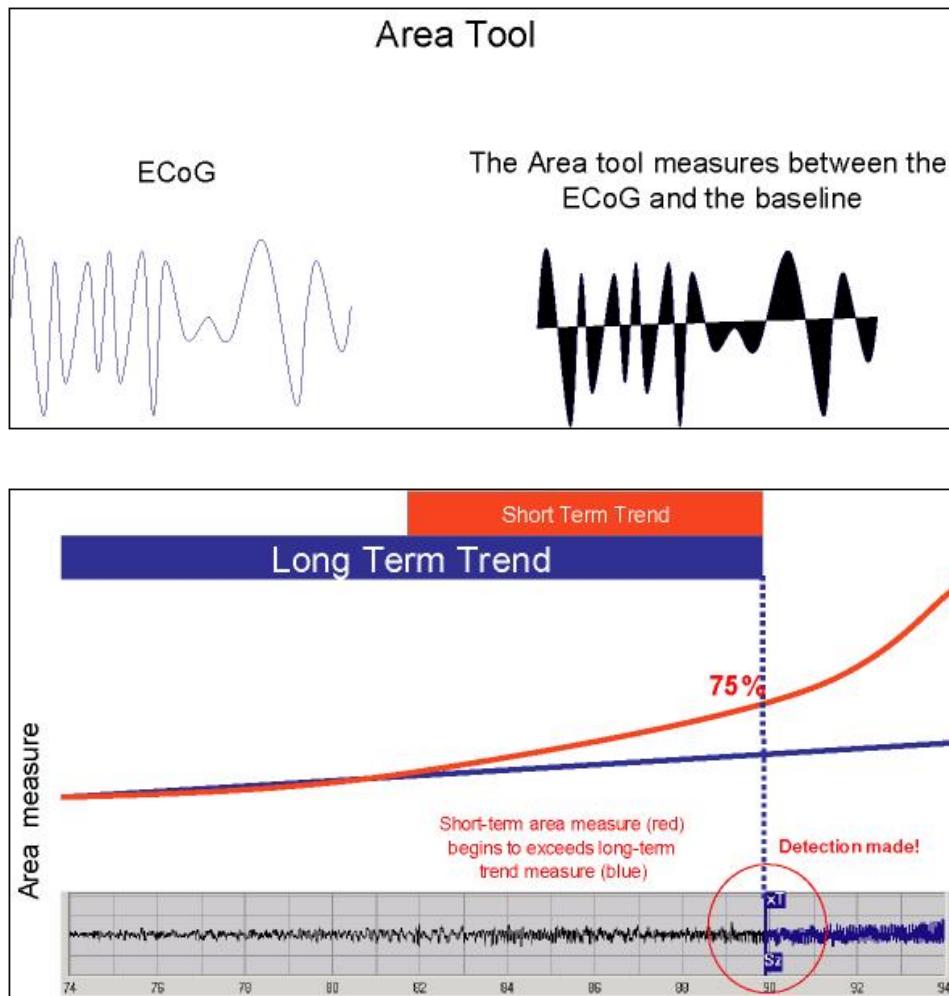


Figure 46 Area detection tool example

The Area detector operates by comparing the normalized ECoG area (area per second) in a short-term trend window (typically 1 to 4 seconds) to the normalized ECoG area in a long-term trend window (typically 1 to 2 minutes). Detection occurs when the difference between the short-term and long-term trend windows exceeds a programmed threshold (typically 37.5% to 100%).

Advanced Controls

The physician has access to advanced controls from the **DEFINE PATTERNS** dialog. Select **MORE CONTROLS > ADV CONTROLS** to open the **ADVANCED SETTINGS** dialog for the pattern you are currently defining, which gives access to additional settings for the detection tools **BANDPASS**, **POWER CHANGE** and **AREA**. (**AREA** can be selected only in this dialog.) Select the down-arrow next to **TECHNICAL PARAMETERS** to access even more settings. You can select more than one tool for use at the same time in this dialog.

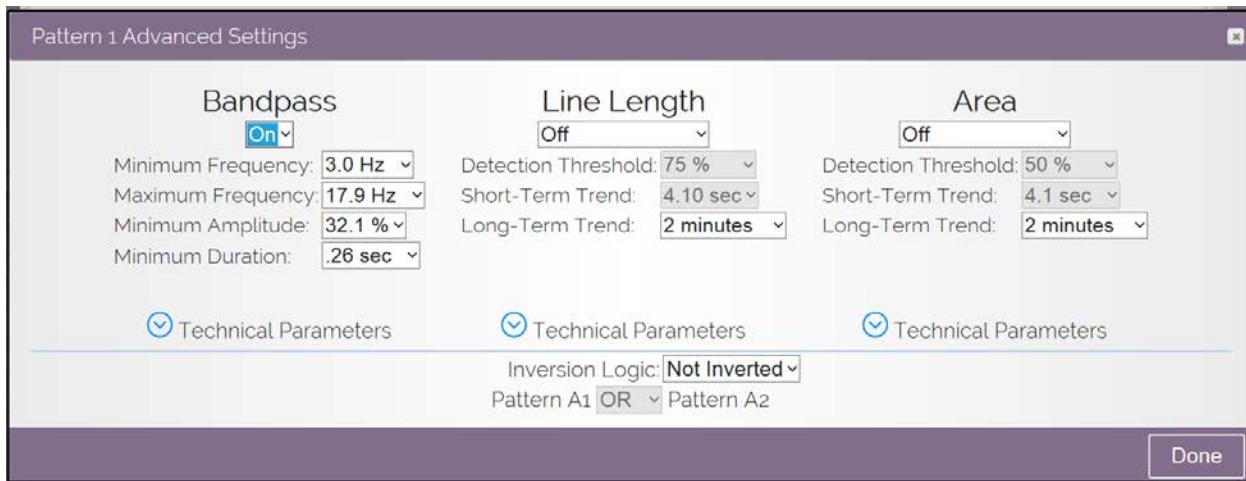


Figure 47 Bandpass, Power Change and Area settings, in Advanced Settings dialog

SHARE AND/OR DOWNLOAD A DETECTION SET

This section explains how to share a detection set with the users or programmers you select, and how to download a shared detection set. For example, use this feature when more than one physician is treating a patient.

You must be connected to the PDMS (via the tablet or a web browser on a personal computer) to define or modify a custom detection set, and by default, these are saved to the PDMS only—not locally. If you are connected to the PDMS during programming, you can program the neurostimulator with any detection set saved under the current username or saved to the programmer you are using. However, if you are not connected to the PDMS (offline) during programming, you can program only the pre-programmed default detection sets or those detection sets previously downloaded from the PDMS to the programmer.

Steps to Share a Detection Set

1. Whenever connected to the PDMS, either on a programmer or using a web browser on a personal computer, you can share the currently selected detection set by selecting the **SHARE** button on the **NEUROSTIM SETTINGS > PATTERN DETECTION** screen. The **SHARE DETECTION SET** dialog opens.

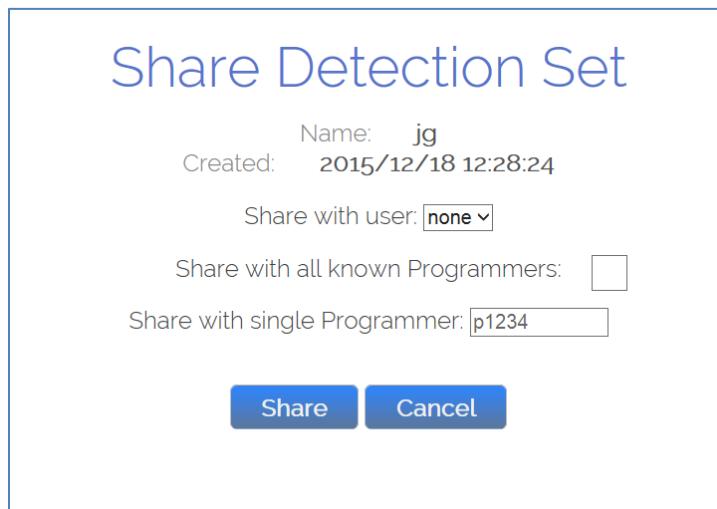


Figure 48 The **SHARE DETECTION SET** dialog

2. You can select to share with:
 - a. A specific user, who could then access it on any programmer to which he is logged in. Touch the **SHARE WITH USER** field to select from a list of available usernames, if any. (Available usernames would include those users who have been logged in to a programmer while interrogating the neurostimulator associated with the patient to which the detection set applies.)
 - b. All known programmers, which means any programmer that has interrogated the neurostimulator associated with the patient to which the detection set applies
 - c. One specific programmer, which you must identify by its serial number. The serial number is on the label affixed to the back of the programmer.
3. When you have selected a user or programmer(s) to share with, select the **SHARE** button in the **SHARE DETECTION SET** dialog. A message appears to tell you that the detection set was saved and will be sent to the indicated user and/or programmer(s) the next time they are synchronized with the PDMS.

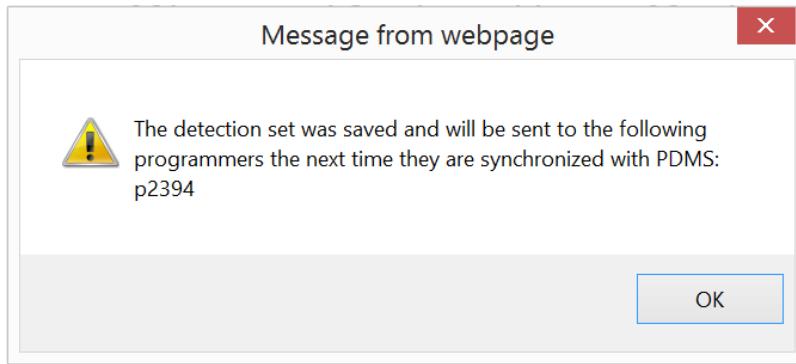


Figure 49 Example message reporting successful sharing

What Happens to a Shared Detection Set

Sharing a detection set using the foregoing method merely makes it available to the targeted user or programmer when connected to the PDMS. Shared detection sets are not automatically downloaded and saved locally to the programmer.

What happens is that the next time the targeted user or programmer connects to the PDMS, an icon with two wavy red lines will appear on the **NEUROSTIM INFO** tile. When they select the icon, a **NEUROSTIM INFO** dialog opens, where the alert icon is shown again with the text "**A pending Detection Set is available on PDMS. Click the icon to being loading it.**" When the user clicks this icon, an overlay appears that describes the detection set: the name, the time it was created, and the name of the person who shared it. They can **ACCEPT** the detection set, which will bring the user to the **PATTERN DETECTION** screen and automatically load the detection set, or they can **SKIP** the detection set, which clears the alert and allows the user to continue with their session. Whether accepted or skipped, the detection set will still be available in the list of detection sets in **NEUROSTIM SETTINGS**—as long as the programmer is connected to the PDMS. However, whether you accept or skip the detection set, the programmer does **not** download the detection set to the programmer for offline use. There is only one way to download a detection set to the programmer for offline use, described next.

Steps to Download a Detection Set to the Tablet

To download a detection set and save it locally to the tablet for offline use, follow these steps:

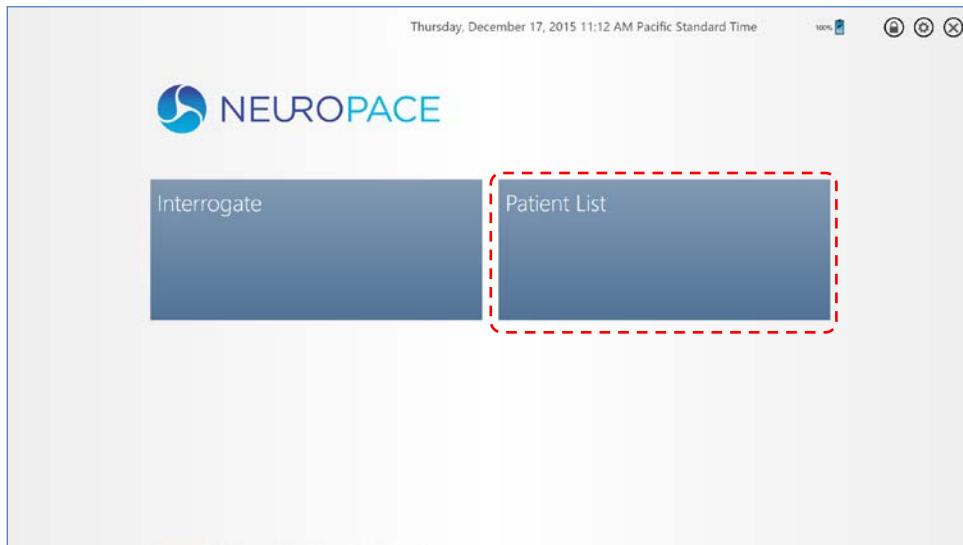


Figure 50 Start screen

1. Select the **PATIENT LIST** tile on the tablet Start screen. The patient list opens.
2. Select a patient from the patient list. The **PATIENT HOME** screen opens. Note this is different than the tablet **HOME** screen (shown in **Figure 6** on page 21).

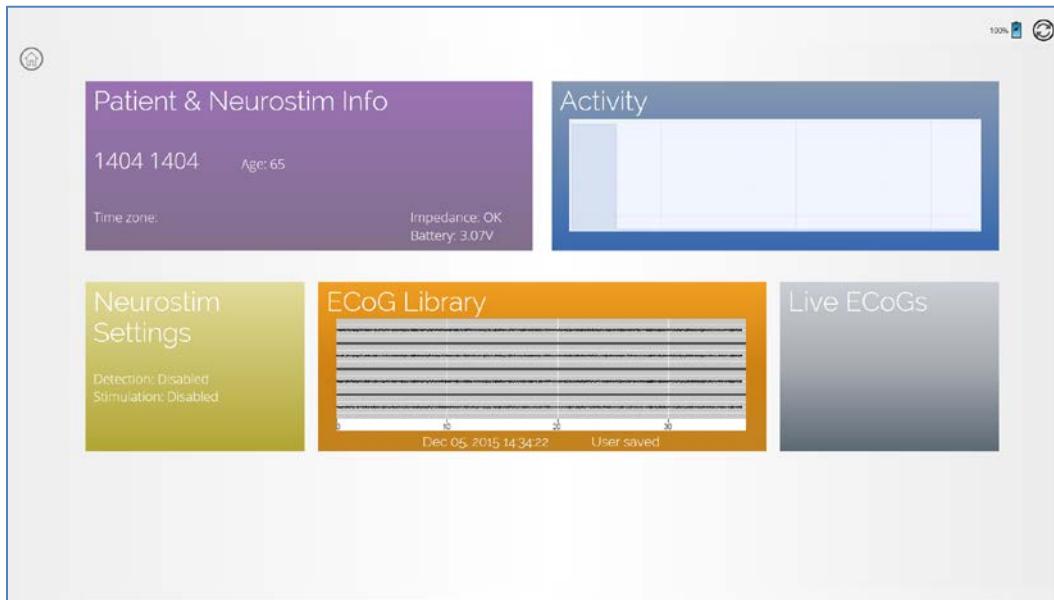


Figure 51 Patient Home screen

3. Select the Neurostim Settings tile. The **NEUROSTIM SETTINGS > PATTERN DETECTION** screen opens.

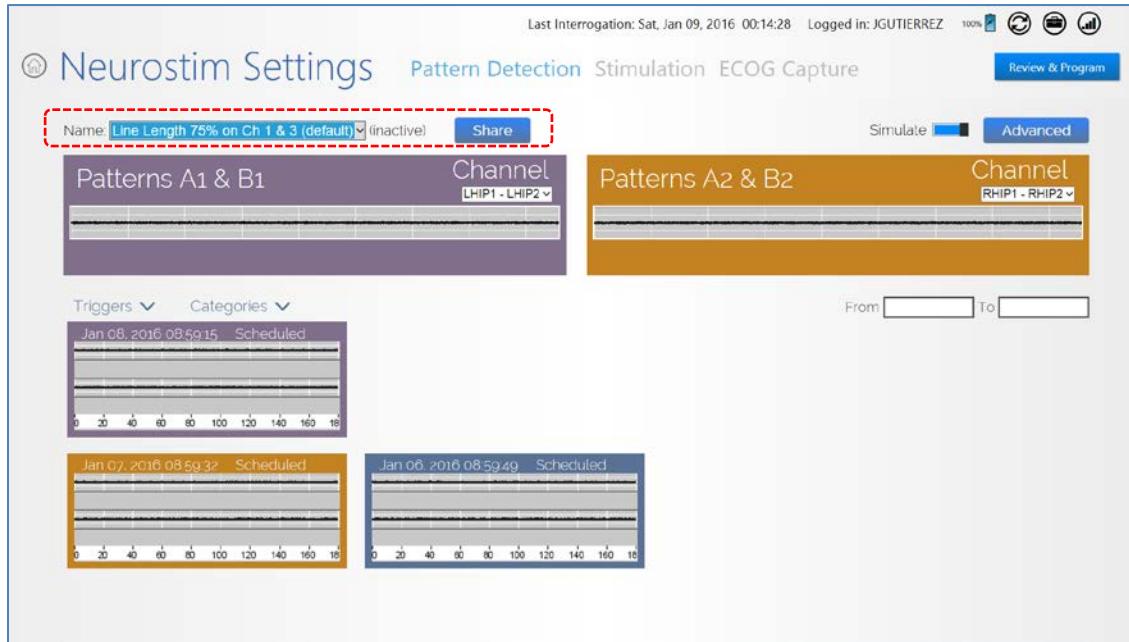


Figure 52 Pattern Detection screen

4. Select the **NAME** field and select the name of detection set you want to download to the tablet.
5. Select the **SHARE** button. The **SHARE DETECTION SET** dialog opens, but in this case, it includes a checkbox named “**Download to this Programmer now**” along with the other sharing options.

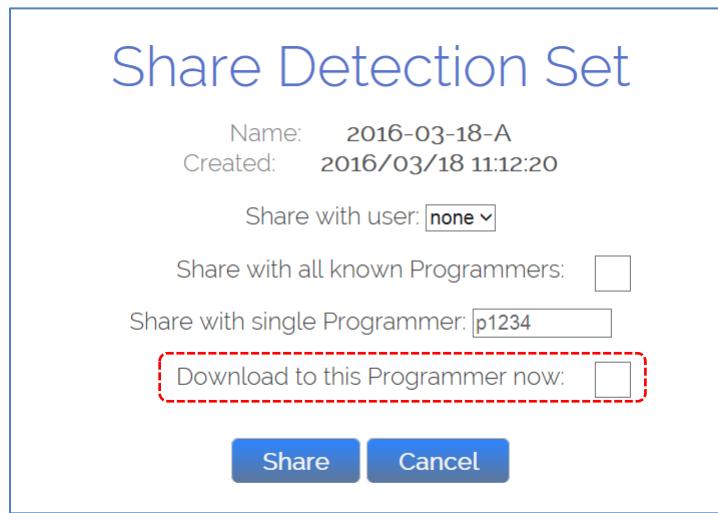


Figure 53 The SHARE DETECTION SET dialog with download now option

6. Select the “**Download to this Programmer now**” checkbox.

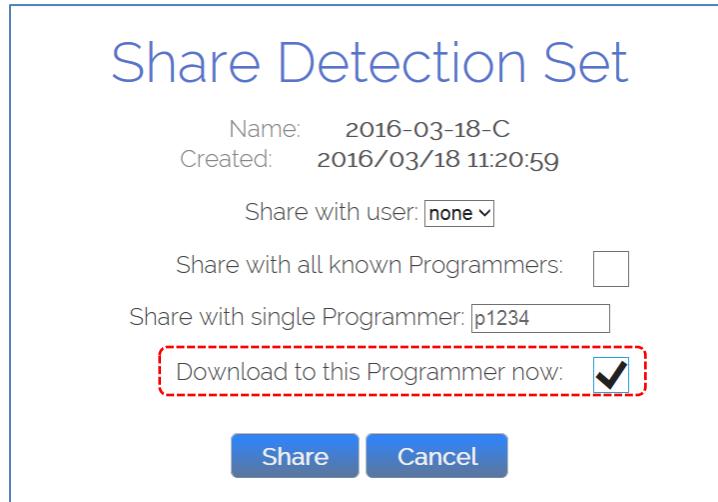


Figure 54 The SHARE DETECTION SET dialog with download now option selected

7. Select the **SHARE** button. A download progress dialog appears briefly, and when download is complete, a PDMS Transfer message reports “**The detection set download has completed successfully and will be available in offline mode.**”

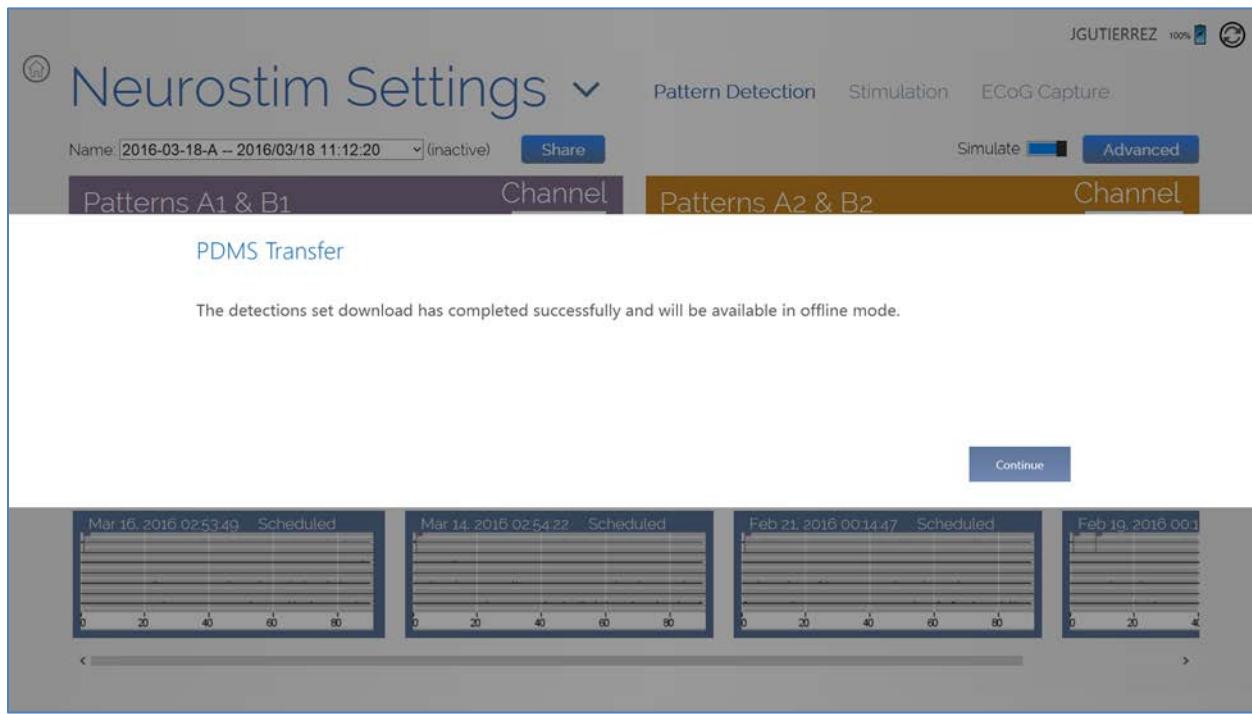


Figure 55 PDMS TRANSFER message reports a successful detection set download

OVERVIEW OF RESPONSIVE THERAPY

The RNS® Neurostimulator can deliver a therapy sequence of up to 5 individually configured therapies in response to detected activity. The neurostimulator resumes sensing after each therapy is delivered. If the activity is still detected, the next (sequential) therapy is delivered. If the activity is no longer detected, the remaining therapies are not delivered and the episode ends. Therapy sequence refreshes with the detection of each new episode.

The system enables you to configure behavior for the first therapy in therapy sequence in one of two ways:

- You can select different bursts to be delivered in response to a specific event detector (Pattern A or Pattern B), as depicted in **Figure 53**; or
- You can select up to two bursts to be delivered in response to any event detector, as depicted in **Figure 54**.

All subsequent therapies deliver a maximum of 2 bursts and are independent of the event detector.

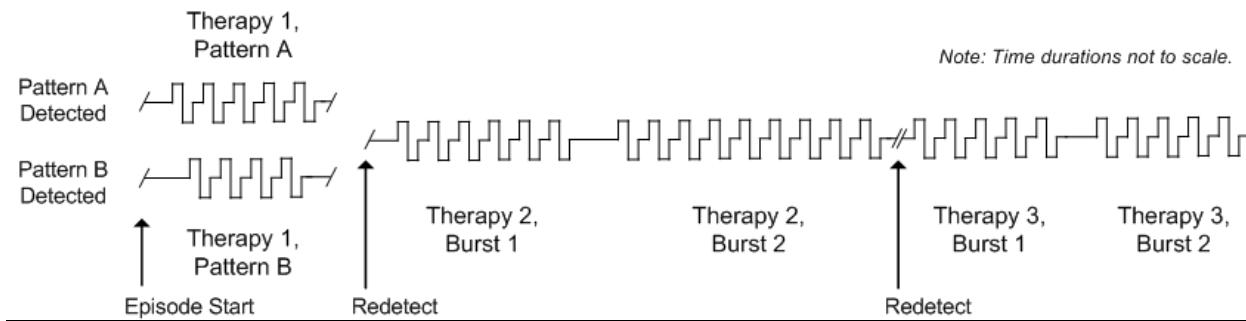


Figure 56 A THERAPY SEQUENCE DELIVERED IN RESPONSE TO A SPECIFIC EVENT DETECTOR

Therapy 1 is configured to deliver a single burst in response to a specific detector. In this example, Pattern A is detected at the episode start; Therapy Pattern A burst is delivered. If Pattern B is detected at the episode start, Therapy Pattern B burst is delivered. Subsequent therapies, Therapy 2 and Therapy 3, are configured to deliver 2 bursts each, no matter which pattern is subsequently detected.

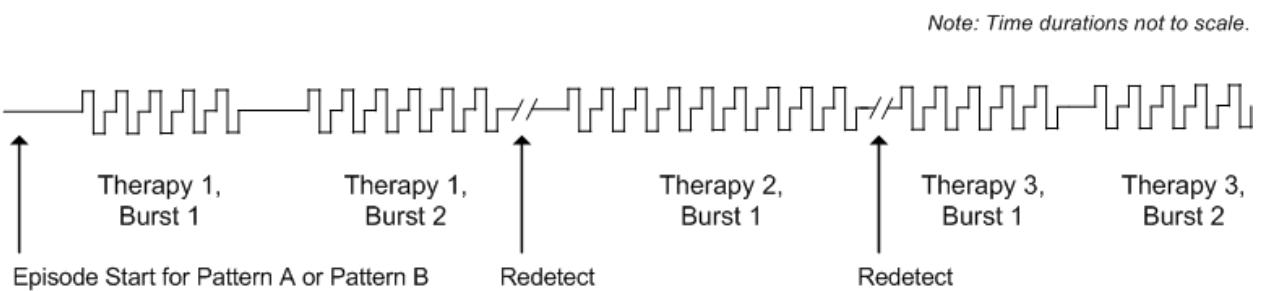


Figure 57 A THERAPY SEQUENCE DELIVERED IN RESPONSE TO ANY EVENT DETECTOR

Therapy 1 is configured to deliver Burst 1 and Burst 2 in response to any event detector. In subsequent therapies, Therapy 2 is configured to deliver Burst 1 only and Therapy 3 is configured to deliver Burst 1 and Burst 2.

The following figure identifies the parameters of a burst of stimulation therapy.

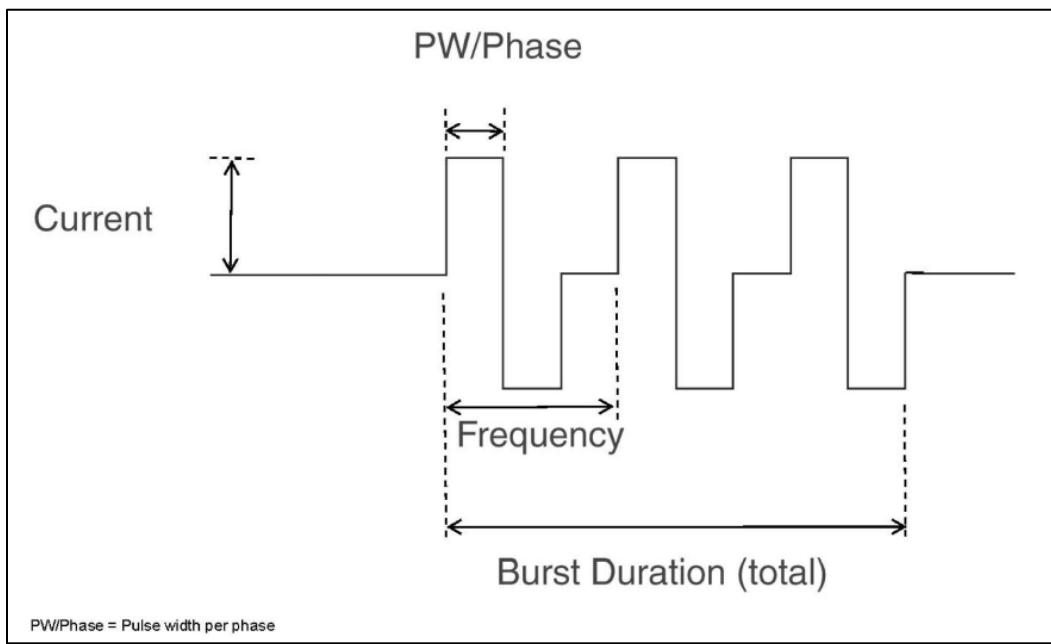



Figure 58 Parameters of a single burst of stimulation

CONFIGURE RESPONSIVE THERAPY

Detection must be enabled before responsive therapy can be enabled. Before programming responsive therapy, at least one detector must be defined in a detection set. (If you select therapy specific to a detection pattern, at least one detector must be defined for each pattern.) To set up a therapy sequence, determine the number of therapies desired for delivery in response to an episode. Select each therapy and configure each burst, which is a collection of one to many biphasic pulses.

Note: For all settings, changed settings are indicated by blue buttons, with the new setting above and the old setting below. Gray buttons indicate settings have not changed.

Note: It is recommended to review and perform detection analysis prior to enabling responsive therapy.

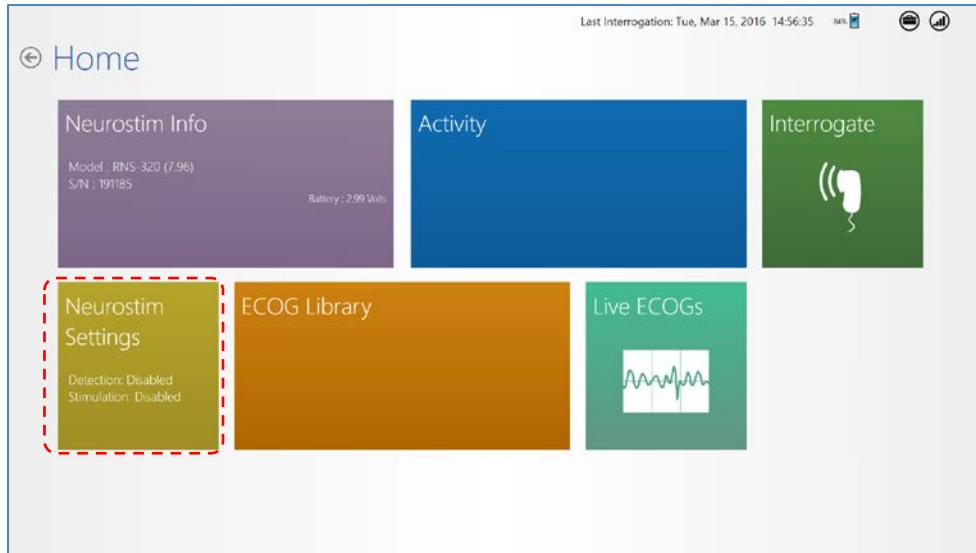

1. Select **NEUROSTIM SETTINGS** from the **HOME** screen.
2. Select **STIMULATION** at upper center of the **NEUROSTIM SETTINGS** screen. The **STIMULATION** screen opens. (Swipe left to see the right side.)

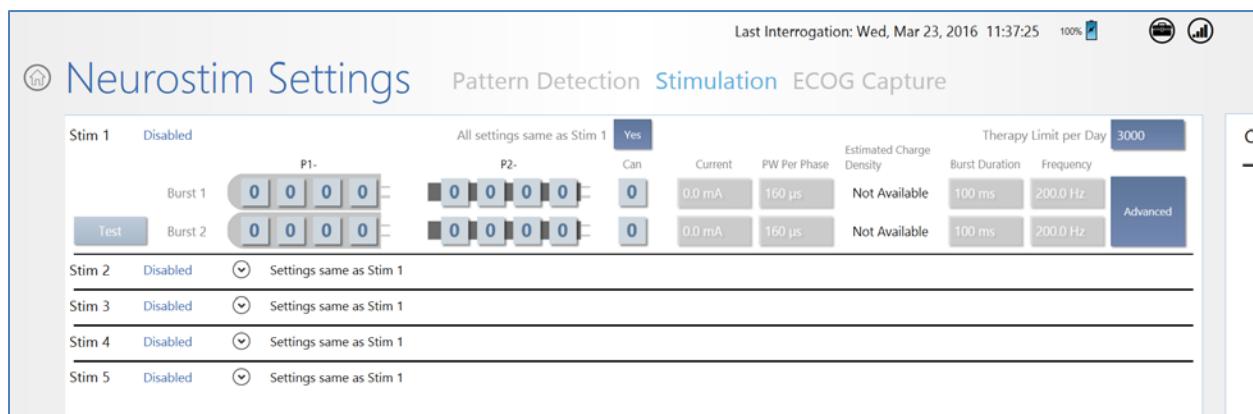


Figure 59 Stimulation screen, left side. Swipe left to see right side.

The screen has rows **STIM 1** through **STIM 5** to program up to five stimulation therapies. Each therapy is either configured or disabled. When one of the five therapies is disabled, all successive therapies are disabled also—you cannot skip a therapy. Recall that detection must be enabled before any therapy can be configured.

3. **Select Stimulation Pathway:** The stimulation pathway consists of the electrodes on each lead across which stimulation is delivered, as configured by the user. For each desired stimulation therapy and burst, starting with **STIM 1 BURST 1**, select the electrodes on each lead across which to deliver therapy, by making at least one electrode positive [+] and another negative [-]. Note that the neurostimulator can is an optional electrode.

In general, stimulation should be delivered to the leads and electrodes from which electrographic patterns of interest are observed. For example, if electrographic activity of interest is observed on all channels, then the stimulation pathway should be configured to

stimulate across all electrodes. However, if electrographic activity of interest is observed on only 2 channels, then the stimulation pathway should be configured such that current is delivered through only those electrodes with electrographic activity of interest.

4. Select the level of **CURRENT** to be delivered during pulses.
5. Select the **PW PER PHASE**, which is the pulse width or duration (in μs) of each phase of the biphasic current pulse.

*Note: The **ESTIMATED CHARGE DENSITY** reported on screen is a function of **CURRENT**, **PW PER PHASE**, and the number of electrodes in the stimulation pathway. The tablet limits the available choices for **CURRENT** and **PW PER PHASE** to those which yield an **ESTIMATED CHARGE DENSITY** less than 25 $\mu\text{C}/\text{sq cm}$.*

6. Select the **BURST DURATION**, from 10 ms to 5000 ms.
7. Select the **FREQUENCY** of the biphasic current pulses, from 1 Hz to 333.3 Hz. Increasing frequency decreases the time between each biphasic pulse of current, and vice versa.
8. If desired, apply the same settings to all five stimulation therapies. To do this, select the **YES** button next to **ALL SETTINGS SAME AS STIM 1**, near top center of the screen.

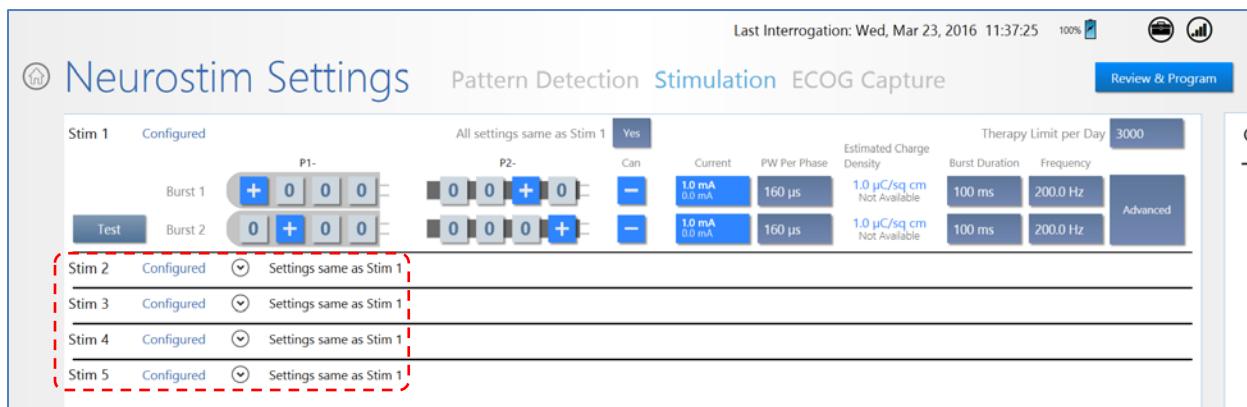


Figure 60 Stim 1 settings copied to all therapies

9. If desired, select **THERAPY LIMIT PER DAY** at upper right to adjust the maximum number of episodes the neurostimulator will treat per day.
10. If desired, swipe left to access the **OPTIONS** area on the right side of the **NEUROSTIM SETTINGS** screen to find the following options:



Figure 61 Options on the right side of the Stimulation settings screen

- a. **BEHAVIOR FOR THERAPY 1:** The default is **EITHER DETECTION PATTERN RESULTS IN SAME THERAPY**. Select the field to select the alternate **DETECTION PATTERN SPECIFIC THERAPY**, to apply only Therapy 1 in response to detection of one specific pattern you have programmed.
- b. **THERAPY LIMIT RESET TIME:** Select the time of day (in 24-hour time) when the daily therapy limit count resets to zero.
- c. **MAGNET THERAPY WITHHOLDING:** The default is **ENABLED**, which enables the magnet to withhold therapy as long as the magnet is over the neurostimulator. Once the magnet is removed, therapy can resume. To disable this functionality, select the field to toggle to **DISABLED**.

Skip to **Test Stimulation Before Enabling** on page **83** if you do not want to use the advanced settings.

Optional Advanced Settings

If desired, select **ADVANCED** at far right to open the **ADVANCED SETTINGS** screen for the selected therapy (**STIM** row). By default, all advanced settings are disabled.

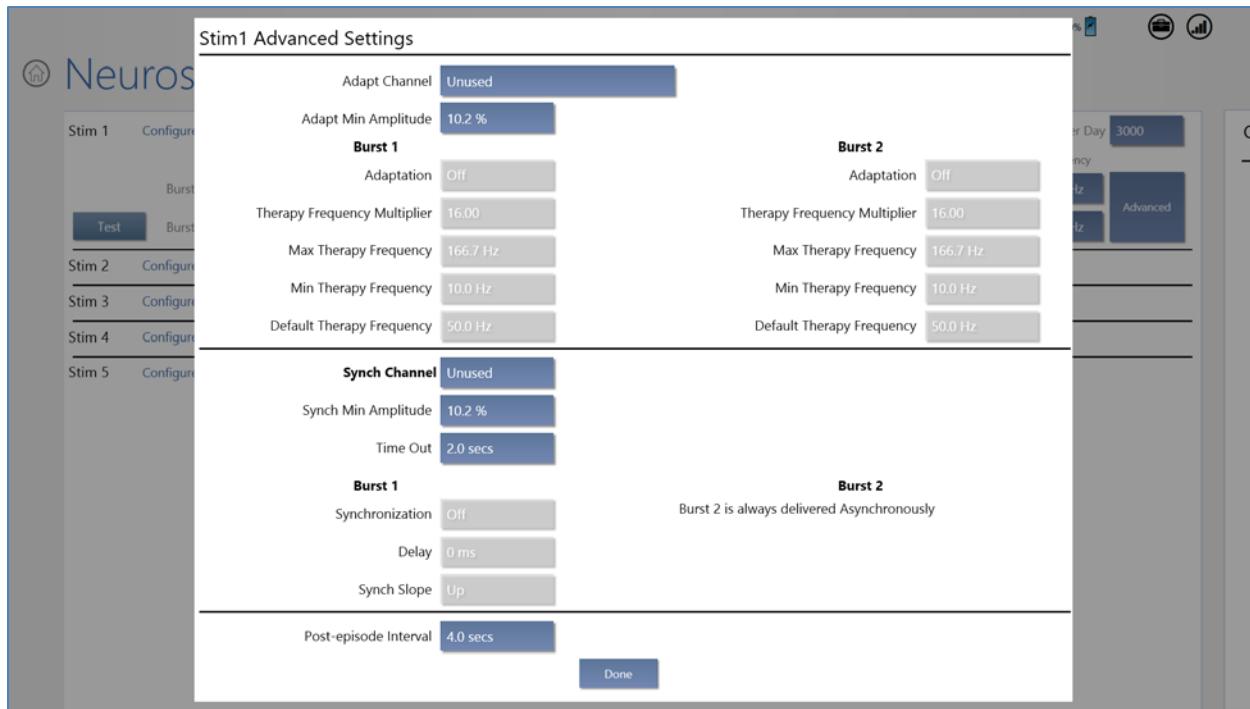


Figure 62 Stim 1 Advanced Settings

The **ADVANCED SETTINGS** screen has upper, lower and bottom portions that let you select the following:

Upper portion—Adaptive Therapy: Adapts the frequency and the minimum amplitude of the stimulation therapy as a percentage of the detected signal amplitude over a chosen channel, for one or both bursts. The available options are:

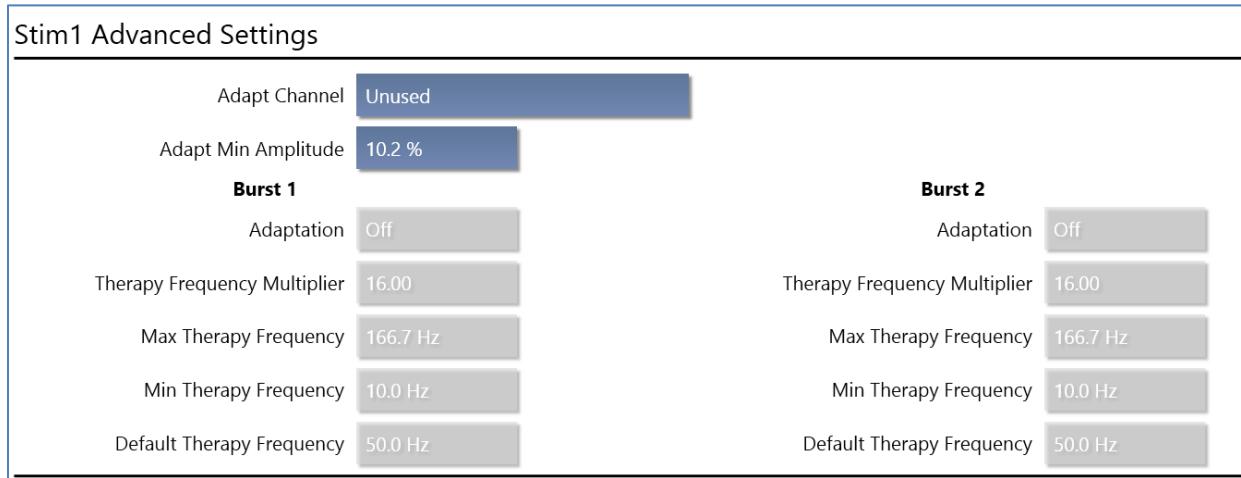


Figure 63 Adaptive therapy settings—upper portion of Advanced Settings dialog

- **ADAPT CHANNEL:** Select one of four channels over which to enable adaptive therapy. By default, this is **UNUSED**, which means adaptive therapy is disabled.
- **ADAPT MIN AMPLITUDE:** Select a percentage of the detected signal amplitude to be used as the minimum therapy amplitude.

Then for **BURST 1** and/or **BURST 2**, select settings for the following:

- **ADAPTATION:** Turn **ON** or **OFF** adaptive therapy for this burst.
- **THERAPY FREQUENCY MULTIPLIER:** A greater multiplier increases frequency of the burst.
- **MAX THERAPY FREQUENCY:** Sets maximum frequency for this burst.
- **MIN THERAPY FREQUENCY:** Sets minimum frequency for this burst.
- **DEFAULT THERAPY FREQUENCY:** Sets default frequency for this burst.

Illustration: If **ADAPT MIN AMPLITUDE** = 100% and the detected signal frequency = 20 Hz (50 ms period), the stimulation frequency is delivered at 20 Hz (50 ms p-p interval). If **ADAPT MIN AMPLITUDE** = 50% and the sensed signal frequency = 20 Hz (50 ms period), the stimulation frequency is delivered at 40 Hz (25 ms p-p interval).

If the neurostimulator cannot detect a stable frequency in the ECoG signal, or if the calculated stimulation frequency is not between the **MAX THERAPY FREQUENCY** and **MIN THERAPY FREQUENCY**, then the burst is delivered at the **DEFAULT THERAPY FREQUENCY**.

Lower portion—Synchronized Therapy: For **BURST 1** only, synchronize delivery of therapy into a specific part of the detected ECoG waveform over a chosen channel. The available options are:

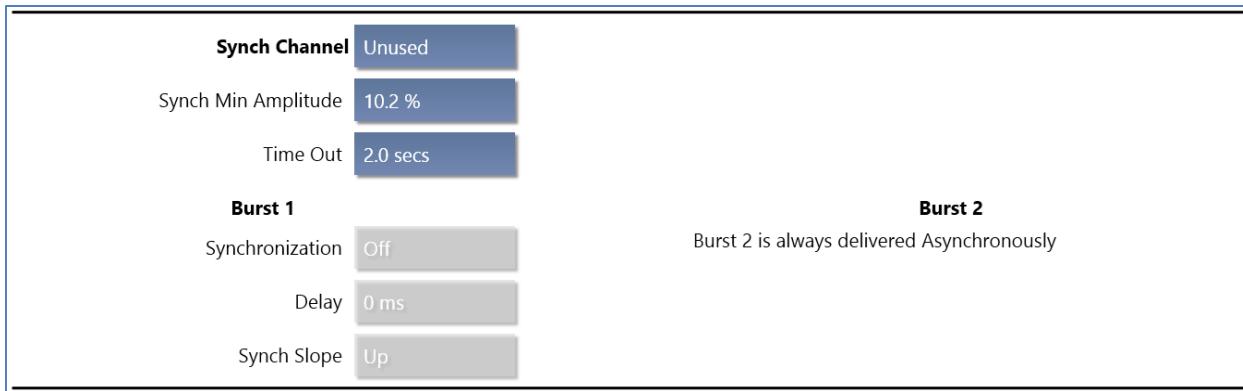


Figure 64 Synchronized therapy settings—lower portion of Advanced Settings dialog

- **SYNCH CHANNEL:** Select one of four channels over which to deliver synchronized therapy.
- **SYNCH MIN AMPLITUDE:** Select a percentage of the detected signal amplitude to be used as the minimum amplitude for synchronized therapy.

Then for **BURST 1** and/or **BURST 2**, select settings for the following:

- **SYNCHRONIZATION:** Turn **ON** or **OFF** synchronized therapy for **BURST 1**. (Burst 2 is always delivered asynchronously.)
- **DELAY:** Set the delay (in ms) for delivery of **BURST 1** therapy after detection of the target ECoG waveform.
- **SYNCH SLOPE:** Select **UP** or **DOWN**.

Note: If the neurostimulator is unable to synchronize therapy into a detected waveform, it reverts to asynchronous therapy delivery after 2 seconds.

Bottom portion—Post-Episode Interval: Select the **POST-EPIISODE INTERVAL**, from **OFF** (default) to **24 HRS**, to specify the time the neurostimulator waits after the end of an episode before treating a new episode. During the post-episode interval, monitoring continues but the system will not deliver therapy for a new episode.



Figure 65 Post-episode interval settings—bottom portion of Advanced Settings dialog

Save Advanced Settings: Select **DONE** to save the new advanced settings and return to the **STIMULATION** screen.

Test Stimulation Before Enabling

Before enabling responsive therapy, it is important to test stimulation to make sure it is well tolerated by the patient and there are no undesired changes in the ECoG, such as afterdischarges.

Note: The system cannot test stimulation during a detected episode or during therapy delivery. To test stimulation, either wait for the episode or therapy to complete, or disable detection and therapy.

1. Select **TEST** on the far left of the row to open the **TEST STIMULATION** screen and test the new stimulation settings. Note that the **TEST** button is not available if adaptive therapy is on (see **Optional Advanced Settings** above). An ECoG is running when the screen opens. Select **DELIVER** to deliver the test burst. Select **STOP/STORE** when desired to stop recording.

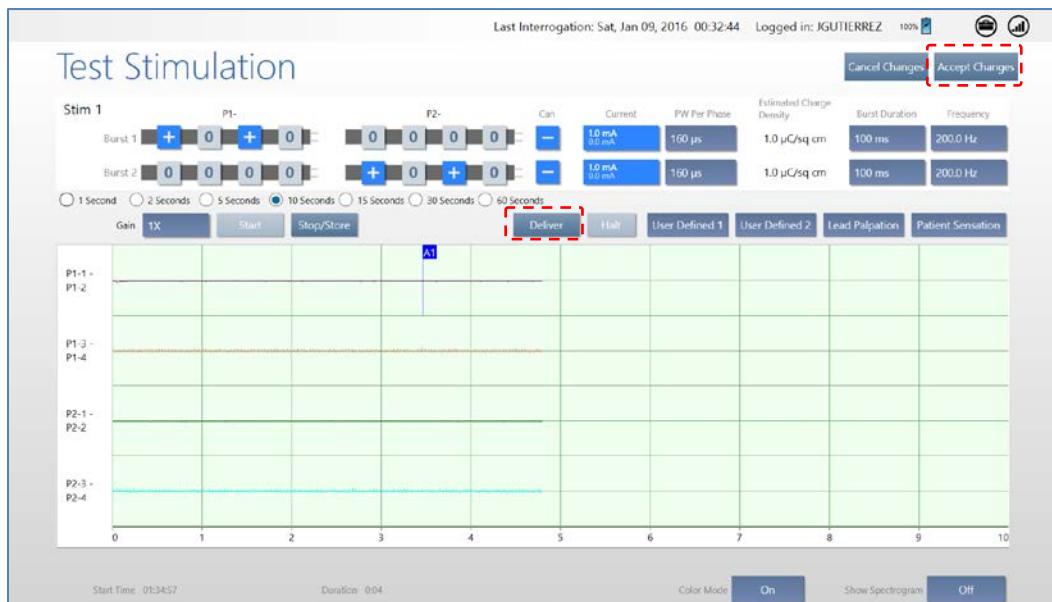


Figure 66 Test Stimulation screen

- You can adjust settings and retest while this screen is open. To exit test stimulation, select **ACCEPT CHANGES** or **CANCEL CHANGES** at upper right.
- Stimulation may evoke visual, olfactory, auditory, taste, sensory or motor responses.

- Test stimulations should be delivered to each electrode individually (electrode to neurostimulator) to determine whether there are acute side effects from stimulation.
- Habituation to side effects resulting from stimulation has been observed; however it is not advised to program settings that cause patient discomfort.
- Avoid stimulation settings that cause an electrographic afterdischarge.

Caution: Afterdischarge Activity

If evidence of afterdischarge activity resulting from stimulation is seen either on stored ECoGs or during test stimulation delivery, stimulation parameters should be adjusted to prevent such occurrence.

Enable Therapy: Review & Program

The new settings are not sent to the neurostimulator until you complete these steps.

1. Select **REVIEW & PROGRAM** at upper right. The **REVIEW** screen opens, showing a list of all changed settings. If you are satisfied with the new or changed settings, proceed to the next step. If not satisfied, select **CANCEL** at upper right to return to the previous screen.



Figure 67 Review screen after configuring therapy. Blue items have changed.

2. To activate the configured therapy, make sure **STIMULATION** is set to **ENABLED** at upper center.
3. On the **REVIEW** screen, if you are satisfied with the new or changed settings, place the wand within approximately 1 inch of the neurostimulator and select **CONFIRM PROGRAMMING**. If not satisfied, select **CANCEL** to return to the previous screen.

The **PROGRAMMING NEUROSTIMULATOR** dialog shows programming progress until complete.



Figure 68 Programming neurostimulator dialog shows progress of programming

Recommended Initial Responsive Therapy Settings

Responsive stimulation therapy should be enabled once the physician has determined the electrographic activity of interest is being detected.

- Program **BURST 1** for electrodes on one lead and **BURST 2** for electrodes on the other lead.
- **FREQUENCY** = 200 Hz (default)
- **BURST DURATION** = 100 ms (default)
- **PW PER PHASE** = 160 μ s (default)
- **CURRENT** = Adjust current to achieve Charge Density of 0.5 μ C/sq cm.
 - The RNS Tablet estimates charge density based on current, PW per phase, and the number of electrodes in the stimulation pathway. After selecting PW per phase and stimulation pathway, the physician should adjust current to achieve targeted charge density. For example, for a stimulation pathway with all electrodes on one lead programmed cathodal and the Can programmed anodal, using the default PW per phase value of 160 μ s and setting current to 1.0 mA results in charge density of 0.5 μ C/sq cm.
- Apply the same settings to all five stimulation therapies: select the **YES** button next to **ALL SETTINGS SAME AS STIM 1**, near top center of the screen.

In general, stimulation should be delivered to the leads and electrodes from which electrographic patterns of interest are observed. For example, if electrographic activity of interest is observed on all channels, then the stimulation pathway should be configured to stimulate across all electrodes. However, if electrographic activity of interest is observed on only 2 channels, then the stimulation pathway should be configured such that current is delivered through only those electrodes with electrographic activity of interest.

Recommended Therapy Modifications After Observing Clinical Response

Once the neurostimulator is programmed to provide responsive stimulation, the patient should be seen at approximately 3-month intervals to determine whether detection and/or responsive therapy settings should be modified. At each patient visit, the physician reviews ECoGs that are stored on the PDMS or programmer. Detection settings should promptly detect ECoG activity and not detect on baseline activity.

Based on the clinical response, the physician may choose to modify stimulation. After selecting PW per phase and stimulation pathway, current can be used to adjust charge density. It is recommended that the charge density be increased in 0.5 µC/sq cm increments. Longer pulse widths or fewer electrodes result in larger charge density changes for each mA current change. It is recommended that the new settings be tested (see **Test Stimulation Before Enabling** on page 83) to ensure that it is well tolerated by the patient and there are no undesired changes in the ECoG such as afterdischarges.

- Review data such as lead impedance, battery voltage, number of detections and stimulations, and stored ECoGs.
- Adjust stimulation settings if the patient reports an acute worsening or increase in clinical seizures that is not typical.
- Consider making minor adjustments if the patient reports only a mild or no improvement in clinical seizures.
- Allow sufficient time between stimulation adjustments to assess the patient's clinical seizure response.

In the RNS® System Pivotal study very few subjects used any of the additional responsive stimulation therapy options (see **Table 1**). These were: pattern specific therapy (each detector triggers a different stimulation setting); adaptive therapy (the stimulation frequency adjusts with the ECoG frequency); synchronized stimulation (stimulation is delivered into a specific part of the ECoG waveform); and post episode interval (responsive therapies are disabled for a specified period of time after detecting the end of the episode). **Table 1** provides the number and percentage of subjects who were treated with any of the additional stimulation therapy options.

Table 1: Additional responsive stimulation therapy options

	Number and (%) of Subjects Programmed
Pattern Specific Therapy	29 (15%)
Adaptive Therapy	4 (2%)
Synchronized Stimulation	3 (2%)
Post-Episode Interval	4 (2%)

EXPORT AND IMPORT DEVICE SETTINGS FOR NEUROSTIMULATOR REPLACEMENT

The tablet supports the transfer of all device settings from an old to a new neurostimulator as part of a neurostimulator replacement procedure. With the **EXPORT DEVICE SETTINGS** and **IMPORT DEVICE** functions, you can transfer all current neurostimulator settings (recording montage, detection settings, stimulation settings, etc.) to the new neurostimulator.

Export and import are done directly through the tablet without using the PDMS—exported settings are kept locally on the tablet—so no Internet connection is required. This also means that the programmer to which the settings have been exported must be physically present with the new neurostimulator to import settings to the new neurostimulator. This is not usually an issue because export and import are done on the same day: before, during or after the procedure to replace the neurostimulator.

Transfer Compatibility

You can transfer settings from an RNS-300M to an RNS-300M or to an RNS-320, and you can transfer settings between RNS-320 devices.

Note: *SCHEDULED IMPEDANCE* is a new feature in the RNS-320. It allows you to select which electrodes to test for impedance. Since the feature is not present in the RNS-300M, you must select the electrodes to test (and confirm programming) on the new RNS-320 before scheduled impedance works. Its settings are created as part of the recording montage—see **Create the Recording Montage** on page 38 for instructions.

Export Neurostimulator Settings

1. Interrogate the old neurostimulator, the one **from** which you will transfer settings.
2. Select the **UTILITIES** (toolbox) icon at upper right of the screen. The **UTILITIES** dialog opens.



Figure 69 Utilities dialog

3. Select **EXPORT DEVICE SETTINGS**. You are done. All current neurostimulator settings are stored on the tablet ready for import.

Note: The tablet retains only the last device settings exported to it—from any neurostimulator. Every export operation from any neurostimulator overwrites these device settings, if any. Successfully importing and programming the settings from the tablet to a new neurostimulator also deletes these device settings from the tablet.

Import Neurostimulator Settings

1. Interrogate the new neurostimulator, the one **to** which you will transfer settings.
2. Select the **UTILITIES** (toolbox) icon at upper right. The **UTILITIES** dialog opens (see figure above).
3. Select **IMPORT DEVICE**. The **IMPORT DEVICE SETTINGS** dialog appears. It shows the device settings available for import by device serial number and the date and time of last interrogation before the settings were exported to the tablet. It also notifies you if the tablet has no settings that came from a compatible neurostimulator (for example, if you exported from an RNS-320 for use on an RNS-300M).

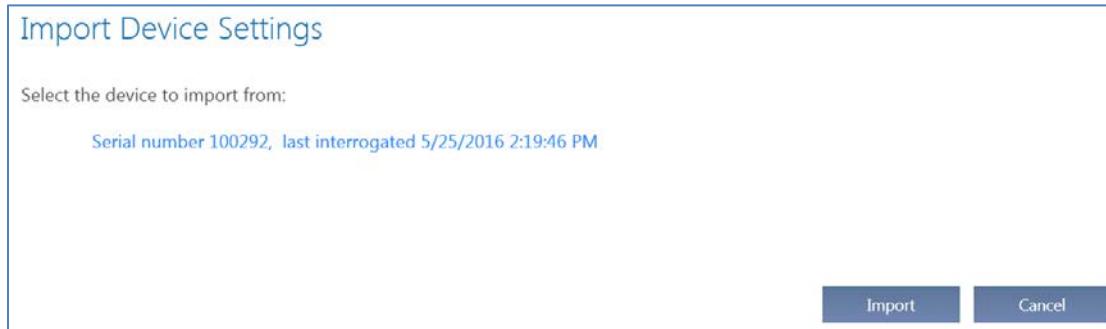


Figure 70 Import Device Settings dialog

Note: Device settings can be imported only once because they are deleted from the tablet once you successfully import and select **CONFIRM PROGRAMMING** from the **REVIEW** screen, as described below.

4. Select the device settings—the selected settings will be highlighted—and then select **IMPORT**. The dialog notifies you that it is importing device settings.

When complete, the tablet Home screen opens. It will include a red **REVIEW & PROGRAM** tile at lower right.


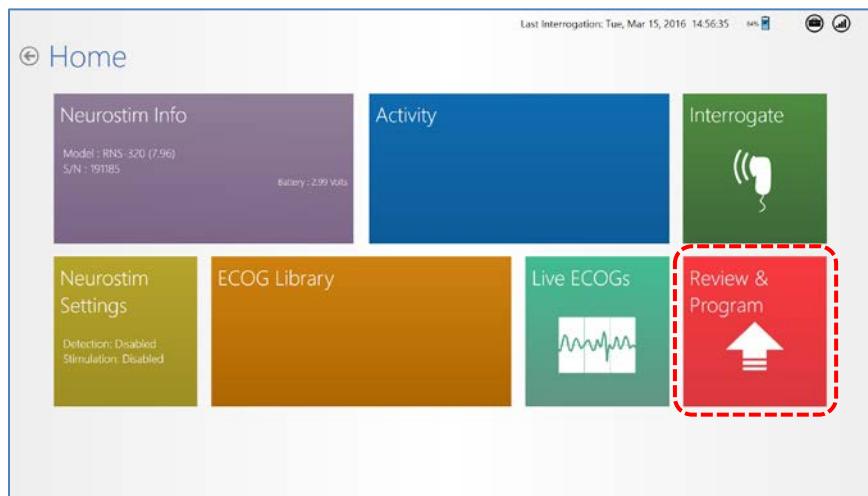


Figure 71 Home screen with the red Review & Program tile at lower right

5. Select **REVIEW & PROGRAM**. The **REVIEW** screen opens, showing all the settings that will be programmed.

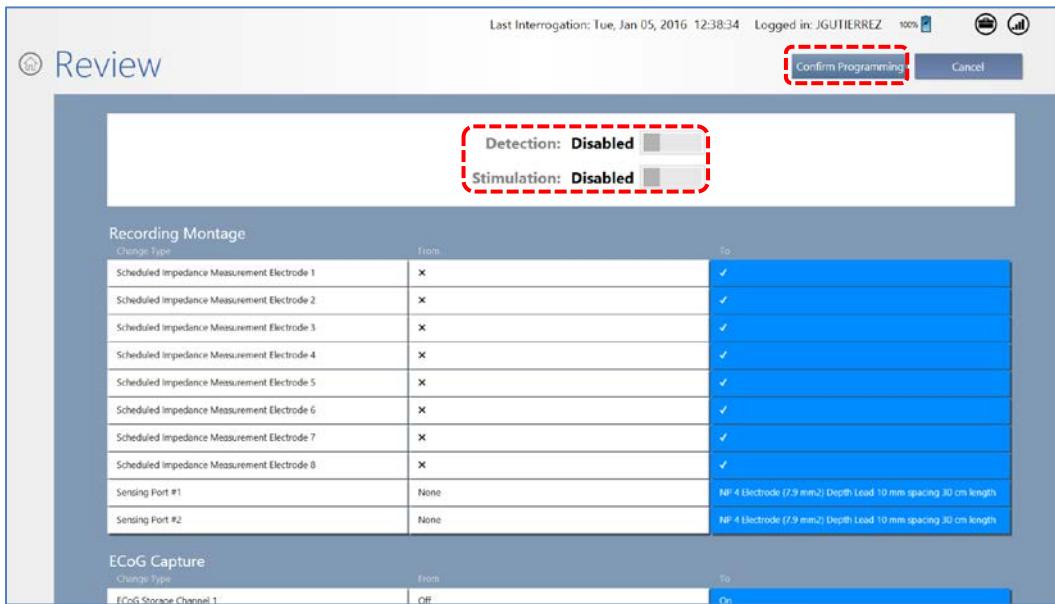

- If you want to enable detection and stimulation at this time, make sure **DETECTION** and **STIMULATION** are set to **ENABLED** at upper center.
6. Place the wand within approximately 1 inch of the neurostimulator and select **CONFIRM PROGRAMMING**. The **PROGRAMMING NEUROSTIMULATOR** dialog shows programming progress until complete. The imported settings are deleted from the tablet when programming is complete. If you do not complete this step, the settings will remain available on the tablet.

*Note: After implanting a replacement neurostimulator, monitor the detection behavior. It may be necessary to adjust the detection settings to ensure desired performance. If a bandpass detection tool (**RHYTHMIC ACTIVITY** or **SPIKE ACTIVITY**) is enabled and the settings need to be adjusted, begin by modifying the amplitude settings. See **Recommended Therapy Modifications After Observing Clinical Response** on page 85.*

PATIENT LIST ON TABLET AND THE PDMS ON A PERSONAL COMPUTER

The **PATIENT LIST** on the tablet and the PDMS on a personal computer give access to stored data for a selected patient. Their appearance and behavior is the same. The **PATIENT LIST** tile on the tablet Start screen works only when the tablet is connected to the PDMS (through the Internet).

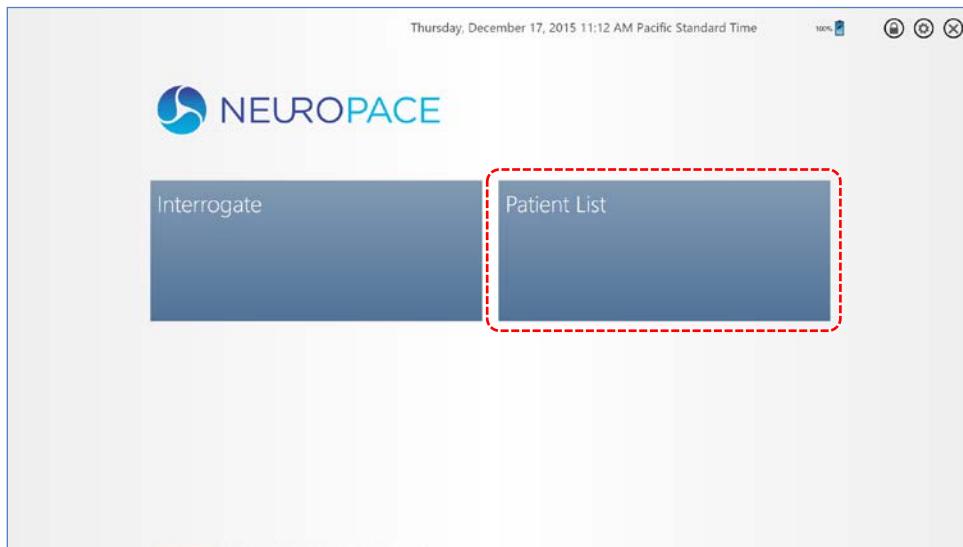


Figure 72 Tablet Start screen

Whether you select **PATIENT LIST** on the tablet Start screen or log in to the PDMS on a personal computer, the information is displayed in a similar way, and on the same-named screens, as when you start an interactive programming session by interrogating the neurostimulator. When you select **PATIENT LIST** on the tablet, the difference is that you are not interacting with a neurostimulator; you are only accessing stored neurostimulator information like detection and stimulation settings, and stored activity reports.

PDMS ACCESS ON A PERSONAL COMPUTER

Access to the secure PDMS database on a personal computer requires a username and password. If NeuroPace has not already provided these to you, contact NeuroPace to obtain a username and temporary password. You are responsible for maintaining the confidentiality of your username and password.

Note: The PDMS requires you to change your password after logging in for the first time.

1. Using any computer with Internet access, open a web browser and select the PDMS link on www.neuropace.com
2. Enter your **USERNAME** and **PASSWORD**.
3. Click **Log In**.

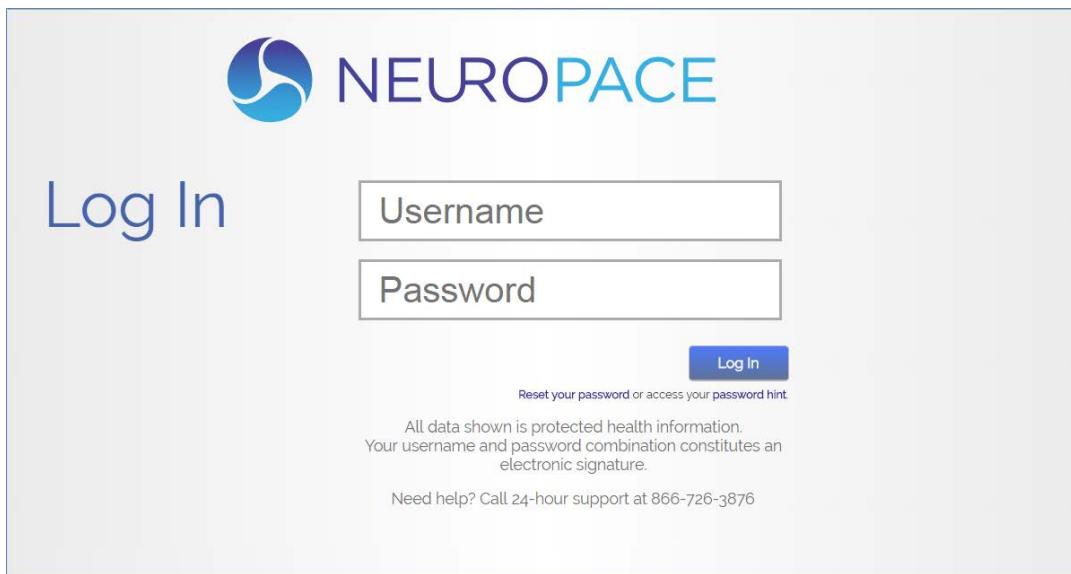


Figure 73 PDMS Login on a personal computer

*Note: If you have forgotten your password select **RESET YOUR PASSWORD OR ACCESS YOUR PASSWORD HINT**. The password hint you entered will be sent to the email address listed in your profile.*

The **PATIENT LIST** opens immediately when you log in to the PDMS using a web browser, just as when you select the **PATIENT LIST** tile on the tablet. Follow the continuing instructions below to access PDMS functions, which apply both on the tablet and on a personal computer.

PDMS PROFILE AND SETTINGS

The toolbox icon at upper right of every screen gives access to view and change your:

- Profile
- Password
- Time Zone



- To reach the **PATIENT LIST** on the tablet, select **PATIENT LIST** from the Start screen. The **PATIENT LIST** screen opens. When you log in to the PDMS via a web browser, it opens directly on the **PATIENT LIST** screen.

Alerts	Center	Last	First and Middle	Patient ID	Gender	Most Recent Implant	First Implant
! 2	NPCLIN	1303	1303	3182673	Female	03/01/2015	03/01/2015
Ω	NPCLIN	1344	1344	3368613	Female	03/01/2015	03/01/2015
Ω	NPCLIN	1365	1365	3264280	Female	03/01/2015	03/01/2015
	NPCLIN	1380	1380	3346920	Female	03/01/2015	03/01/2015
	NPCLIN	1404	1404	3335557	Female	03/01/2015	03/01/2015
Ω	NPCLIN	189142	189142	3408900	Female	03/01/2015	03/01/2015
Ω	NPCLIN	191185	191185	3367580	Female	03/01/2015	03/01/2015
Ω	NPCLIN	191392	191392	3284940	Female	03/01/2015	03/01/2015

Figure 74 Patient List screen

- Select a patient from the **PATIENT LIST**. The **PATIENT HOME** screen opens.

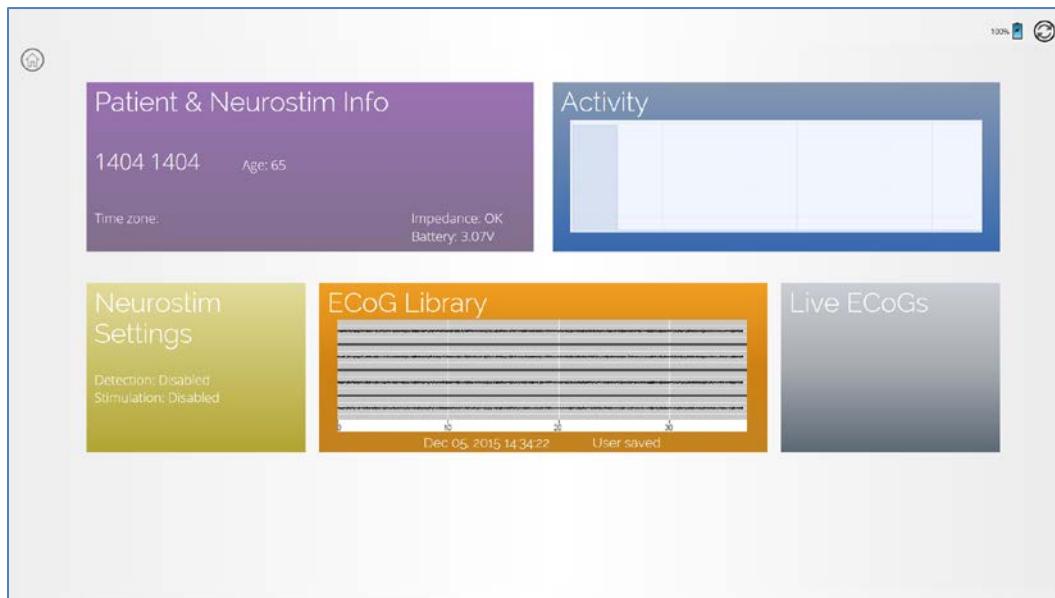


Figure 75 Patient Home screen

3. Select a tile to view its associated data for this patient.

- **PATIENT & NEUROSTIM INFO**
- **NEUROSTIM SETTINGS**
- **ACTIVITY**
- **ECoG LIBRARY**
- **LIVE ECoGs**

You can switch between patient views or return to the **PATIENT HOME** screen by touching the screen title at upper left (or the down arrow next to it).

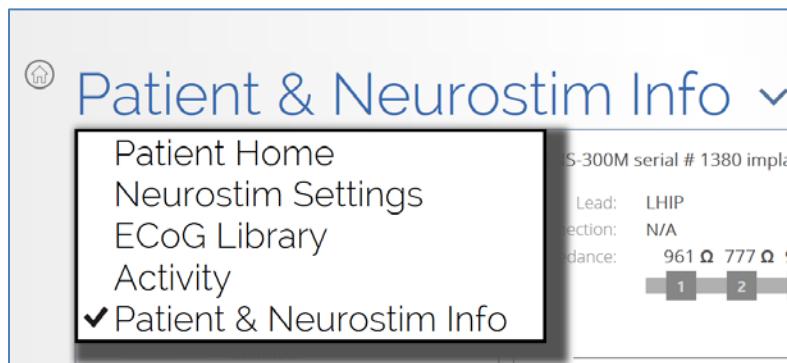


Figure 76 Select to switch between patient views

The following sections provide details about each view.

PATIENT & NEUROSTIM INFO

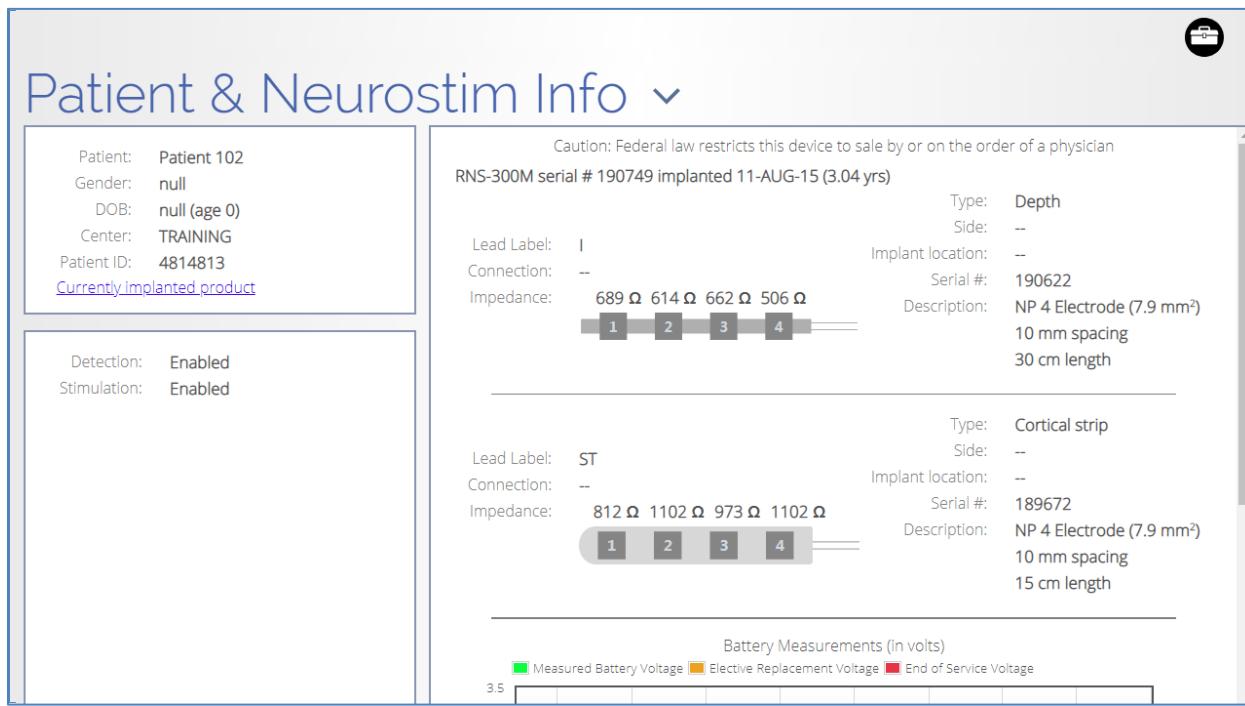


Figure 77 Patient & Neurostim Info screen

This screen provides:

- At upper left, patient information. Select **CURRENTLY IMPLANTED PRODUCT** to see details about the implanted products.
- At lower left, whether **DETECTION** and **STIMULATION** are **ENABLED** or **DISABLED**.
- At upper right, information about the neurostimulator (serial number, implant date), leads, and electrodes.
- At lower right (swipe up to see), graphs of battery and impedance measurements over time.

NEUROSTIM SETTINGS

- This **NEUROSTIM SETTINGS** screen gives access to the **PATTERN DETECTION**, **STIMULATION** and **ECoG CAPTURE** screens. These have the same names as found through the **NEUROSTIM SETTINGS** screen after you interrogate the neurostimulator. Here are the similarities and differences:
 - PATTERN DETECTION** screen: It looks exactly the same as seen in **Figure 32 Pattern Detection screen** (page 54). It behaves the same in that you can define, modify, and save pattern detection settings for the current patient, with one difference: you cannot program the new settings immediately. You can program the settings later when the programmer is in communication with the patient's neurostimulator (after interrogation). For instructions to use this screen, see **Configure Pattern Detection: Define and Modify Detection Settings** (page 50), starting with **Phase 2: Define a Pattern Detection Set** (page 54).

- When accessing the PDMS via a web browser (not via the tablet), you can define or edit and save a detection set to the PDMS, and later download the detection set to the programmer. Having saved a detection set to the programmer, you can program the neurostimulator while the programmer is offline, as is required, for example, if an operating room lacks network connectivity. Note that these instructions apply whether the target programmer is a model PGM-300 (laptop) or model 5000 (RNS® Tablet).
- **STIMULATION** screen: This screen enables you only to review the patient's history of stimulation settings. You can adjust stimulation settings only using the tablet during an interactive programming session. See **Configure Responsive Therapy** on page 76.



Figure 78 Stimulation screen through the patient list—review settings history only

- **ECoG CAPTURE screen:** This screen reports the current capture settings. You cannot adjust them here. You can adjust ECoG capture settings only using the tablet during an interactive programming session. See **Set up ECoG Capture** on page 44.

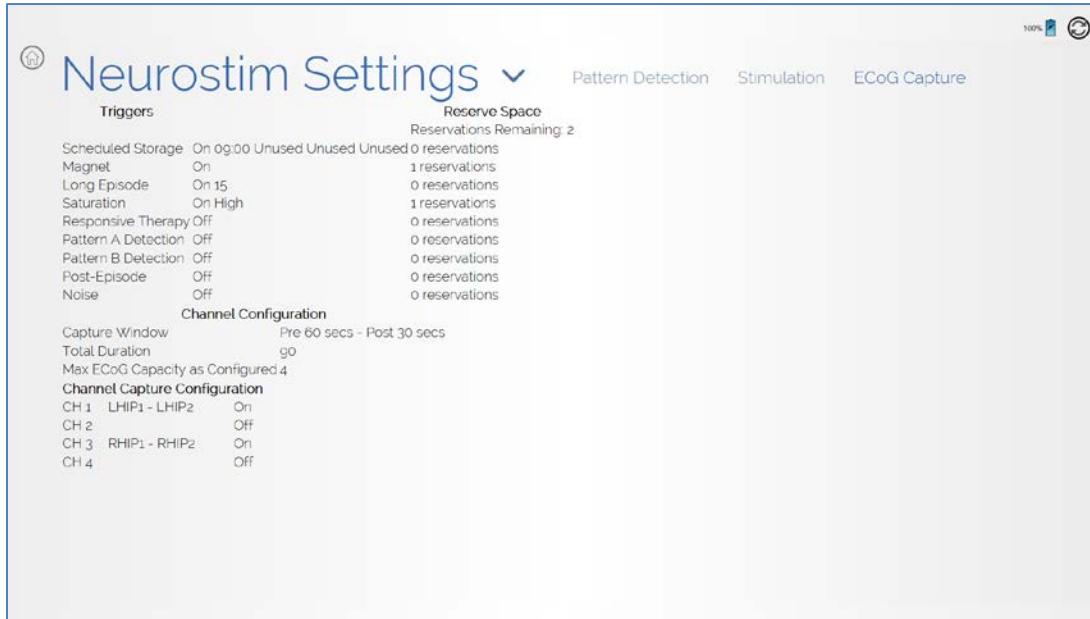


Figure 79 ECoG Capture screen through the patient list—review current settings only

ACTIVITY

The **ACTIVITY** screen is a reporting screen only. It looks and behaves as shown and explained under **Activity Screen** starting on page 28.

ECoG LIBRARY

The **ECoG LIBRARY** screen is a reporting screen only. It looks and behaves as shown and explained under **ECoG Library Screen** starting on page 34.

LIVE ECoGs

The **Live ECoGs** tile is inactive through the patient list and the PDMS via the Internet. It can be active only on the tablet, provided it is in communication with a neurostimulator, as described under **Live ECoGs Screen** starting on page 37.

PATIENT FOLLOW-UP ACTIVITIES

The following activities are recommended at a follow-up visit:

- Interrogate the RNS® Neurostimulator. See **Interrogate the RNS® Neurostimulator** on page 18.
- Review measured battery voltage and electrode impedance on the **Neurostim Info Screen**—see page 23. These are measured routinely during interrogation.

*Note: Battery voltage and electrode impedance are measured during interrogation. If these two values were not measured (e.g., because therapy was being delivered), measure them manually using the applicable **MEASURE** button on the **NEUROSTIM INFO** screen.*

- Review stored ECoGs on the **ECoG Library Screen**—see page 34—to review and consider adjustments to:
 - ECoG capture settings—see **Set up ECoG Capture** on page 44
 - Detection settings—see **Configure Pattern Detection: Define and Modify Detection Settings** on page 50
 - Responsive therapy settings—see **Configure Responsive Therapy** on page 76

*Note: You must be connected to the PDMS via the Internet to review stored ECoGs, neurostimulator activity and all other neurostimulator information in a patient's history. See **Connect and Login Overview** on page 10.*

- View live ECoGs—see **Live ECoGs Screen** on page 37.
- If you adjust detection settings, review the simulated detection on an ensemble of stored ECoGs before programming. See **Phase 2: Define a Pattern Detection Set** on page 54 or **Phase 3: Customize Pattern Detection** on page 58.
- If you adjust stimulation settings, test the patient's tolerance for the new settings before programming. See **Test Stimulation Before Enabling** on page 83.
- If you change any settings, remember to program the new settings into the neurostimulator, as described in each applicable section.

Caution: Afterdischarge Activity

If evidence of afterdischarge activity resulting from stimulation is seen either on stored ECoGs or during test stimulation delivery, stimulation parameters should be adjusted to prevent such occurrence.

TROUBLESHOOTING

Note: This section addresses only problems specific to the RNS® Tablet. For troubleshooting information regarding other components of the RNS® System, see the RNS® System manual. For troubleshooting information regarding the NeuroPace® Remote Monitor, see the NeuroPace® Remote Monitor Manual.

If the information in this section does not resolve the problem, contact NeuroPace for further assistance (see **Contacting NeuroPace** on page 5).

Network Connectivity Problems

For an overview and step-by-step instructions to establish full network connectivity with the tablet, see **Connect and Login Overview** on page 10 and following.

If you are already logged in to the tablet and have lost network connectivity, first log off to return to the login screen: return to the Start screen (see **Figure 3 Start screen** on page 18) and touch the circle-X (exit) button at upper right, and then select **EXIT** on the **PROGRAMMER EXIT** dialog. Follow the instructions under **Connect to the Internet for PDMS Access** on page 11 to establish a connection.

- If a PDMS screen appears to get “stuck” or is not responding, select the **REFRESH** button at upper right.

Note: If you attempt to connect to a Wi-Fi network that requires you to go to an Internet website to log in, as in some hotels and coffee shops, you will not be able to log in because the tablet does not have an Internet browser program. As a medical device, its Internet connectivity is restricted to the PDMS only.

Abnormal Lead Impedance (greater than 3500 Ohms or less than 250 Ohms)

Troubleshooting During an Intraoperative Procedure

The goal of the following procedures is to ascertain if the abnormal lead impedance measurement is associated with the lead or with the RNS® Neurostimulator port to which the lead is connected.

1. Repeat the impedance measurements.
2. If the impedance remains abnormal, note the electrode label associated with the abnormal impedance
3. Reconnect leads in the same neurostimulator port.
 - Disconnect and reconnect the leads to the connector cover and neurostimulator according to the instructions in the relevant section of the system manual.
 - View the live ECoGs to verify that the signal is sufficient and without noise or artifacts.
 - Repeat impedance measurements.
4. If the impedance remains abnormal, note the electrode label associated with the abnormal impedance.
5. Place leads in the opposite neurostimulator port.

- Disconnect leads, place leads in the opposite connector cover ports, and reconnect the connector cover and neurostimulator according to the instructions in the relevant section of the system manual.
- Program the new lead information and montage and use the same electrode labels for the new lead configuration.
- View the live ECoGs to verify that the signal is sufficient and without noise or artifacts.
- Repeat impedance measurements.

Note: If the abnormal impedance is associated with the same electrode label, then the problem may be with the lead.

Depth lead - Impedance > 3500 Ohms

- The depth lead may be damaged and should be replaced. Refer to instructions in the relevant section of the system manual.

Depth lead - Impedance < 250 Ohms

- The depth lead may be damaged and should be replaced. Refer to the instructions in the relevant section of the system manual.

Cortical strip lead - Impedance > 3500 Ohms

- The high impedance may be due to air.
- Irrigate the area around the cortical strip lead and examine the geometry.
- Repeat impedance measurements.
- If the impedance remains abnormal, the cortical strip lead may be damaged and should be replaced. Refer to the instructions in the relevant section of the system manual.

Cortical strip lead - Impedance < 250 Ohms

- The cortical strip lead may be damaged and should be replaced. Refer to the instructions in the relevant section of the system manual.

Note: If the abnormal impedance is not associated with the same electrode label, the problem may be with the neurostimulator or connector cover.

6. Replace the connector cover.

- Disconnect the leads from the connector cover and reconnect the leads to a new connector cover. Connect to the neurostimulator and troubleshoot according to the instructions in the relevant section of the system manual.
- View the live ECoGs to verify that the signal is sufficient and without noise or artifacts.
- Repeat impedance measurements.

7. If the impedance remains abnormal, replace the neurostimulator. Refer to instructions in the relevant section of the system manual.

Troubleshooting During a Routine Follow-up


1. Repeat lead impedance measurements to verify abnormal lead impedance.
2. Take an image (X-ray or CT scan) of the implanted lead(s) and connector cover to discern the integrity of the lead.
 - **Lead Damage** – If the lead is damaged or fractured, it is recommended to abandon and replace the lead. Refer to the instructions in the relevant section of the system manual.
 - **Poor Connection** – If the lead is not completely installed in the connector cover, it is recommended to surgically revise the connection between the connector cover and lead. Refer to the instructions in the relevant section of the system manual.

Note: It is normal for lead impedance measurements to fluctuate following the implant procedure due to the acute injury and healing processes of the brain tissue.

*Note: Lead damage may not be immediately apparent from the live ECoGs. Using the **NOISE** ECoG capture trigger may be helpful in future troubleshooting techniques.*

Insufficient Charge

Insufficient charge occurs when the neurostimulator is unable to fully deliver the programmed therapy. It may be due to high lead impedance or the neurostimulator configuration settings. If you receive an insufficient charge alert (see **Figure 9** and **Figure 10** on page 27) that cannot be explained by observed lead impedance changes, or if the first entries in the PDMS interrogation report event list show the “Programmed current” is higher than the “Lowest delivered” current, contact your NeuroPace Field Clinical Engineer. The NeuroPace Field Clinical Engineer will help you review the neurostimulator configuration settings to determine the cause. An example interrogation report is shown below in **Figure 77**.

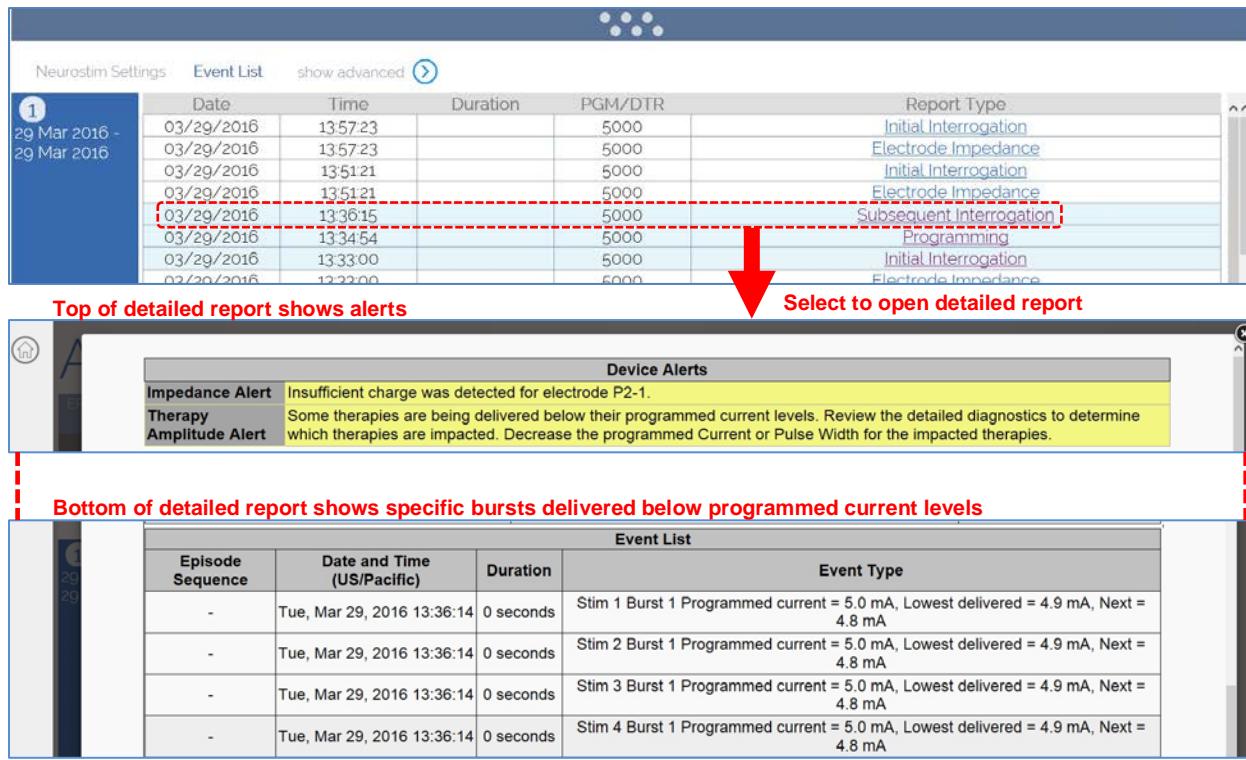

The screenshot shows a software interface for Neurostim Settings. At the top, there's a header with 'Neurostim Settings' and 'Event List'. Below this is a date range selector '29 Mar 2016 - 29 Mar 2016'. A red dashed box highlights the first few rows of the event list table. A red arrow points from the text 'Select to open detailed report' to the bottom of the event list table.

Date	Time	Duration	PGM/DTR	Report Type
03/29/2016	13:57:23		5000	Initial Interrogation
03/29/2016	13:57:23		5000	Electrode Impedance
03/29/2016	13:51:21		5000	Initial Interrogation
03/29/2016	13:51:21		5000	Electrode Impedance
03/29/2016	13:36:15		5000	Subsequent Interrogation
03/29/2016	13:34:54		5000	Programming
03/29/2016	13:33:00		5000	Initial Interrogation
03/29/2016	12:22:00		5000	Electrode Impedance

Top of detailed report shows alerts

Select to open detailed report

Bottom of detailed report shows specific bursts delivered below programmed current levels

Figure 80 Accessing an interrogation report that shows insufficient charge

Noise, Artifacts, Poor Signal Displayed, or No Signal Displayed in Live ECoG

Troubleshooting During an Intraoperative Procedure

1. If no ECoG signal is seen, make sure the recording montage has been programmed.
2. If the signal is too high in amplitude or saturating, decrease the gain setting.
3. If the leads are positioned over the RNS® Neurostimulator, move them away from the neurostimulator.
4. Remove the connector cover from the RNS® Neurostimulator using the torque driver.
5. Flush the connector cover (with leads inserted) and the RNS® Neurostimulator connector area with saline.

6. Remove the leads from the connector cover. Clean the proximal end of the lead with saline and visually inspect that the lead end is free of excess blood or other tissue.
7. Reinsert the leads. Verify that they are fully inserted into the connector cover ports.
8. Place the connector cover back onto the RNS® Neurostimulator. Verify, once again, that the leads are fully inserted in the connector cover ports.
9. Secure the connector cover onto the RNS® Neurostimulator using the torque driver.
10. View the live ECoGs to verify that the signal is sufficient and free of noise or artifacts. Isolated artifacts may be the result of electrostatic discharge (ESD) / static electricity shocks to the tablet.
11. If the signal contains noise or artifacts, disconnect the AC adapter from the tablet and run on battery power. (Make sure the battery is charged before doing so.) Disconnect the network cable as well if it is connected.
12. If the signal is not sufficient or still contains noise or artifacts, measure the impedances. If the impedances are out of range, follow the instructions in the relevant section of the system manual.
13. If impedance measurements are normal, but the problem remains, try using a different tablet, if available.
14. If the problem persists, replace the neurostimulator. Refer to the relevant section of the system manual.

Troubleshooting During a Routine Follow-up

1. Reposition the wand over the RNS® Neurostimulator and interrogate again.
2. If the signal is not sufficient, follow the **Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet** instructions in this section.
3. If the signal is not sufficient, measure the impedances. If the impedances are out of range, follow the instructions in the relevant section of the system manual.

Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet

1. Confirm that the wand cable is inserted securely into the USB port on the tablet.
2. Confirm that the concave side of the wand faces the neurostimulator and that the wand is within approximately 1 inch of the neurostimulator.
3. Disconnect the AC adapter from the tablet and run on battery power. (Make sure the battery is charged before doing so.) Disconnect the network cable as well if it is connected.
4. Move electronic equipment that may be a source of interference or move away from such equipment if possible. (Turning electronic equipment on and off may help identify the source of any interference.)
 - Reposition the wand over the neurostimulator and interrogate.
 - View the live ECoGs to verify that the signal is sufficient and void of noise or artifacts.
5. Restart (turn OFF and back ON) the tablet. Attempt the interrogation again.
6. Try a different tablet, if available.

7. Try interrogating a different neurostimulator, if available, to see if the problem is with the neurostimulator.

Tablet Does Not Turn On

1. Charge the tablet for at least four (4) hours.
2. Press and hold the power button for four (4) seconds. After 2-3 seconds, the tablet beeps to indicate it is powering on. If the tablet does not turn on, follow these instructions again. If the tablet still does not start, contact NeuroPace for further assistance (see **Contacting NeuroPace** on page 5).

Tablet Freezes or Does Not Turn Off

1. Plug the tablet into a grounded power outlet.
2. Shut down the tablet: Press and hold down the power button for 11 seconds to shut down. Then release the power button and wait 10 seconds.
3. Restart the tablet: Press and hold down the power button for four (4) seconds. After 2-3 seconds, the tablet beeps to indicate it is powering on. If the tablet does not restart, follow these instructions again. If the tablet still does not start, contact NeuroPace for further assistance (see **Contacting NeuroPace** on page 5).

Tablet Shuts Down

1. Plug the tablet into a grounded power outlet.
2. Restart the tablet: Press and hold the tablet power button. After 2-3 seconds, the tablet beeps to indicate it is powering on. If the tablet does not restart, follow these instructions again. If the tablet still does not start, contact NeuroPace for further assistance (see **Contacting NeuroPace** on page 5).

Impedance Measurement Was Rejected / Test Request Was Rejected

Message

"Impedance measurement was rejected by the neurostimulator"

"Test request was rejected by the neurostimulator"

Issue

The neurostimulator cannot measure impedance or test stimulation (test therapy).

Possible Cause

- The neurostimulator may be delivering therapy. The system cannot measure impedance or test stimulation during therapy delivery.
- The neurostimulator may be detecting an episode. The system cannot measure impedance or test stimulation during episode detection.

Actions To Be Taken

1. Wait until therapy delivery completes or temporarily disable therapy delivery.
2. Wait until episode detection completes or temporarily disable detection.

SPECIFICATIONS AND CHARACTERISTICS

RNS® TABLET

DIMENSIONS*	Width: 29 cm Height: 21 cm Depth: 1 cm
MATERIAL	Commercial Materials
WEIGHT*	0.86 kg
POWER SOURCE	100 – 240 VAC, 1.7 A or internal rechargeable battery
USB POWER SUPPLY	5V, 500 mA
OPERATING TEMPERATURE	0 °C to 40 °C

* Typical, approximate values. The values for your tablet may be different.

WAND (MODEL W-02)

DIMENSIONS (LENGTH X WIDTH X DEPTH)	7" x 3.5" x 1.3" (18 cm x 9 cm x 3 cm)
WEIGHT	0.4 pounds (181 g)
POWER SOURCE	USB port of the remote monitor laptop
OPERATING CONDITIONS	Temperature: 32 to 95 °F (0 to 35 °C) Humidity: 15 to 90%, non-condensing Atmospheric pressure: 700 to 1060 hPa
MATERIAL	ABS copolymer
LEAST FAVORABLE WORKING CONDITIONS	Wand output power and data rate vary with communication distance. Communication at a far distance (3 cm), indicated by a low signal level when using the wand, results in the slowest rate of transmission at the highest output power.
STORAGE AND TRANSPORT TEMPERATURE	32 to 140 °F (0 to 60 °C)
EXPECTED SERVICE LIFE	5 years

TABLET WIRELESS

WIRELESS FUNCTION	Transfer data between the Neurostimulator and the Patient Data Management System (PDMS)
WIRELESS TECHNOLOGY TYPE	Wi-Fi dual band 802.11n and 802.11ac 4G (fallback to 3G)
OPERATING DISTANCE	Varies - Wi-Fi and cellular range similar to other smartphones and tablets
FREQUENCY BAND OF OPERATION	Wi-Fi: 2.4 GHz or 5 GHz Cellular: 700, 1900, 1700, 2100 MHz
INTENDED USE ENVIRONMENT	Clinical setting only
REQUIRED WIRELESS QUALITY OF SERVICE	The requirements for quality of service vary depending on the clinical function. The tablet login screen displays the network status indicator on the lower right hand corner of the screen. The indicator is a signal strength symbol consistent with commercial tablet computers.
RECOMMENDED WIRELESS SECURITY	There are no requirements for timely data transmission. The tablet locally stores data until a connection to PDMS is available. Data are resent until they are successfully received.
WIRELESS FUNCTION	Wi-Fi: WPA2

RNS® SYSTEM WIRELESS

WIRELESS FUNCTION	Transfer data between the neurostimulator and the programmer or remote monitor and wand: interrogation, programming and live ECoGs
WIRELESS TECHNOLOGY TYPE	Short range, low power inductive coil to coil telemetry
INTENDED USE ENVIRONMENT	Clinical setting and home environment
OPERATING RANGE	0 – 3 cm
FREQUENCY BAND	20 kHz – 50 kHz
RECEIVE BANDWIDTH OF THE NEUROSTIMULATOR	100 kHz
RECEIVE BANDWIDTH OF THE WAND	≥ 50 kHz
NUMBER OF CHANNELS	Single Channel
MODULATION TYPE	On/Off pulse amplitude modulation
RF DATA FLOW CHARACTERISTICS	Half duplex
EFFECTIVE RADIATED POWER	37.25 nW or less

ELECTROMAGNETIC EMISSIONS AND IMMUNITY

Electromagnetic interference (EMI) is a field of energy generated by equipment found in the home, work, medical, or public environments that is strong enough to interfere with neurostimulator function. The RNS® System is designed to be immune from common sources of electromagnetic interference. The most common sources of EMI are discussed in this section.

The neurostimulator is intended to sense electrical activity of the brain, detect activity of interest, deliver stimulation therapy, and communicate with a programmer. In addition, the neurostimulator is intended to operate autonomously, conserve battery life, check for proper connectivity to leads, and automatically recover from certain conditions.

Electromagnetic disturbances may be detected by the neurostimulator as noise or saturation and therefore suppress the ability of the neurostimulator to detect activity of interest during the exposure to the electromagnetic field. Conversely, electromagnetic disturbances may be detected as activity of interest and result in the delivery of programmed stimulation to the patient. Electromagnetic disturbances will not cause the neurostimulator to change modes or operating parameters.

Electromagnetic disturbances will not cause the neurostimulator to exceed prescribed therapy limits or charge density limits. Electromagnetic disturbances can be expected to reduce the ability of the neurostimulator to communicate with the programmer. If this happens, refer to **Poor or No Communication Between the RNS® Neurostimulator and the RNS® Tablet** on page 102. The neurostimulator is expected to recover all intended functions following the removal of the electromagnetic disturbance.

Hospital or Medical Environments

Patients should always inform healthcare personnel that they have an implanted RNS® System (and show their medical implant identification card) before any procedure is performed. Most diagnostic procedures, such as x-rays and ultrasounds, may be performed without affecting the RNS® System. However other diagnostic and therapeutic equipment with higher energy levels may interfere with the RNS® System. Refer to **Warnings and Cautions – Programmer** on page 6 for specific information.

Home, Work or Public Environments

The patient should avoid or exercise caution when in the presence of the following potential sources of EMI that may affect the operation of the neurostimulator:

- Radiofrequency identification (RFID) sources
- Airport security and other surveillance systems
- Power lines and transmission towers
- Electric substations, power generators and large transformers
- Portable and mobile RF communications equipment
- Electric arc welding equipment
- Electric steel furnaces
- Electric induction heaters
- Electric fences
- Body fat measurement scales
- Jackhammers
- Stun guns

The following commonly used items should not affect the operation of the neurostimulator:

- Cell phones and Bluetooth devices

- Electric toothbrushes, electric shavers, and hair trimmers
- Microwave ovens
- Appliances such as washing machines, dryers, electric stoves, toasters, blenders, electric can openers, and food processors
- Electric blankets and heating pads
- Personal computers, electric typewriters, copiers, and fax machines
- Televisions, AM/FM radios, stereos, personal music players
- Vacuum cleaners and electric brooms

For additional information about devices that generate electromagnetic interference contact NeuroPace (see **Contacting NeuroPace** on page 5). If a patient suspects EMI is disrupting the operations of their neurostimulator, advise the patient to move away from the source of the EMI.

Guidance and Manufacturer's Declaration

Medical electrical equipment needs special precautions regarding EMC and needs to be installed and put into service according to EMC information provided. This declaration applies for the following devices:

- RNS® Neurostimulator, models RNS-300M and RNS-320
- RNS® Tablet, model 5000
- Wand, model W-02
- Remote Monitor, models DTR-300 and DTR-300-E, and 5100

The devices comply with IEC 60601-1-2, ISO 14708-3, and FCC 47 CFR Parts 2 and 15.

- Portable and mobile RF communications equipment can affect the device.
- This equipment/system may cause radio interference or may disrupt the operation of nearby equipment. The devices may be interfered with by other equipment, even if that other equipment complies with CISPR emission requirements. It may be necessary to take mitigation measures, such as re-orienting or relocating the devices or shielding their location.

Emissions and Immunity Information

The programmer with wand is designed for use in a medical environment by a trained healthcare professional. The remote monitor with wand is designed for use in the home by a patient. The devices are intended for use in the electromagnetic environment specified below. The customer or user of the system should assure they are used in such an environment.

Note: Unless otherwise indicated in the table footnotes, emissions testing information in the tables below apply to all of the devices addressed in this section as listed above.

**Table 1. Guidance and manufacturer's declaration –
electromagnetic emissions – for all equipment and systems**

Emissions test	Compliance	Electromagnetic environment – guidance
Conducted emissions (CISPR 11)	Class B, Group 1 150 kHz to 30 MHz	The RNS® System uses RF energy only for its internal function. Nearby electronic equipment may be affected.
RF emissions (CISPR 11)	Class B, Group 1 30 MHz to 1 GHz	
Harmonic emissions (IEC 61000-3-2)	Class A Device ^a	The RNS® System is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations / flicker emissions (IEC 61000-3-3)	Limits per Clause 5 of the Standard ^a	

- a. Tablet (model 5000), Remote Monitor (model 5100), with Wand (model W-02) tested for harmonic emissions and flicker. Other products excluded from harmonic emissions and flicker testing.

Table 2. Guidance and manufacturer's declaration – electromagnetic immunity – for all equipment and systems

Immunity test	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) (IEC 61000-4-2)	±2, 4, 6, 8kV contact discharge ^a ± 2, 4, 8, 15kV air discharge ^a	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient / burst (IEC 61000-4-4)	± 2 kV for power supply lines ± 1 kV for input / output lines ^b	Mains power quality should be that of a typical commercial or hospital environment.
Surge (IEC 61000-4-5)	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth ^c	
Voltage dips, short interruptions and voltage variations on power supply input lines (IEC 61000-4-11)	0% U _T 0.5 cycle 0% U _T 1 cycle ^d 40% U _T 5 cycle ^d 70% U _T 25 cycles 0% U _T 5 Sec	If the user of the programmer and wand requires continued operation during power mains interruptions, it is recommended that the programmer and wand be powered from an uninterruptible power supply or a battery.
Magnetic immunity (IEC 61000-4-8)	3, 30 ^e A/m	

^a Immunity to ESD tested to ±8 kV contact and ±15 kV air with Tablet (model 5000), Remote Monitor (model 5100), and Wand (model W-02).

^b Immunity to electrical fast transients on I/O lines tested with Tablet (model 5000), Remote Monitor (model 5100), and Wand (model W-02).

^c Immunity to surge line to earth tested with Programmer (model PGM-300) and Remote Monitor (model DTR-300-E).

^d Immunity to voltage dips with compliance to 0% U_T 1 cycle for Tablet (model 5000), Remote Monitor (model 5100), and Wand (model W-02) and to 40% U_T 5 cycle for Programmer (model PGM-300) and Remote Monitor (model DTR-300-E).

^e Magnetic immunity tested with Tablet (model 5000), Remote Monitor (model 5100), and Wand (model W-02) to 30 A/m.

Table 3. Guidance and manufacturer's declaration – electromagnetic immunity – for equipment and systems that are not life supporting

Immunity Test	Compliance Level	Electromagnetic environment – guidance
Conducted RF (IEC 61000-4-6)	3 Vrms 150 KHz to 80 MHz 6 Vrms ISM bands ^a	Portable and mobile RF communications equipment should be used no closer to any part of the RNS® System, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance ^{b, c} $d = 1.17\sqrt{P}$ (80 MHz to 800MHz) $d = 2.33\sqrt{P}$ (800 MHz to 2.7 GHz) Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).
Radiated RF (IEC 61000-4-3)	3, 10, 20 V/m 80 MHz to 2.7 GHz ^e	Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ^d , should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 

^a Conducted immunity of the Tablet (model 5000), Remote Monitor (model 5100), and Wand (model W-02) compliant to 6 V/m in the ISM bands.

^b Separation distance relevant to Programmer (model PGM-300) and Remote Monitor (model DTR-300). At 80 MHz and 800 MHz, the higher frequency range applies.

^c These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^d Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitter, an electromagnetic site survey should be considered. If the measured field strength in the location in which the RNS® System is used exceeds the applicable RF compliance level above, the RNS® System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the RNS® System.

^e Radiated immunity of the Programmer (model PGM-300), Remote Monitor (model DTR-300-E) compliant to 3 V/m up to 2.5 GHz. Tablet (model 5000), Remote Monitor (model 5100) and Wand (model W-02) compliant to 10 V/m to 2.7 GHz and spot frequencies with pulse modulation and to 20 V/m in the communication passband. Spot frequencies tested are: 385 MHz, 450 MHz, 710 MHz, 745 MHz, 780 MHz, 810 MHz, 870 MHz, 930 MHz, 1.720 GHz, 1.845 GHz, 1.970 GHz, 2.450 GHz, 5.240 GHz, 5.500 GHz, 5.785 GHz.

Table 4. Recommended separation distances between portable and mobile RF communications equipment and the RNS® System^{a, b, c}

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = 1.17\sqrt{P}$	$d = 1.17\sqrt{P}$	$d = 2.33\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
1	1.17	1.17	2.33
10	3.69	3.69	7.38
100	11.67	11.67	23.33

^a All separation distances relevant to Programmer (model PGM-300) and Remote Monitor (model DTR-300). At 80 MHz and 800 MHz, the higher frequency range applies.

^b These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^c For transmitter rated at a maximum output power not listed above, the recommended separation distance of d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

GLOSSARY

Adapt Channel. The adapt channel is an advanced responsive therapy setting. The neurostimulator will monitor this channel from which to calculate the adaptive stimulation frequency. This channel is used for all adaptive bursts.

Adapt Min Amplitude. The percentage of the detected signal amplitude to be used as the minimum therapy amplitude.

Adaptation. The neurostimulator may be programmed to deliver the pulse interval for a responsive therapy burst by adapting to the electrographic signal by a selected percentage.

Anode. The positive pole of an electrode selected for stimulation.

Area Detection Tool. The area detection tool measures the area under the curve of an ECoG signal.

Asynchronous Pulse. Delivery of the first pulse in a burst independent of the ECoG activity.

Bandpass Detection Tool. The bandpass detection tool measures the approximate frequency of the signal based on the duration and amplitude of waveform segments.

Battery Measurement. The neurostimulator battery voltage may be measured interactively through the programmer.

Burr Hole. Hole made in the skull used for the insertion and securing of implanted cortical and/or depth leads.

Burr Hole Cap. The burr hole cap is intended to secure the leads in relation to the cranium.

Burr Hole Cover. The burr hole cover is intended to secure the leads to the burr hole cap.

Burst. An individually programmed group of one to many pulse(s) of current delivered by the neurostimulator.

Burst Duration. The amount of time during which a burst will be delivered by the neurostimulator.

Burst 1. The first burst of stimulation delivered within an individual responsive therapy. This burst may be delivered asynchronously or synchronously.

Burst 2. The second burst of stimulation delivered within an individual responsive therapy. This burst may be programmed ON or OFF for each therapy and will always be delivered asynchronously following Burst 1 of the therapy if programmed ON.

Can. “Can” is used to describe the neurostimulator titanium body. The can may be used as either an anode (+) or cathode (−) in a stimulation pathway.

Cathode. The negative pole of an electrode selected for stimulation.

Channel. The amplifier channel of electrodes selected to detect ECoG activity.

Connector Cover. The connector cover is intended to secure the proximal lead contacts and lead strain relief to the body of the neurostimulator.

Connector Cover Port. The connector cover ports contain the point of contact between the proximal end of the lead(s) and the neurostimulator. The connector cover contains two (2) ports, one for each lead.

Connector Plug. The connector plugs are used to fill any ports on the connector head of the neurostimulator that are not in use.

Cortical Strip Lead. The NeuroPace® Cortical Strip Lead lies on the surface of the brain (cortex) and can detect the electrical activity of the brain and deliver stimulation.

Cranial Prosthesis. The cranial prosthesis is intended to occupy the vacant craniectomy after the neurostimulator has been explanted.

Craniectomy. Defect formed in the skull in order to implant the neurostimulator.

Craniectomy Template. The craniectomy template is intended to be used as a pattern to mark and delineate the shape of the ferrule on the skull prior to making the craniectomy.

Current. The amplitude of current that will be delivered in each phase of a current pulse.

Delay. A programmed delay in the delivery of the first pulse of a burst, measured from the time of synchronizing to a signal peak.

Default Therapy Frequency. The default therapy frequency is an advanced setting available for each responsive therapy burst. For an adaptive therapy burst (Adaptation = ON), if the calculated stimulation frequency does not fall between the programmed min and max frequencies, then the stimulation is delivered at the default therapy frequency.

Delay. A programmed delay in the delivery of the first pulse of a burst, measured from the time of synchronizing to a signal peak.

Depth Lead. The NeuroPace® Depth Lead is implanted into the brain and can detect the electrical activity of the brain and deliver stimulation.

Detection Settings. Settings that combine patterns and detection tool settings in order to detect activity.

Detection Tools. There are three (3) types of detection tools selectable for Pattern A or Pattern B detectors. The detection tools are: the bandpass tool, power change tool (also known as the line length tool), and area tool.

ECoG. Electrocorticogram. Electrical activity derived directly from the cerebral cortex. Also used to describe the neurostimulator or programmer stored record of this activity (e.g. "ECoG record").

ECoG Activity. ECoG activity refers to activity which is detected by the neurostimulator using one of the event detectors.

ECoG Trigger. The type of neurostimulator activity selected by the user to result in the capture of an ECoG record if the activity occurs.

Elective Replacement Indicator (ERI). The ERI indicator is displayed on the programmer to notify the user that the neurostimulator is approximately three months prior to end-of-service (EOS). Neurostimulator elective replacement is suggested at a battery voltage equal to the ERI.

End-Of-Service (EOS). End of service is the point in service at which the battery level is still sufficient to allow operation, however some features may be lost.

Estimated Charge Density. The charge density of a selected and/or programmed stimulation is estimated by the programmer. Charge densities greater than 25 µC/sq cm are not allowed to be selected or programmed into the neurostimulator.

Event. Activity detected by the neurostimulator. Events are considered any of the following: long episodes, pattern A, pattern B, responsive therapy, noise or saturation.

Ferrule. The ferrule is intended to be installed in a craniectomy to secure and mechanically support the neurostimulator in the cranium.

Ferrule Clamp. The ferrule clamp is used to secure the neurostimulator to the ferrule.

Frequency. The frequency is the rate at which pulses are delivered which is expressed in Hz (pulses/second).

Gain. ECoG signal amplification on one channel of the neurostimulator.

Impedance Measurement. The opposition of current flow between electrodes as determined by measuring the voltage (Ohms) resulting from the delivery of a current pulse. A current pulse must be applied to the patient in order to obtain this measurement.

Interrogate. The act of retrieving stored data and settings from the neurostimulator using the programmer.

Lead Cap. The lead cap is intended to physically protect the proximal ends of implanted leads that are not connected to the neurostimulator.

Lead ID. Numerical and/or alpha characters assigned to a lead to describe its anatomical location.

Lead Strain Relief. The lead strain relief is intended to mechanically support the proximal ends of the lead at their exit point from the neurostimulator.

Line Length Detection Tool: See Power Change Detection Tool.

Live ECoG. ECoG data viewed in real-time through the programmer.

Long Episode. Initiate ECoG storage when a detected episode continues beyond a preset duration.

Magnet. The magnet can be used by a patient to mark a clinical seizure, and depending on how the physician programs the neurostimulator, to trigger electrocorticographic (ECoG) storage and withhold therapy.

Max Therapy Frequency. The max therapy frequency is an advanced setting available for each responsive therapy burst. For an adaptive therapy burst (Adaptation = ON), the max therapy frequency specifies the maximum stimulation frequency. If the calculated stimulation frequency does not fall between the programmed min and max frequencies, then the stimulation is delivered at the default therapy frequency.

Min Therapy Frequency. The min therapy frequency is an advanced setting available for each responsive therapy burst. For an adaptive therapy burst (Adaptation = ON), the min therapy frequency specifies the minimum stimulation frequency. If the calculated stimulation frequency does not fall between the programmed min and max frequencies, then the stimulation is delivered at the default therapy frequency.

Montage. Collective term used to describe the assignment of electrode labels and gain settings to the amplifier channels.

Neurostim Info. Programmer display of measurements of the neurostimulator lead impedances and battery voltage.

Neurostimulator. The RNS® Neurostimulator is an implantable, battery powered, microprocessor controlled device that can amplify and analyze the patient's electrocorticographic activity, detect activity from intracranial electrodes and deliver a short train of current pulses to the brain to attempt to interrupt the detected activity.

Noise. 1. An event detector that identifies 60 Hz noise from selected detection channels. 2. Environmental interference in the telemetry between the neurostimulator and wand/programmer.

Overwritten. Once the memory for ECoG or diagnostics is full, the neurostimulator will store new data by overwriting the oldest ECoG records or diagnostics. Reserved ECoGs are only overwritten by newer ECoGs of the same ECoG capture trigger type.

Patient Data Management System (PDMS). The NeuroPace® Patient Data Management System is a secure database that provides a means for a clinician to review information that has been transmitted by the programmer and the remote monitor.

Pattern Detection. Analysis of detection settings using the programmer detection tools and retrieved neurostimulator ECoG records.

Patterns. Patterns (pattern A and pattern B) use a collection of detection tools to classify electrographic activity. Patterns A1 and B1 apply to the first chosen channel and Patterns A2 and B2 to the second chosen channel.

Post-Episode Interval. Responsive therapies will be DISABLED for this period of time after detecting the end of an episode.

Power Change Detection Tool. The power change detection tool (also known as the line length detection tool) performs an estimate of the length dimension of a signal related to the complexity of a signal.

Programmer. The programmer provides the clinician with a user interface to select and download operating parameters to the neurostimulator for detection and responsive stimulation settings, to view live ECoG signals, to test the RNS® System integrity, and to upload data and diagnostic information from the neurostimulator for viewing.

Programming. Using the programmer to program settings into the neurostimulator.

Pulse. A biphasic waveform of electrical stimulation delivered by the neurostimulator.

Pulse-Width Per Phase (PW). Duration of a single phase within a biphasic pulse, measured in milliseconds.

Reset. When the neurostimulator resets, detection and all therapies are disabled, and the neurostimulator stops storing ECoGs and diagnostics. (However, events prior to reset will be saved.) If you interrogate a neurostimulator that has reset, the programmer displays one or more alerts that describe the cause of the reset. When the type of reset is recoverable, the programmer displays instructions to recover, which typically guide you to reprogram the neurostimulator. If the reset is not recoverable, or if recovery is unsuccessful, note the cause of the reset and contact NeuroPace (see **Contacting NeuroPace** on page 5).

Responsive Stimulation. Electrical stimulation output to cortical tissue by the RNS® Neurostimulator in response to a detection.

Saturation Event Detector. Saturation occurs when an input signal exceeds the dynamic range of a particular channel. The neurostimulator incorporates a saturation detector to reduce occurrences of detecting saturation as a neurological event.

Scheduled Storage. A selectable neurostimulator function that causes the neurostimulator to capture an ECoG at the desired time(s) of day.

SimpleStart Feature. This feature allows the physician to easily begin the process of defining detection parameters by touching or clicking on the location in the ECoG where the physician would like detection to occur.

Simulate Performance. Simulation uses the programmer to test a proposed set of detection settings with retrieved ECoGs. The simulation will generate markers, simulating how the neurostimulator may detect if those settings are programmed into the neurostimulator.

Start Time. Scheduled storage setting controlling the time(s) at which ECoG capture will be triggered.

Stereotactic Frame. A medical device that is attached to the patient's head and provides a three-dimensional frame of reference for accurate implantation of depth leads.

Stereotactic Surgery. Surgery in which a system of three-dimensional coordinates is used to locate the site to be operated on.

Stim Pathway. Collective term describing the anodes [+] and cathodes [-] selected for a burst.

Stop Gauge. The stop gauge is placed on a depth lead prior to implantation to indicate the appropriate depth of its insertion.

Stylet. Support wire contained in both the cortical strip and depth leads. The stylet is removed after proper lead placement has been achieved.

Suture Sleeve. The suture sleeve is intended to protect the lead body at the point at which the lead is sutured to subcutaneous tissue to immobilize it and prevent its dislodgement while chronically implanted.

Synch Channel. The neurostimulator will monitor this channel for ECoG signal peaks to calculate synchronous delivery timing. One channel will be used for all bursts selected to have synchronous pulse delivery.

Synchronization. Delivering the first pulse in a burst or after a first pulse delay (synchronous) from a peak in an ECoG signal.

Technical Parameters. Option to configure detection using all detection parameters.

Telemetry. Communication between the neurostimulator and the programmer.

Test Stimulation. A stimulation that is configured and delivered in real-time using a programmer.

Therapy Frequency Multiplier. Therapy frequency multiplier is an advanced responsive therapy setting. The therapy frequency multiplier specifies the ratio of the stimulation frequency to the sensed signal frequency.

Therapy Limit Per Day. Responsive therapy delivery will be DISABLED for the remainder of the day after the programmed number of responsive therapies has been delivered.

Therapy Sequence. Responsive therapy is delivered as a therapy sequence of up to 5 individually configured sequential therapies (electrical stimulation) in response to each detected episode.

Tunneling. Part of the implant procedure during which the proximal end of an implanted lead is tunneled from the burr hole location to the neurostimulator implant location to achieve a lead/neurostimulator connection.

Tunneling Tool. The cranial tunneling tool is intended to tunnel implanted leads from their point of exit from the cranium through a sub-galeal pathway to the implanted neurostimulator location.

INDEX

- 2-channel detection set, 57
Abnormal Lead Impedance, 27, 98
About this manual, ii
ACCEPT CHANGES, 83
Access PDMS, 10
Activity histogram, 29
Activity screen, 28, 96
ADAPT CHANNEL, 81
ADAPT MIN AMPLITUDE, 81
ADAPTATION, 81
Adaptive Therapy, 81
ADV (Advanced) column, 32
ADV CONTROLS, 68, 69
Advanced Controls, 69
ADVANCED SETTINGS dialog, 69
Advanced Settings—therapy, 80
Afterdischarge Activity—caution, 8, 84, 97
AMPLITUDE, 65, 66
AMPLITUDE (%) slider, 57, 61
AREA, 69
Area Detection Tool, 68
Artifacts, 41, 98, 99, 101, 102
ASSIGN CATEGORY field, 36
Assign Channels, 41
BACK button, 17
BANDPASS, 69
Bandpass detection tool, 66
Bandpass Detection Tool, 67
Battery Measurement, 27
Battery Voltage, 26
burst, 77
BURST DURATION, 78
bursts or burst, 75
Capture Duration, 46
Capture ECoGs, 44
Capture Triggers, 46
CAPTURE WINDOW field, 36, 46
CATEGORIES, 35
CHANGE CHANNEL button, 57
CHANGE RECORDING MONTAGE, 24, 40
CHANNEL(S) (CH 1 through CH 4), 46
CLEAR button on DEFINE PATTERNS dialog, 57, 61
Clinical Response and recommended therapy modifications, 86
Computerized Tomography (CT), 100
CONFIGURE CHANNELS area, 46
Configure Pattern Detection, 37, 50, 94, 97
Configure Responsive Therapy, 76
CONFIRM PROGRAMMING, 43, 49, 52, 58, 84, 89
CONFIRM PROGRAMMING button, 61
Connect to the Internet, 10, 11, 28, 34, 54, 98
Connectivity Problems, 98
Contacting NeuroPace, 5
CONTINUE button, 53, 61
Create the Recording Montage, 38
Creating the Montage, 39
CURRENT, 41, 78
Customize Pattern Detection, 58
date range, 35
Default Detection Set, 50, 51
DEFAULT THERAPY FREQUENCY, 81

- Define a Pattern Detection Set, 54
Define and Modify Detection Settings, 50
DEFINE PATTERNS dialog, 56, 59
DELAY, 82
DETECT slider, 57, 61
Detection parameters, 57, 61
detection parameters affect sensitivity, 57
DETECTION set, 31
DETECTION SETS, 35
DETECTION SETTINGS—share, 70
Detection Tools and Detectors, 62
Detector Power Change, 63
Detector Spike Activity, 66
Device Description, 4
DURATION, 65, 66
DURATION area, 36, 46
DURATION slider, 57, 61
ECoG
 Artifact, 41
 Real-time, 97
 Records, 97
ECoG Capture, 44
ECoG CAPTURE screen, 36, 96
ECoG CAPTURE screen, 45
ECoG CAPTURE screen, 45
ECoG Capture Settings—Initial, 47
ECoG Capture Settings—Review, 49
ECoG LIBRARY screen, 96
ECoG Library Screen, 34
ECoG Records, 34
ECoG, Artifacts, 101
ECoG, Live, 101
ECoGs—Live, 37, 41
Electrode 1 is most distal, 39
Electrode Impedance, 26
electrodes numbered 1 through 4, 39
Electromagnetic Emissions and Immunity, 105
Emissions and Immunity Information, 108
Enable Therapy: Review & Program, 84
EPOCHS - PROGRAMMING, 29
ESTIMATED CHARGE DENSITY, 78
EVENT LIST, 33
FAVORITES, 35
FCC Information, iii
FLAGS, 29
FREQUENCY, 78
FROM and **To** fields, 35
Functionality when Not Connected to PDMS, 11
GAIN, 36, 41
Glossary, 112
Heating—caution, 8
HIDE button on DEFINE PATTERNS dialog, 57, 61
hide the simulation, 61
Home Screen, 21
Identifying Ports, Leads and Electrodes on the Neurostimulator, 38
Impedance, 27, 98
Impedance measurement, 26, 103
impedance significance, 27
IMPLANTS, 29
Incompatibility of Programmer with Other Medical Devices—caution, 8
Initial ECoG Capture Settings, 47
Initial pattern detection settings, 51
Initial Responsive Therapy Settings, 85
Insufficient Charge, 101
Interference, iii, 19, 101, 102, 107, 110

- Interference, RFID—warning, 7
Interference—caution, 8
Internet connection, 11, 28, 34, 54, 98
Interrogation, 18
INTERVAL, 41
LAST PROGRAMMED date, 23
Lead damage, 100
Lead Type, 40
Live ECoG, 101
Live ECoG artifact—caution, 7
Live ECoGs Screen, 37
Log In screen, 17
Log In to the PDMS, 17
Login overview, 10
Login to the Patient Data Management System (PDMS), 10
Login to the RNS® Tablet, 14
LONG EPISODE trigger, 46
lower parameter values result in greater sensitivity, 57
MAGNET, 46, 49
MAGNET SWIPE, 29
MAGNET THERAPY WITHHOLDING, 79
MAGNET trigger, 46
MAX THERAPY FREQUENCY, 81
MAXIMUM FREQUENCY, 65, 66
MEASURED BATTERY VOLTAGE, 24
MIN THERAPY FREQUENCY, 81
MINIMUM FREQUENCY, 65, 66
Mobile Broadband Works Automatically, 11
MODEL, 23
MODIFY buttons (**A1 & A2**, **B1 & B2**), 55
Modify Detection Settings, 50
Montage—Create, 38
MORE CONTROLS button, 57, 61
NAME field, 54, 59
Network Connectivity Problems, 98
NETWORKS panel, 12
NeuroPace components—warning, 7
NeuroPace Update Center, 15
NEUROSTIM EVENTS, 29
Neurostim Info Screen, 23
Neurostim Settings, 94
NEUROSTIM SETTINGS screen, 52
Neurostimulator Interrogate, 18
Neurostimulator Wireless Specifications, 105
Neurostimulator, Configure Pattern Detection, 50
Neurostimulator—Capture ECoGs, 44
Neurostimulator—Initial Montage Settings, 41
Noise, 19, 101
NOISE trigger, 47
Noise—caution, 8
OFFICE INTERROGATION, 29
ON/OFF switch, 46
Optional Advanced Settings—therapy, 80
OPTIONS area, 79
Overview of Responsive Therapy, 75
parameter values to the right of sliders, 57
parameters—detection, 57, 61
Patient & Neurostim Info, 94
Patient Follow-Up, 97
PATIENT HOME screen, 92
Patient List on Tablet and PDMS, 10, 90
PATIENT LIST screen, 92
PATTERN A DETECTION trigger, 47
PATTERN B DETECTION trigger, 47
pattern detection has three phases, 50

- Pattern Detection Phase 1**, 50, 51
Pattern Detection Phase 2, 54
Pattern Detection Phase 3, 58
PATTERN DETECTION screen, 59, 94
PATTERNS A1 & B1, 54, 59
PATTERNS A2 & B2, 54, 59
PDMS login, 10, 17
PDMS, access Internet, 10
PDMS, Internet connection, 10
Phase 1: Select a Default Detection Set upon First Use, 51
Phase 2: Define a Pattern Detection Set, 54
Phase 3: Customize Pattern Detection, 58
Placement of the Programmer Power Cords—caution, 8
Poor Connection of a lead, 100
port #1 and #2 locations, 38
PORT #1 and PORT #2, 24
PORT #1 and PORT #2—IDENTIFYING, 38
Portable Radio Frequency (RF) Communications Equipment—warning, 7
Post-Episode Interval, 82
POST-EPIISODE trigger, 47
Potential Shock—warning, 6
POWER CHANGE, 57, 61, 69
Power Change Detector, 63
Power On the Tablet, 10
Products and Accessories, 4
Program Neurostimulator, 43, 49
Programmer Failure—caution, 8
Programmer login, 14
programming epoch, 28
PW PER PHASE, 78
Recommendation for Review of Initial ECoG Capture Settings, 49
Recommended Initial ECoG Capture Settings, 47
Recommended Initial Montage Settings, 41
Recommended Initial Responsive Therapy Settings, 85
Recommended Therapy Modifications After Observing Clinical Response, 86
Recording Montage, 38
REPORT TYPES, 33
Request rejected message, 103
RESERVE SPACE column, 47
Responsive Therapy, 75
RESPONSIVE THERAPY trigger, 47
REVERT button, 62
REVIEW & PROGRAM, 43, 58
Review of Initial ECoG Capture Settings, 49
REVIEW screen, 43, 84
RFID Interference—warning, 7
RHYTHMIC ACTIVITY, 57, 61, 67
Rhythmic Activity Detector, 65
RNS® Neurostimulator Programming, 10
RNS® Tablet Instructions, 10
RNS® Tablet Kit Model 5001, 4
S/N, serial number, 23
SATURATION trigger, 47
Save Advanced Settings, 83
Save and Program, 58, 61
SAVE DETECTION SET dialog, 61
SCHEDULED STORAGE trigger, 46
Security and Electronic Tracking Systems—warning, 7
Selecting ECoG Capture Settings, 44
sensitivity, 57
SENSITIVITY (%) slider, 57, 61
serial number of tablet, 23

- Set Capture Duration**, 46
Set Capture Triggers, 46
Set Gain, 41
SHARE button, 70
signal level, 19
signal quality, 19
Significance of Battery Voltage, 27
Significance of Impedance, 27
SIMULATE button, 55, 59
Specifications, 104
SPIKE ACTIVITY, 57, 61, 67
Spike Activity Detector, 66
Start screen, 18, 72
Steps to Connect to a Wi-Fi Network, 12
STIM 1 through **STIM 5**, 77
STIMULATION button, 77
STIMULATION screen, 77, 95
Stimulation testing, 83
STOP/STORE, 37
SYNCH CHANNEL, 82
SYNCH MIN AMPLITUDE, 82
SYNCH SLOPE, 82
SYNCHRONIZATION, 82
Synchronized Therapy, 82
Tablet Does Not Turn On, 103
Tablet Freezes, 103
Tablet Instructions, 10
Tablet login, 14
Tablet Shuts Down, 103
Tablet software updates, 15
Tablet Specifications, 104
Tablet Wireless Specifications, 105
TECHNICAL PARAMETERS, 69
Test Stimulation Before Enabling, 83
TEST STIMULATION screen, 83
Therapy Configuration, 76
THERAPY FREQUENCY MULTIPLIER, 81
Therapy limit per day, 79
THERAPY LIMIT PER DAY, 79
THERAPY LIMIT RESET TIME, 79
Therapy Modifications, 86, 89
Therapy or Stimulation, 75
Therapy Settings—initial, 85
three phases of pattern detection, 50
TIME, 41
TRIGGERS, 35, 46
Troubleshooting, 98
Troubleshooting During a Routine Follow-up, 100
Troubleshooting During an Intraoperative Procedure, 98
Troubleshooting Neurostimulator to Tablet Communication, 102
Troubleshooting, ECoGs, 101
Typographic conventions, 5
Update Center, 15
Update Options, 15
Update Types, 15
USB port, 18
Using the Define Patterns Dialog, 56
wand, 18
Warnings and Cautions – Programmer, 6
Wi-Fi connection, 11, 28, 34, 54, 98
WIRELESS SIGNAL button, 22
X-ray, 100
ZOOM in, 36

NeuroPace, Inc.
455 N. Bernardo Ave.
Mountain View, CA 94043

Customer Support: 1-866-726-3876
(Toll Free in the US)
Fax: 650-237-2855
Website: www.NeuroPace.com

